

CE 0123



## USER MANUAL

# Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators

Rev. 6.0  
Date revised 20.09.2022  
Rx only



Esco Medical Technologies, UAB

Draugystes g. 19 • Kaunas, Lithuania

Tel +370 37 470 000

www.esco-medical.com • support-medical@escolifesciences.com

For Technical Service, contact  
Europe

Esco Medical Technologies, UAB

Draugystes g. 19 • Kaunas, Lithuania

Tel +370 37 470 000

www.esco-medical.com • support-medical@escolifesciences.com

North America

Esco Technologies, Inc.

903 Sheehy Drive, Suite F, Horsham, PA 19044, USA

Tel 215-441-9661 • Fax 484-698-7757

www.escolifesciences.us • eti.admin@escoglobal.com

Rest of the World

Esco Micro Pte. Ltd.

21 Changi South Street 1 • Singapore 486 777

Tel +65 6542 0833 • Fax +65 6542 6920

www.escolifesciences.com • mail@escolifesciences.com

## Copyright Information

© Copyright 2014 Esco Micro Pte Ltd. All rights reserved.

The information in this manual and the accompanying product is copyrighted and all rights are reserved by Esco.

Esco reserves the right to make periodic minor design changes without obligation to notify any person or entity of such change.

Sentinel™ is a registered trademark of Esco.

Caution: Federal law restricts this device to sale by or on a licensed healthcare practitioner's order.

Only to be used by a trained and qualified professional. The device is sold under exemption 21 CFR 801 Subpart D.

*"Material in this manual is provided for informational purposes only. The contents and the product described in this manual (including any appendix, addendum, attachment or inclusion) are subject to change without notice. Esco makes no representations or warranties as to the accuracy of the information contained in this manual. In no event shall Esco be held liable for any damages, direct or consequential, arising out of or related to the use of this manual."*

## **Unpacking and Inspection**

Follow standard receiving practices upon receipt of the instrument. Check the shipping carton for damage. If the damage is found, stop unpacking the instrument. Notify the freight carrier and ask for an agent to be present while the instrument is unpacked. There are no special unpacking instructions, but be careful not to damage the instrument when unpacking it. Inspect the instrument for physical damage such as bent or broken parts, dents, or scratches.

## **Claims**

Our routine method of shipment is via common carrier. If physical damage is found, upon delivery, retain all packing materials in their original condition and immediately contact the carrier to file a claim.

If the instrument is delivered in good physical condition but does not operate within specifications, or if there are any other problems not caused by shipping damage, please contact your local sales representative or Esco Medical immediately.

## **Standard Terms and Conditions**

### **Refunds & Credits**

Please note only serialized products (products labeled with a distinct serial number) and accessories are eligible for a partial refund and/or credit. Non-serialized parts and accessory items (cables, carrying cases, auxiliary modules, etc.) are not eligible for return or refund. In order to receive a partial refund/credit, the product must not have been damaged. It must be returned complete (meaning all manuals, cables, accessories, etc.) within 30 days of original purchase, in "as new" and resalable condition. The *Return Procedure* must be followed.

### **Return Procedure**

Every product returned for refund/credit must be accompanied by a Return Material Authorization (RMA) number obtained from Esco Medical Customer Service. All items being returned must be sent *prepaid* (freight, duty, brokerage and taxes) to our factory location.

### **Restocking Charges**

Products returned within 30 days of original purchase are subject to a minimum restocking fee of 20% of the list price. Additional charges for damage and/or missing parts and accessories will be applied to all returns. Products that are not in "as new" and resalable condition are not eligible for credit return and will be returned to the customer at their own expense.

### **Certification**

This instrument has been thoroughly tested/inspected and found to meet Esco Medical's manufacturing specifications when shipped from the factory. Calibration measurements and testing are traceable and done according to Esco Medical's ISO certification.

### **Warranty and Product Support**

Esco Medical warrants this instrument to be free from defects in materials and workmanship under regular use and service for two (2) years from the original purchase date, provided the instrument is calibrated and maintained following this manual. During the warranty period, Esco Medical will, at our option, either repair or replace a product

that proves to be defective at no charge, provided you return the product (shipping, duty, brokerage and taxes prepaid) to Esco Medical. Any transportation charges incurred are the responsibility of the purchaser and are not included within this warranty. This warranty extends only to the original purchaser. It does not cover damage from abuse, neglect, accident or misuse, or as the result of service or modification by parties other than Esco Medical.

IN NO EVENT SHALL ESCO MEDICAL LTD. BE LIABLE FOR CONSEQUENTIAL DAMAGES. No warranty shall apply when any of the following causes damage:

- Power failure, surges, or spikes
- Damage in transit or when moving the instrument
- An improper power supply such as low voltage, incorrect voltage, defective wiring or inadequate fuses
- Accident, alteration, abuse or misuse of the instrument
- Fire, water damage, theft, war, riot, hostility, acts of God such as hurricanes, floods, etc.

Only serialized products (those items bearing a distinct serial number tag) and their accessory items are covered under this warranty.

PHYSICAL DAMAGE CAUSED BY MISUSE OR PHYSICAL ABUSE IS NOT COVERED UNDER THE WARRANTY. Items such as cables and non-serialized modules are not covered under this warranty.

This warranty gives you specific legal rights and you may have other rights, which vary from province to province, state to state, or country to country. This warranty is limited to repairing the instrument per Esco Medical's specifications.

When you return an instrument to Esco Medical for service, repair or calibration, we recommend shipment using the original shipping foam and container. If the original packing materials are not available, we recommend the following guide for repackaging:

- Use a double-walled carton of sufficient strength for the weight being shipped
- Use heavy paper or cardboard to protect all instrument surfaces. Use non-abrasive material around all projecting parts
- Use at least four inches of tightly packed, industrial-approved, shock-absorbent material all around the instrument

Esco Medical will not be responsible for lost shipments or instruments received in damaged condition due to improper packaging or handling. All warranty claim shipments must be made on a prepaid basis (freight, duty, brokerage, and taxes). No returns will be accepted without a Return Materials Authorization ("RMA") number. Please contact Esco Medical to obtain an RMA number and receive help with shipping/customs documentation.

Re-calibration of instruments, which have a recommended annual calibration frequency, is not covered under warranty.

### **Warranty Disclaimer**

If your instrument is serviced and/or calibrated by someone other than Esco Medical Ltd. and their representatives, please be advised that the original warranty covering your product becomes void when the tamper-resistant Quality Seal is removed or broken without proper factory authorization.

In all cases, breaking the tamper-resistant Quality Seal should be avoided at all cost, as this seal is key to your original instrument warranty. In an event where the seal must be broken to gain internal access to the instrument, you must first contact Esco Medical Ltd.

You will be required to provide us with the serial number for your instrument, as well as a valid reason for breaking the Quality Seal. You should break this seal only after you have received factory authorization. Do not break the Quality Seal before you have contacted us! Following these steps will help ensure that you will retain the original warranty on your instrument without interruption.

### **WARNING**

Unauthorized user modifications or applications beyond the published specifications may result in an electrical shock hazard or improper operation. Esco Medical will not be responsible for any injury sustained due to unauthorized equipment modifications.

ESCO MEDICAL LTD. DISCLAIMS ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE OR APPLICATION.

THIS PRODUCT CONTAINS NO USER-SERVICEABLE COMPONENTS.

UNAUTHORIZED REMOVAL OF THE INSTRUMENT COVER SHALL VOID THIS AND ALL OTHER EXPRESSED OR IMPLIED WARRANTIES.

# Table of contents

1 How to use this manual.....	11
2 Safety warning.....	11
3 Indication for use.....	12
4 About the product.....	12
5 Transport, Storage and Disposal .....	14
5.1 Transportation requirements .....	14
5.2 Storage and operation environment requirements.....	14
5.2.1 Storage requirements .....	14
5.2.2 Operation environment requirements.....	15
5.3 Disposal.....	15
6 Supplied service parts and accessories .....	15
7 Safety symbols and labels .....	16
8 Important safety instructions and warnings.....	18
8.1 Before installation .....	18
8.2 During installation.....	18
8.3 Post installation .....	19
9 Getting started.....	20
10 Mains connection.....	20
11 Gas connections.....	21
12 HEPA / VOC filter (applicable only for the Mini MIRI® Dry model).....	22
12.1 Installation of new filter capsule.....	23
13 User interface .....	24
13.1 Activating the heat and gas controls.....	25
13.2 System menu .....	25
13.3 Status.....	26
13.4 Main menu.....	26
13.4.1 Temperature sub-menu .....	27
13.4.2 CO <sub>2</sub> sub-menu .....	28
13.4.3 O <sub>2</sub> sub-menu.....	30
13.4.4 UVC light sub-menu (applicable only for the Mini MIRI® model).....	32
13.4.5 Service sub-menu.....	32

14 Installation with premixed gas.....	33
14.1 Installation procedure at the site.....	33
14.2 User training.....	36
15 Alarms .....	36
15.1 Temperature alarms .....	36
15.2 Gas level alarms .....	37
15.2.1 CO <sub>2</sub> alarms .....	37
15.2.2 O <sub>2</sub> alarms .....	38
15.3 Gas pressure alarms.....	39
15.3.1 CO <sub>2</sub> pressure alarm .....	39
15.3.2 N <sub>2</sub> pressure alarm.....	39
15.4 Multiple alarms .....	40
15.5 Alarm UVC light (applicable only for the Mini MIRI® Dry model) .....	40
15.6 Loss of power alarm.....	40
15.7 Summary of the alarms.....	41
15.8 Alarm verification .....	42
16 Changing the set points .....	42
16.1 The temperature set point.....	42
16.2 The CO <sub>2</sub> gas concentration set point.....	43
16.3 The O <sub>2</sub> gas concentration set point.....	43
16.4 The culture mode .....	44
17 Surface temperatures and measuring temperature.....	45
18 Pressure.....	47
18.1 CO <sub>2</sub> gas pressure.....	47
18.2 N <sub>2</sub> gas pressure .....	47
19 Firmware.....	48
20 Safe sense function.....	48
21 Cleaning instructions .....	48
21.1 Considerations about a sterile device .....	48
21.2 Manufacturer recommended cleaning procedure.....	49
21.3 Manufacturer recommended disinfection procedure.....	50
22 Heat optimization plates.....	50
23 Humidification .....	51

23.1 Mini MIRI® Dry multiroom IVF incubator .....	51
23.2 Mini MIRI® Humidity multiroom IVF incubator.....	52
24 Temperature validation .....	52
25 Gas level validation .....	53
26 Alarm switch for an external system .....	53
27 Writing area on the compartment lids .....	55
28 Maintenance .....	55
29 Emergency Procedures.....	56
30 User Troubleshooting .....	58
31 Specifications.....	60
32 Electromagnetic compatibility .....	61
33 The Validation guide.....	64
33.1 Product release criteria .....	64
33.1.1 Performance.....	64
33.1.2 Electrical safety.....	64
33.1.3 Communication & data logging .....	64
33.1.4 Gas concentration levels and consumption .....	64
33.1.5 Cosmetic .....	65
34 Validation on-site.....	65
34.1 Mandatory equipment.....	66
34.2 Recommended additional equipment.....	66
35 Testing .....	66
35.1 Gas supply CO <sub>2</sub> .....	66
35.1.1 About CO <sub>2</sub> .....	67
35.2 Gas supply N <sub>2</sub> .....	68
35.2.1 About N <sub>2</sub> .....	68
35.3 CO <sub>2</sub> gas pressure check .....	69
35.4 N <sub>2</sub> gas pressure check.....	69
35.5 Voltage supply .....	69
35.6 CO <sub>2</sub> gas concentration check.....	70
35.7 O <sub>2</sub> gas concentration check .....	70
35.8 Temperature check: Compartment bottoms .....	71
35.9 Temperature check: Compartment lids.....	72



35.10 6-hour stability test .....	72
35.11 Cleaning .....	73
35.12 Test documentation form .....	74
35.13 Recommended additional testing .....	74
35.13.1 A VOC meter (applicable only for the Mini MIRI® Dry model) .....	74
35.13.2 A laser particle counter .....	74
36 Clinical use .....	74
36.1 Temperature check .....	75
36.2 CO <sub>2</sub> gas concentration check .....	75
36.3 O <sub>2</sub> gas concentration check .....	76
36.4 CO <sub>2</sub> gas pressure check .....	76
36.5 N <sub>2</sub> gas pressure check .....	77
37 The Maintenance guide .....	77
37.1 VOC/HEPA filter capsule (applicable only for the Mini MIRI® model) .....	78
37.2 Humidification bottle (applicable only for the Mini MIRI® Humidity model) .....	79
37.3 In-line HEPA filter for CO <sub>2</sub> gas .....	79
37.4 In-line HEPA filter for N <sub>2</sub> gas .....	80
37.5 O <sub>2</sub> sensor .....	80
37.6 CO <sub>2</sub> sensor .....	81
37.7 UV light (applicable only for the Mini MIRI® model) .....	81
37.8 Cooling fan .....	82
37.9 Internal gas pump (applicable only for the Mini MIRI® model) .....	82
37.10 Pump module (applicable only for the Mini MIRI® Humidity model) .....	83
37.11 Proportional valves .....	83
37.12 Gas lines .....	84
37.13 Flow sensors .....	84
37.14 Pressure regulators .....	85
37.15 Internal 0.2µ filter for CO <sub>2</sub> gas .....	85
37.16 Internal 0.2µ filter for N <sub>2</sub> gas .....	85
37.17 Firmware update .....	86
38 The Installation guide .....	86
38.1 Responsibilities .....	86
38.2 Before installation .....	87

38.3 Preparing for installation .....	87
38.4 Bring the following to the installation site .....	87
38.5 Installation procedure at the site .....	88
38.6 User training.....	88
38.7 After the installation .....	89

## 1 How to use this manual

The manual is designed to be read by sections and not ideally from cover to cover. It means that if the manual is read from start to finish, there will be some repetition and overlap. We recommend the following method for going through the manual: first, familiarize yourself with the safety instructions; then, proceed to the essential user functions that are needed for operating the equipment on a day-to-day basis; then, review the alarm functions. The menu functions of the user interface detail information that is required only for advanced users. All parts must be read before the device is taken into use. The Validation guide is detailed described in sections 33 – 36. The Maintenance guide is detailed described in section 37. The Installation procedures are detailed described in section 38.

## 2 Safety warning

- Anyone working with, on or around this equipment should read this manual. Failure to read, understand and follow the instructions given in this documentation may result in damage to the unit, injury to the operating personnel and/or poor equipment performance
- Any internal adjustment, modification or maintenance to this equipment must be undertaken by qualified service personnel.
- If the equipment must be relocated, make sure it is appropriately fixed on a support stand or base and move it on a flat surface. When necessary, move the equipment and the support stand/base separately.
- The use of any hazardous materials in this equipment must be monitored by an industrial hygienist, safety officer or other suitably qualified individuals.
- Before you proceed, you should thoroughly understand the installation procedures and note the environmental/electrical requirements.
- In this manual, important safety-related points will be marked with the following symbols:



### **NOTE**

Used to direct attention to a specific item.



### **WARNING**

Use caution.

- If the equipment is used in a manner not specified by this manual, the protection provided by this equipment may be impaired.

### 3 Indication for use

The Esco Medical MIRI® family's multiroom IVF incubators are intended to be used to provide a stable culture environment at or near body temperature and CO<sub>2</sub>/N<sub>2</sub> or premixed gases and humidification for the development of gametes and embryos during in vitro fertilization (IVF) / assisted reproduction technology (ART) treatments.

### 4 About the product

The Esco Medical Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are CO<sub>2</sub>/O<sub>2</sub> gas incubators.

Direct warming of the dishes in the chambers gives superior temperature conditions in comparison to conventional multiroom IVF incubators.

The temperature in the compartment will remain stable up to 1 °C (even when a lid is open for 30s) and will recover within 1 min after the lid is closed.

The Esco Medical Mini MIRI® multiroom IVF incubators have 2 completely separate culture heat chambers. Each chamber has its own heated lid and warming plate for Petri dish. Mini MIRI® Dry and Mini MIRI® Humidity capacity for 35 mm Petri dish is 16 pcs. And for 60 mm and 4-well Petri dishes – 8 pcs.

To ensure maximum performance, the system of Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators have 4 separate PID temperature controllers. They control and regulate temperature in culture chambers and lids. Compartments do not affect each other's temperatures in any way. The top and the bottom of each compartment is separated with a PET layer so that the lid temperature would not affect the bottom. For validation purposes, each compartment has a PT-1000 sensor built in. The circuitry is separated from the unit's electronics so it remains a truly separate validation system.

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator needs 100% CO<sub>2</sub> and 100% N<sub>2</sub> or premixed gas (for instance 5% CO<sub>2</sub>; 5% O<sub>2</sub> and 90% N<sub>2</sub>) in order to be able to control the CO<sub>2</sub> and O<sub>2</sub> concentrations in the culture chambers.

The multiroom IVF incubator needs 100% CO<sub>2</sub> and 100% N<sub>2</sub> or premixed gas (for instance, 5% CO<sub>2</sub>; 5% O<sub>2</sub> and 90% N<sub>2</sub>) in order to be able to control CO<sub>2</sub> and O<sub>2</sub> concentrations in the culture chambers.

A dual beam infrared CO<sub>2</sub> sensor with extremely low drift rates controls the CO<sub>2</sub> level. A chemical medical grade oxygen sensor controls the level of O<sub>2</sub>.

Gas recovery time is less than 3 min. after opening the lid. To validate gas concentration, the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator is fitted with 2 gas sample ports that allow the user to sample gas from the individual compartment.

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator features a recirculated gas system where gas is continuously put into the compartment and taken out at the same rate. Gas is cleaned via 254 nm UVC light with direct gas contact between the bulb and gas, then through a VOC filter and through a HEPA filter. The UVC light has filters that inhibit any 185 nm radiation that would produce dangerous ozone. The VOC filter is located under the UVC light.

UVC light modules and HEPA-VOC filters are not applied in Mini MIRI® Humidity multiroom IVF incubator.

Complete gas repletion in the system takes less than 5 min.

The total gas consumption is very low. Less than 2 l/h CO<sub>2</sub> and 5 l/h N<sub>2</sub> in use.

For safety reasons the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator has a very complete gas control system that consists of: pressure regulator (preventing dangerous gas pressure problems), gas flow sensors (actual consumption can be accumulated), gas pressure sensors (then user knows that the pressure and variation can be logged to avoid dangerous conditions), gas filters (to avoid valve problems).

Petri dish location in a compartment is easy to reach and safe because of the compartment numbering and the ability to write on the white lid with a pen.

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator has been primarily developed and designed for incubation of gametes and embryos with an overlay of either Paraffin or mineral oil.

 **Refer to section “16.4 The culture mode” for more detailed information.**

The upright LED display in Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators is large, clear and easy to read from a distance. The user can tell if the parameters are correct without going near the unit.

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator can be connected to a PC running the Esco Medical Data logger software for long term data logging and data storage.

MIRI® family's multiroom IVF incubators are stationary devices. The term refers to equipment that, once installed and placed into service, is not intended to be moved from one place to another.

The device is manufactured under a full EU certified 13485 ISO quality management system.

This product fulfils the requirements of EN60601-1 3rd edition standards as a Class I type B equivalent device suited for continuous operation. It also conforms to the requirements of the Regulation (EU) 2017/745 concerning medical devices and is classified as a Class IIa device under rule II.

Personal Protective Equipment (89/686/EEC) and Machine Directive (2006/42/EC) is not applicable for the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators. Also, the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators do not contain or incorporate: a medical substance, including a human blood or plasma derivate; tissues or cells, or their derivatives, of human origin; or tissues or cells of animal origin, or their derivatives, as referred to in Regulation (EU) NO. 722/2012.

## 5 Transport, Storage and Disposal

### 5.1 Transportation requirements

The device is packed in a carton box, and it is wrapped in polyethylene. The box is affixed to a pallet with special straps.

A visual inspection should be done if there is any damage. If no damage is found, the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator can be prepared for transport.

These labels should be glued on the box:

- Label with the marked packing date
- Label with the product name and serial number

### 5.2 Storage and operation environment requirements

#### 5.2.1 Storage requirements

The device may only be store under the following conditions:

- The unit can be in storage for one year. If stored longer than one year, the unit must be returned to the manufacturer for a new release test
- The unit can be stored at temperatures between -20 °C and + 50 °C
- Keep away from direct sunlight

- Caution: consult the accompanying documents for important safety-related information such as warnings and precautions that cannot be presented on the device itself for various reasons
- Do not use if the packing material is damaged
- Keep dry

### 5.2.2 Operation environment requirements

The device may only be used under the following conditions:

- Operating humidity: 5 – 95% RH (Non-Condensing)
- Operating altitude – up to 2000 meters (6560 feet or 80kPa – 106kPa)
- Non-operating altitude – more than 2000 meters (6560 feet or more than 80kPa – 106kPa)
- Environmental temperature: 18 – 30 °C
- Away from direct sunlight
- Kept dry
- For indoor use only

### 5.3 Disposal

Information on handling of the unit as per the WEEE Directive (Waste Electrical and Electronic Equipment).



**The device may have been used for treating and processing infectious substances. Therefore, the device and device components may be contaminated. Prior to disposal, the whole device must be disinfected or decontaminated.**

The unit contains reusable materials. All components (except for the VOC/HEPA and HEPA filters) can be discarded as electrical waste after cleaning and disinfection.

Please note that the VOC/HEPA and HEPA filters must be discarded following the applicable national regulations for special solid waste.

## 6 Supplied service parts and accessories

### Service parts:

- 1 VOC/HEPA filter capsule (only for Mini MIRI® Dry model)
- 1 humidity bottle (only for Mini MIRI® Humidity model)
- 2 HEPA filter for input gas supply
- 2 warming blocks
- 4 warranty labels

- 1 pump box calibration tool (only for Mini MIRI® Humidity model)
- 1 USB stick containing Esco Medical Data logger software and a PDF version of the user manual
- 1 medical grade power cord
- 1 3,5 mm external alarm jack connector
- 1 set of fast male connectors with 15 silicone pipes

Accessories do not apply with the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators.

## 7 Safety symbols and labels

There are several user labels on the surfaces of Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators to guide the user. User labels are shown below.

**Table 7.1** Packing box and electrical safety labels

Description	Image
<p><b>Packing box label:</b></p> <ol style="list-style-type: none"> <li>1. If stored longer than the shelf life, the unit must be returned to the manufacturer for a new release test.</li> <li>2. Shipping temperature between -20 °C and +50 °C.</li> <li>3. Keep away from direct sunlight.</li> <li>4. Caution: consult the accompanying documents for important safety-related information such as warnings and precautions that cannot, for a variety of reasons, be presented on the device itself.</li> <li>5. Consult instructions for proper use of the device.</li> <li>6. Do not use if the packing material is damaged.</li> <li>7. Rx Only.</li> <li>8. Keep dry.</li> </ol>	<p>The diagram shows a rectangular label with the following elements:         <ul style="list-style-type: none"> <li>1. A line pointing to the top text: "This device must be installed before _____ else contact the manufacturer"</li> <li>2. A thermometer icon with "-20°C" at the bottom and "+50°C" at the top.</li> <li>3. A sun icon with rays.</li> <li>4. A warning triangle with an exclamation mark.</li> <li>5. An information icon (book with 'i').</li> <li>6. A circular arrow icon.</li> <li>7. The text "Rx Only" with a circle around the 'x'.</li> <li>8. An umbrella icon with raindrops.</li> </ul> </p>
<ol style="list-style-type: none"> <li>1. View the instructions for use.</li> <li>2. Warning on the back of the device indicates that an earth connection is needed, as well as the mains information, and "ON/OFF" push button.</li> <li>3. "Lightning bolt" indicates the potential risk of electrical shock (never remove any cover).</li> </ol>	<p>The diagram shows a rectangular label with the following elements:         <ul style="list-style-type: none"> <li>1. A line pointing to the title "Operating Instructions" next to an information icon.</li> <li>2. A line pointing to the text "Warning: Equipment must be earthed".</li> <li>3. A line pointing to a lightning bolt icon, with the text "230V~, 50Hz, 160W" to its right.</li> </ul> <p>At the bottom of the label, it says: "Fuses: 2x3.15A-250V – CB 1.5KA"</p> </p>






**Table 7.2** Device label

Description	Image
<ol style="list-style-type: none"> <li>Model.</li> <li>Main's power rating.</li> <li>CE mark.</li> <li>Not protected against ingress of water.</li> <li>Manufacturer's address and country of origin.</li> <li>View instruction for use.</li> <li>Temperature limit.</li> <li>Rx only.</li> <li>Serial number.</li> <li>UDI-DI code.</li> <li>Logo.</li> <li>Keep away from direct sunlight.</li> <li>Observe WEEE.</li> <li>Keep dry.</li> <li>Year of manufacture.</li> </ol>	<p>The image shows two identical device labels for 'MULTIROOM IVF INCUBATOR'. The top label is for 'Mini MIRI® Dry' (SN: 2692) and the bottom label is for 'Mini MIRI® Humidity' (SN: 2758). Both labels include technical specifications, safety warnings, and a barcode.</p> <p><b>Label 1 (Dry):</b></p> <ul style="list-style-type: none"> <li>MODEL: Mini MIRI® Dry</li> <li>MAINS: ~230V, 50Hz, 160W</li> <li>Esco Medical Technologies, UAB, Draugystes g.19, Kaunas 51230 Lithuania</li> <li>CE 0123</li> <li>IPXO</li> <li>Consult instruction for use</li> <li>Temperature limit: 18°C to 30°C</li> <li>Rx only</li> <li>SN: 2692</li> <li>Keep away from direct sunlight</li> <li>Not for general waste</li> <li>Keep dry</li> <li>2022</li> <li>Barcode: (01)4779041940151(11)220919(21)2692</li> </ul> <p><b>Label 2 (Humidity):</b></p> <ul style="list-style-type: none"> <li>MODEL: Mini MIRI® Humidity</li> <li>MAINS: ~230V, 50Hz, 160W</li> <li>Esco Medical Technologies, UAB, Draugystes g.19, Kaunas 51230 Lithuania</li> <li>CE 0123</li> <li>IPXO</li> <li>Consult instruction for use</li> <li>Temperature limit: 18°C to 30°C</li> <li>Rx only</li> <li>SN: 2758</li> <li>Keep away from direct sunlight</li> <li>Not for general waste</li> <li>Keep dry</li> <li>2022</li> <li>Barcode: (01)4779041940162(11)2201019(21)2758</li> </ul>

**Table 7.2** Info labels on Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators

Description	Image
USB communication port	USB communication port
CO <sub>2</sub> inlet	CO <sub>2</sub> 100% Inlet
N <sub>2</sub> inlet	N <sub>2</sub> 100% Inlet
Alarm port	Alarm port
Compartments numbers are indicated in the top corner of the lid with a label	1 2

Description	Image
Humidification bottle <sup>1</sup>	
VOC/HEPA filter <sup>2</sup>	

 **The connected external device to signal input/output connections should be compliant with the appropriate safety standard for medical equipment EN 60601-1. It applies to USB connection.**

Compartment numbers are shown in the picture below and also indicated on the top of the lid with labels:



**Figure 7.1** Compartment numbers on Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators

## 8 Important safety instructions and warnings

### 8.1 Before installation

1. Do not use the product if the package is damaged. Contact Esco Medical or the local representative.
2. Read the user manual thoroughly before use.
3. Always keep these instructions easily accessible near the device.

### 8.2 During installation

1. Never place this unit on top of other equipment that gives off heat.
2. Place this unit on a flat, hard and stable surface.

<sup>1</sup> Only for Mini MIRI® Humidity model

<sup>2</sup> Only for Mini MIRI® Dry model

3. Do not place the unit on a carpet or similar surfaces.
4. Do not defeat the safety purpose of the grounding-type (earthing) plug.
5. A grounding-type (earthing) plug with two blades and a third prong are provided for your safety. If the provided plug does not fit into your outlet, consult an electrician to replace the outlet.
6. Always connect the power cord to a properly grounded outlet and only use the cord that came with the device.
7. Do not install near any heat sources such as radiators, heat registers, stoves or other apparatus that produce heat.
8. Do not use this device near water sources.
9. Use only 100% concentration CO<sub>2</sub> and 100% concentration N<sub>2</sub> gases.
10. Always use an external HEPA filter for input CO<sub>2</sub> and N<sub>2</sub> gases.
11. Do not use this product if the room temperature exceeds 30 °C.
12. Place this unit in a location with adequate ventilation to prevent internal heat build-up. Leave at least 10 cm clearance from the rear, 30 cm from the top and 20 cm from left and right to prevent overheating and allow access to the ON/OFF switch in the back.
13. This unit is intended for indoor purposes only.
14. The unit must be connected to a suitable uninterrupted power supply (UPS) source.

### 8.3 Post installation

1. Refer all servicing procedures to qualified service personnel.
2. Servicing is required according to the service manual as well as cases when the device has been damaged in any way, e. g. suppose the apparatus has been dropped, exposed to rain or moisture or does not operate normally. The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators contain high voltage components that may be hazardous.
3. Unplug this device during lightning storms or when unused for an extended period of time.
4. Protect the power cord from being walked on or pinched, particularly at the plug, convenience receptacles and the point where it exits from the apparatus.
5. Perform temperature and gas calibration at the intervals described in the manuals.
6. Never leave the lids open for more than 10 sec while in use.
7. VOC/HEPA filters must be changed every 3 months (do not apply for Mini MIRI® Humidity multiroom IVF incubator).
8. A maintenance plan must be fulfilled to keep the device safe.
9. NEVER block gas supply holes in the compartment.
10. Ensure that CO<sub>2</sub> and N<sub>2</sub> gas supply pressures are kept stable at 0.4 – 0.6 bar (5.80 – 8.70 PSI).

11. Never use any other except Esco Medical filter. Otherwise, the warranty will be void.
12. Do not use the device without a proper Esco Medical VOC/HEPA filter attached (do not apply for Mini MIRI® Humidity multiroom IVF incubator).

## 9 Getting started



**The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators must be installed by authorized and trained personnel only!**

1. Follow the guidelines in the safety instructions and warnings section.
2. Connect the mains cable to the UPS.
3. Connect the mains cable to the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator.
4. Connect gas lines.
5. Set the gas pressure on the external gas regulator at 0.4 – 0.6 bar (5.80 – 8.70 PSI).
6. Switch on the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator in the back.
7. Observe for standard functionality.
8. Let the unit warm up and stabilize for 20 min.
9. Follow the guidelines in the Validation guide (see section “33 The validation guide”)
10. Complete user training and finish reading instructions.
11. After a burn-in phase of 24 hours, the unit is ready for use IF the testing is successful.



**Clean and disinfect the device before use. It is not delivered sterile or in a clinically acceptable cleanliness state. Consult the cleaning instructions section in this manual for the manufacturer's recommended guidelines!**

## 10 Mains connection

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators come with a detachable mains power cord. The power cord is prepared for the country in which the unit is intended to be used.

The ON/OFF switch provides the user with a means to isolate the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator from the mains.



**Do not defeat the safety purpose of the grounding-type plug! A grounding-type plug has two blades and a prong, which is provided for your safety. If the provided plug does not fit into your outlet, consult an electrician to replace the outlet.**

The power requirement is 230V 50Hz OR 115V 60Hz. The built-in power supply has a switch mode that automatically adjusts to the correct mains power between 100V-240V AC 50-60 Hz.



**Figure 10.1** Power supply

## 11 Gas connections

There are two gas inlets on the back of the unit. These ports are marked "CO<sub>2</sub> 100% Inlet" and "N<sub>2</sub> 100% Inlet".



**Figure 11.1** Gas inlets on the back of Mini MIRI® Dry and MIRI® Humidity multiroom IVF incubators

CO<sub>2</sub> inlet should be connected to a 100% concentration of CO<sub>2</sub>. CO<sub>2</sub> control in the compartment is available in the range from 1.9% to 9.9%.

The N<sub>2</sub> inlet should be connected to 100% concentration N<sub>2</sub> if low oxygen conditions are required. The O<sub>2</sub> control in the compartments is available in the range from 3.9% - 19.9% by infusing N<sub>2</sub>.

The premixed gas inlet should be connected to the CO<sub>2</sub> inlet.



**The inlet's gas pressure should be between 0.4 – 0.6 bar (5.80 – 8.70 PSI) and it must be kept stable!**

Always use a high-quality pressure regulator that can be set with the required precision for both gases.



**Figure 11.2** Pressure regulator


Connect the CO<sub>2</sub> gas to the CO<sub>2</sub> inlet with a suitable silicone tube. Ensure that the tube is fastened with a clip so that it does not accidentally loosen itself during sudden pressure fluctuation. Use the supplied 0.2μ HEPA filter on the gas line just before the inlet on the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator. Notice the direction.

Connect the N<sub>2</sub> inlet to the Nitrogen Bottle in a similar way.



**Figure 11.3** Gas filter

Connect the N<sub>2</sub> inlet to the Nitrogen Bottle in a similar way.

 **The Mini MIRI® and Mini MIRI® Humidity multiroom IVF incubators can also run-on premixed gas. It is a more expensive option for gas consumption. It also means the user cannot adjust the CO<sub>2</sub> and O<sub>2</sub> levels without changing the gas supply. Please read the "13 Installation with premixed gas" section below for more detailed information about using the device on premixed gas.**

## 12 HEPA / VOC filter (applicable only for the Mini MIRI® Dry model)

VOCs are hydrocarbon-based compounds that are found in fuel, solvents, adhesives and other compounds. Examples of VOCs include isopropanol, benzene, hexane, formaldehyde, vinyl chloride.

VOCs can also occur in medical gases, such as CO<sub>2</sub> and N<sub>2</sub>. It is essential to use in-line VOC filters to prevent these fumes from entering your Mini MIRI® Dry multiroom IVF incubators for your medical gasses.

Unexpected sources of VOCs are commonly found in IVF labs. These can include cleaning agents, perfumes, cabinetry, grease on the wheels of equipment and sources in HVAC equipment.

VOCs are typically measured in parts per million (ppm.) They can also be reported in parts per billion (ppb.) For IVF, the recommended count below 0.5 ppm; the total quantity of VOCs should be below <0.2 ppm or preferably zero.

High levels of VOCs (over 1 ppm) are toxic to embryos, resulting in poor embryo development and even probable failure to reach the blastocyst stage.


VOC levels in the 0.5 ppm range will typically allow an acceptable blastocyst development and reasonable pregnancy rates but will likely result in a high percentage of miscarriages.

A combined HEPA and VOC filter (carbon filter) are integrated into the construction of the Mini MIRI® Dry multiroom IVF incubator. Before entering the Mini MIRI® Dry multiroom IVF incubator, the gas is sent through the filter in a single pass. Then, upon return from the compartment, the gas is filtered again. The recirculation system constantly filters gas in the Mini MIRI® Dry multiroom IVF incubator.

The combined HEPA and VOC filter are mounted on the device's back to ease access and replacement.

## 12.1 Installation of new filter capsule

Two blue caps that are installed on the filter can be discarded during unwrapping. Correct filter performance is crucial for system performance.

 **Filter element must be changed every 3 months. Mark the date when it is put on and make sure to keep this interval!**

Start by putting the blue fittings on the filter into the filter holder sockets. The flow arrow on the Mini MIRI® Dry multiroom IVF incubator and the filter should point in the same direction.



**Figure 12.1** The flow arrow on the Mini MIRI® multiroom IVF incubator



**Figure 12.2** The way of pulling filter



**Figure 12.3** Filter in place

Then, simultaneously press both angle fittings (using both hands) into the holes till they snap into place. The last 4 mm step should feel stiff.

**⚠ A filter element that has been fitted incorrectly will cause the unit to work not as intended. This is dangerous!**

The filter is removed by gently pulling it straight out using both hands.

**⚠ Never run the Mini MIRI® Dry multiroom IVF incubator if the filter element is missing! Dangerous particle contamination could occur!**





## 13 User interface

In the following chapters, the functions associated with keys and menu items are going to be explained.

User interface handles daily used functions and more advanced adjustments that might be made to the device. The main keys and their purpose are presented in table 13.1.



**Table 13.1** Main keys and their purpose

Description	Image
<b>Rotary button</b> Use to toggle and select items on the menu to change their status. It is also used to change the temperature and gas setpoints values	
<b>ON/OFF button</b> It is located in the REAR of the unit	
<b>Alarm key</b> It mutes an audible alarm and visually indicates the alarm condition by a flashing red circle of light. The audio alarm will come back on automatically after 5 min. It can be muted again	
<b>Display panel</b> Shows the information on the current status of the unit. The display consists of 7 x high brightness 16 segment LEDs. The first one is red to indicate a user warning. The other 6 are blue and used to display normal running conditions	

### 13.1 Activating the heat and gas controls

Heat and gas controls are activated using the "ON/OFF" switch in the rear.

Soon after system activation, the main display will alternate the reading between the following 4 parameters:

Temperature	= Temperature in °C
CO <sub>2</sub>	= CO <sub>2</sub> concentration in %
O <sub>2</sub>	= O <sub>2</sub> concentration in %
Mode	= Open/Oil Culture

### 13.2 System menu

Press and hold the rotatory button for 3 seconds to access the menu.


Navigate in the menu by:

- Rotating rotatory button clockwise (↻) or anticlockwise (↺) = previous OR next
- Pressing rotatory button = enter, change OR accept

Rotate the rotatory button (↺) to exit the menu entirely.

### 13.3 Status

Alternating between the 4 values under normal running conditions:

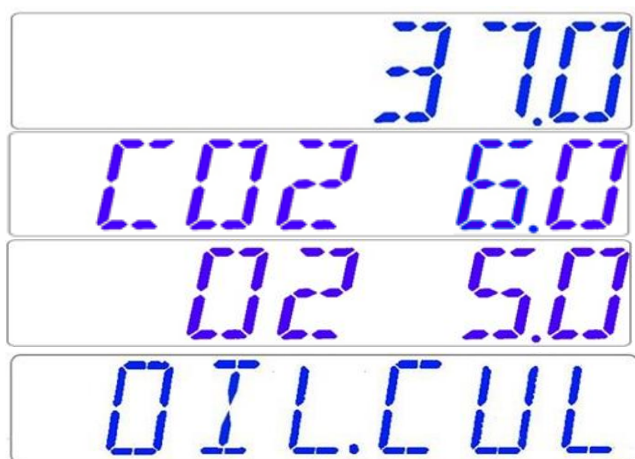

Force scroll between parameters with rotating rotatory button (↻) or (↺).

👉 If the O<sub>2</sub> regulator is deactivated, the system will display “O2 OFF”.



👉 If the use mode is “Open Culture” (no oil or Paraffin overlay culture), the device shall be set for that and will display:



### 13.4 Main menu

Press the rotatory button – enter the menu.

The user can exit the menu by rotating the rotatory button (↺).



Temperature is the first category when the user enters the menu.

Press the rotatory button to enter the Temperature sub-menu.



Rotate the rotatory button (↶) to scroll further down in the menu.  
Press the rotatory button to enter the CO<sub>2</sub> sub-menu.



Rotate the rotatory button (↶) to scroll further down in the menu.  
Press the rotatory button to enter the O<sub>2</sub> sub-menu.



Rotate the rotatory button (↶) to scroll further down in the menu.  
Press the rotatory button to enter the UVC light sub-menu (not available in the Mini MIRI® Humidity multiroom IVF incubator).



Rotate the rotatory button (↶) to scroll further down in the menu.  
Press the rotatory button to enter the Service sub-menu.



#### 13.4.1 Temperature sub-menu

Press the rotatory button on the temperature menu to enter the temperature sub-menu.  
The first item in the temperature sub-menu is T1 sensor calibration:



Calibrate T1 temperature by pressing the rotatory button and rotating it (↶) or (↷) to adjust.

Move to the next sub-menu item with (↶) rotation or one step up with (↷) rotation.

👉 Each compartment has two internal temperature sensors. One in the compartment lid and another in the compartment bottom.

***Example – how to calibrate the temperature:***

The temperature has to be measured with a suitable and calibrated device. With a quality thermometer, it has been estimated that T1 is 37.4 °C. Locate "T1 CAL" in the sub-menu and press and hold the SP key. The display should show:

A digital display showing the text "T1 37.0" in a blue, segmented font. The "T1" is on the left and "37.0" is on the right, with a small decimal point between the 7 and 0.

Rotate rotatory button (↶) to adjust the temperature calibration to the desired level by (↶) or (↷) rotation. The display will show the steps from 37.1, 37.2, 37.3 and 37.4. When temperature equals the measured temperature, press the rotatory button. The new value is stored and the temperature sensor for the T1 area has been modified.

👉 Calibration procedure is the same for T1 – T4.

### 13.4.2 CO<sub>2</sub> sub-menu


Press the rotatory button on the CO<sub>2</sub> menu to enter the CO<sub>2</sub> sub-menu.  
The first item in the CO<sub>2</sub> sub-menu is CO<sub>2</sub> sensor calibration:

A digital display showing the text "CO2CAL" in a blue, segmented font. The "CO2" is on the left and "CAL" is on the right, with a small decimal point between the 2 and C.

Calibrate CO<sub>2</sub> by pressing the rotatory button and rotating it (↶) or (↷) to adjust.  
Move to the next CO<sub>2</sub> sub-menu item with (↶) rotation or one step up with (↷) rotation.

A digital display showing the text "CO2REG" in a blue, segmented font. The "CO2" is on the left and "REG" is on the right, with a small decimal point between the 2 and R.

Toggle CO<sub>2</sub> regulation on/off by pressing the rotatory button and rotating it (↶) or (↷).

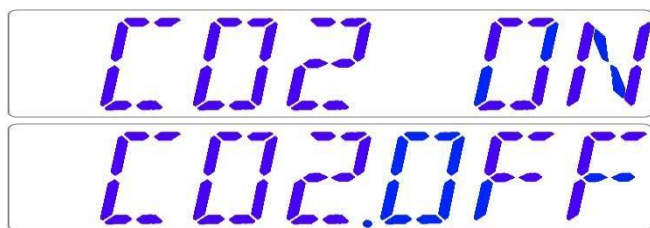

 The default status for the CO<sub>2</sub> control is OFF.

Move to the next CO<sub>2</sub> sub-menu item with (↶) rotation or one step up with (↷) rotation. CO<sub>2</sub> flow rate is shown (it cannot be adjusted):



It shows the amount of CO<sub>2</sub> gas put into the system while regulating. The volume is shown in liters/hour. It usually will fluctuate along with the CO<sub>2</sub> regulation.

Move to the next CO<sub>2</sub> sub-menu item with (↶) rotation or one step up with (↷) rotation. CO<sub>2</sub> internal pressure rate is shown (it cannot be adjusted on the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator. It is adjusted on the external gas regulator):



The value is in bar and it must be 0.4 – 0.6 bar (5.80 – 8.70 PSI) at all times.

#### ***Example – how to calibrate CO<sub>2</sub>:***

CO<sub>2</sub> gas concentration has to be measured with a suitable and calibrated device. The real CO<sub>2</sub> concentration has been estimated to be 6.4% on one of the gas sample ports. Each port is suitable for this purpose.

Locate “CO<sub>2</sub> CAL” in the CO<sub>2</sub> sub-menu and press the rotatory button. The display should show:



Rotate rotator button (↺) to adjust the CO<sub>2</sub> calibration to the desired level by (↺) or (↻) rotation. In this case, we want to change to 6.4%. Rotate (↺) till the display shows 6.0, 6.1, 6.2, 6.3 and 6.4. When CO<sub>2</sub> equals measured CO<sub>2</sub>, press the rotatory button. The new value is stored and the CO<sub>2</sub> sensor calibration has been modified.

👉 **Pure CO<sub>2</sub> 100% gas recovery till 5% is less than 4 minutes.**

👉 **Calibration is performed by adjusting the CO<sub>2</sub> level according to the measurement taken from the gas sampling outlet by an external reliable CO<sub>2</sub> measurement device.**

⚠ **Calibration values should only be changed by a trained user or the technician, according to specific measurements. Done with the calibrated device.**

### 13.4.3 O<sub>2</sub> sub-menu

Press the rotatory button on the O<sub>2</sub> menu to enter the O<sub>2</sub> sub-menu.  
The first item in the O<sub>2</sub> sub-menu is O<sub>2</sub> sensor calibration:



Calibrate O<sub>2</sub> by pressing the rotatory button and rotating it (↺) or (↻) to adjust.  
Move to the next O<sub>2</sub> sub-menu item with (↺) rotation or one step up with (↻) rotation.



Toggle O<sub>2</sub> regulation on/off by pressing the rotatory button and rotating it (↺) or (↻).



👉 **The Default status for the O<sub>2</sub> control is OFF.**

Move to the next O<sub>2</sub> sub-menu item with (↺) rotation or one step up with (↻) rotation.

N<sub>2</sub> flow rate is shown (it cannot be adjusted):



It shows the amount of N<sub>2</sub> gas put into the system while regulating. The volume is shown in liters/hour. It usually will fluctuate along with the O<sub>2</sub> regulation.

Move to the next O<sub>2</sub> sub-menu item with (↶) rotation or one step up with (↷) rotation. N<sub>2</sub> internal pressure rate is shown (it cannot be adjusted on the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator. It is adjusted on the external gas regulator):



The value is in bar and it must be 0.4 – 0.6 bar (5.80 – 8.70 PSI) at all times.

#### ***Example - How to calibrate O<sub>2</sub>:***

O<sub>2</sub> gas concentration has to be measured with a suitable and calibrated device. The real O<sub>2</sub> concentration has been estimated to be 5.3% on one of the gas sample ports. Each port is suitable for this purpose.

Locate “O<sub>2</sub> CAL” in the O<sub>2</sub> sub-menu and press the rotatory button. The display should show:



Rotate rotator button (↶) to adjust the O<sub>2</sub> calibration to the desired level by (↶) or (↷) rotation. In this case, we want to change to 5.3%. Rotate (↶) till the display shows 5.0, 5.1, 5.2 and 5.3. When O<sub>2</sub> equals measured O<sub>2</sub>, press the rotatory button. The new value is stored and the O<sub>2</sub> sensor calibration has been modified.

👉 Calibration is performed by adjusting the O<sub>2</sub> level according to the gas sampling outlet's measurement by an external reliable O<sub>2</sub> measurement device.

👉 Calibration values should only be changed by a trained user or the technician, according to specific measurements. Done with the calibrated device.



#### 13.4.4 UVC light sub-menu (applicable only for the Mini MIRI® model)

Press the rotatory button on the UV-C menu to enter the UV-C light sub-menu.



Toggle UV-C light regulation on/off by pressing the rotatory button and rotating it (↶) or (↷).



 The default status for the UV-C light is “ON”.

The UV light will automatically switch off when the unit is turned off.

 For optimal air cleaning, it is recommended to have the UV-C light set to “ON” when the unit is used.

#### 13.4.5 Service sub-menu

Press the rotatory button on the service menu to enter the service sub-menu.

The display will show the currently installed firmware version:



Ver 2.0 is only shown as an example. Consult Esco Medical or the local representative for the number of the latest version.

Move to the next service sub-menu item with (↶) rotation or one step up with (↷) rotation.



The display will show the “GAS” function:



Press the rotatory button to enter and rotate (↶) or (↷) to choose “PREMIX” or “CO<sub>2</sub>/O<sub>2</sub>”. Press the rotatory button to choose premixed or CO<sub>2</sub>/O<sub>2</sub> gas mode.

**👉 When using the premixed gas mode, it is necessary to use a premixed gas with higher gradation than the setpoint. For example, if you need to achieve a 5% CO<sub>2</sub> gas setpoint, premixed gas should have 6% CO<sub>2</sub> in its mixture.**

Then exit the menu by (↷) or press the rotatory button and hold it until the main menu does not appear.

## 14 Installation with premixed gas

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators have primarily been designed to run on 100% CO<sub>2</sub> and 100% N<sub>2</sub>. It can also run with a premixed gas. Running on 100% CO<sub>2</sub> and 100% N<sub>2</sub> gases, the device accuracy will be significantly higher (< 0.2% from the selected setpoint) compared to using the device on premixed gas. A premixed gas is usually used for simpler incubation systems that do not contain any CO<sub>2</sub> and O<sub>2</sub> sensors and have no gas mixing capabilities.

This section describes how to install the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator at an IVF clinic running with premixed gas.

**👉 The Premixed gas concentration must be chosen specifically to match the requirement of the culture medium. As the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators cannot alter the concentration, the media's resulting pH will depend on the correct concentration choice.**


**👉 Be advised that premixed gas consumption will be significantly higher compared to pure gas. Recover to the setpoint will be longer.**

### 14.1 Installation procedure at the site

Follow all the instructions in the installation manual and the guidelines in the user manual's safety instructions and warnings section.

Instead of connecting Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators to either only 100% CO<sub>2</sub> or both 100% CO<sub>2</sub> and 100% N<sub>2</sub>, the incubator is attached to only a premixed source.

Premixed gas should only be connected to the CO<sub>2</sub> gas port (a 4 mm diameter hose barb).

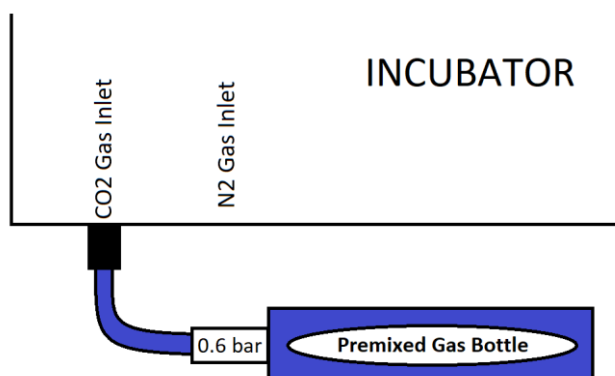

**Figure 14.1** Premixed gas connections to the incubator

👉 Please read the “11 Gas connection” section in this manual above for more detailed gas connection requirements.

Measure the gas concentration from the premixed gas bottle with a calibrated gas analyzer. The result of the measurement is significant for the set-up of the device and the correct operation.

CO<sub>2</sub> regulation must be “ON” in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators menu. CO<sub>2</sub> is generally as a default set to “ON” and O<sub>2</sub> to “OFF”.

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators must be set to premix gas work mode.


**Please follow these instructions:**

Press the rotatory button – enter the menu.

The user can exit the menu by rotating the rotatory button (↺).



Press the rotatory button on the service menu to enter the service sub-menu.

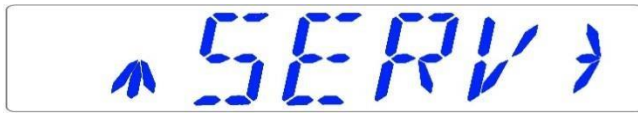


The display will show the currently installed firmware version. Move to the next service sub-menu item with (⏮) rotation or one step up with (⏭) rotation.


The display will show the “GAS” function:




Press the rotatory button to enter the service sub-menu and rotate (⏮) or (⏭) to choose “PREMIX”. Press the rotatory button to select the premixed gas work mode.


Then exit the menu by (⏭) or press the rotatory button and hold it until the main menu does not appear.

 **The CO<sub>2</sub> setpoint must be 0.1% lower than the premixed gas measured value (i.e., 4.9% if 5.0% measured).**

 **The O<sub>2</sub> setpoint must be 1% higher than the premixed gas measured value gas (i.e., 5.0% if 6.0% measured).**

For changing the CO<sub>2</sub> and O<sub>2</sub> setpoints, please read the 16.2 and 16.3 sections in this manual.

 **If the setpoints are not set up correctly, a continual gas flow may occur, which will lead to high gas consumption and incorrect recovery times.**

 **The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators contain a high-grade CO<sub>2</sub> and O<sub>2</sub> sensor. They will measure the incoming premixed gas. Make sure that sensors are reading the anticipated gas percentage in the display of the device. That is a percentage that is the proximity of the values on the certificate of the gas bottle. If this is not the case, it must be established if the bottle's concentration per the certificate is correct. If so, the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators sensors must be calibrated. Refer to the user manual for gas calibration. If the gas bottle does not contain the expected mixture, contact the gas bottle supplier.**

## 14.2 User training

Explain the user:

1. As seen in the display, the gas concentration values must be 0.1% (CO<sub>2</sub> lower and O<sub>2</sub> higher) from the values they expect. If they try to change the setpoint or the calibration to get rid of the offset, the regulation will not work.
2. They cannot set the setpoints they would typically do when using 100% CO<sub>2</sub> and 100% N<sub>2</sub> as the source gas. It is an inherent compromise of using premixed gas. The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators cannot change the gas composition of the premixed gas.
3. If the media's pH is not correct, they must get a new mixture of premixed gas. They cannot adjust anything on the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator.
4. If they change to another concentration, the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators setpoints must be adjusted accordingly, as described above. They should also check the flow rates when they change to a new bottle if it does not precisely contain the same gas mixture.

## 15 Alarms

The display will show a red "A" and the affected parameter's status message on a single fault condition. An audio signal can be muted by pressing the alarm key once (toggled on/off for 5-minutes). There will be a red arrow that indicates if the alarm is triggered due to too high or too low values, and the audio on/off key will blink red:



**Figure 15.1** Alarm key which indicate the alarm condition

The audio pattern is 3 short beeps followed by a 3-second pause. All alarms have the same audio pattern.

### 15.1 Temperature alarms

All 2 compartments can trigger a temperature alarm if their temperature varies over  $\pm 0.5$  °C from the setpoint.

 **Remember that changing the setpoint more than  $\pm 0.5$  °C from the current temperature will result in an alarm. The same goes for all calibration adjustments.**

The number will indicate the zone triggering the alarm following "A".

Temperature is too high in compartment 2:




Temperature is too low in compartment 1:



The display will lock on the alarm condition and will stop alternating between the standard status messages. If the mute key is pressed, the display will shift to normal status and show the parameters for 5 minutes until the audio alarm comes back on again. The mute alarm key will still show the alarm condition by blinking red while the alarm is muted.

The zone layout and sensor placement are described in the section "16 Surface temperatures and measuring temperature".

If a temperature sensor malfunctions, it will be indicated by the following warning:

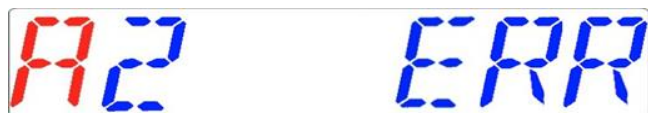


It denotes that the sensor in compartment 2 has failed. As a safety precaution, the heating of the affected area will be switched off.

## 15.2 Gas level alarms

### 15.2.1 CO<sub>2</sub> alarms

The CO<sub>2</sub> gas level alarm is activated if the concentration of the CO<sub>2</sub> gas deviates more than  $\pm 1\%$  from the setpoint.

 **Remember that changing the setpoint more than  $\pm 1\%$  from the current gas level will result in a gas level alarm. The same goes for all calibration adjustments.**

CO<sub>2</sub> gas % is too low:




CO<sub>2</sub> gas % is too high:



The display will lock on the alarm condition and will stop alternating between the standard status messages. If the mute key is pressed, the display will shift to normal status and show the parameters for 5 minutes until the audio alarm comes back on again. The mute alarm key will still show the alarm condition by blinking red while the alarm is muted.

### 15.2.2 O<sub>2</sub> alarms

The O<sub>2</sub> gas level alarm is activated if the concentration of the O<sub>2</sub> gas deviates more than  $\pm 1\%$  from the setpoint.

 **Remember that changing the setpoint more than  $\pm 1\%$  from the current gas level will result in a gas level alarm. The same goes for all calibration adjustments.**

O<sub>2</sub> % is too low:



O<sub>2</sub> % is too high:



The display will lock on the alarm condition and will stop alternating between the standard status messages. If the mute key is pressed, the display will shift to normal status and show the parameters for 5 minutes until the audio alarm comes back on again. The mute alarm key will still show the alarm condition by blinking red while the alarm is muted.

## 15.3 Gas pressure alarms

### 15.3.1 CO<sub>2</sub> pressure alarm

If the CO<sub>2</sub> gas supply is not attached correctly or incorrect CO<sub>2</sub> gas pressure is applied to the system, the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators will go into CO<sub>2</sub> pressure alarm mode. The display will show "CO<sub>2</sub> P", which indicates an incorrect incoming gas pressure. If the pressure falls below 0.3 bar (4.40 PSI) or rises above 0.7 bar (10.20 PSI), it will trigger the alarm.



 "P" stands for pressure.

The display will lock on the alarm condition and will stop alternating between the standard status messages. If the mute key is pressed, the display will shift to normal status and show the parameters for 5 minutes until the audio alarm comes back on again. The mute alarm key will still show the alarm condition by blinking red while the alarm is muted.

### 15.3.2 N<sub>2</sub> pressure alarm

If the N<sub>2</sub> gas supply is not attached correctly or incorrect N<sub>2</sub> gas pressure is applied to the system, the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators will go into N<sub>2</sub> pressure alarm mode. The display will show "N<sub>2</sub> P", which indicates an incorrect incoming gas pressure. If the pressure falls below 0.3 bar (4.40 PSI) or rises above 0.7 bar (10.20 PSI), it will trigger the alarm.



 "P" stands for pressure.

The display will lock on the alarm condition and will stop alternating between the standard status messages. If the mute key is pressed, the display will shift to normal status and show the parameters for 5 minutes until the audio alarm comes back on again. The mute alarm key will still show the alarm condition by blinking red while the alarm is muted.

## 15.4 Multiple alarms

When there are two or more alarms, the display will indicate this by showing first "A MULTI" and then the alarm conditions:



The flow will be forced according to the alarms. The temperature alarms have 1<sup>st</sup> priority, gas level alarms 2<sup>nd</sup>, and gas pressure 3<sup>rd</sup>.

## 15.5 Alarm UVC light (applicable only for the Mini MIRI® Dry model)

Alarms on UV-C light will show only as a warning message during the normal status. A red "S" will appear. **There will be no audio alarm.**



The user should consult the distributor for further guidance or service inspection. Only when the UV-C light works again will the "S" disappear.

## 15.6 Loss of power alarm

If the power is disconnected, the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator will give an audio alarm for approximately 4 seconds, and the LED in the mute alarm key will flash.



**Figure 15.2** Alarm key which indicates the alarm condition



## 15.7 Summary of the alarms

In the table below, there is a list of every possible alarm in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators.

**Table 15.1** Every possible alarm in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators

Alarm name	Conditions	How it is determined	Alarm group	Alarm priority
Low-temperature alarm	If the temperature falls below 0.5 °C from the SP. It is applicable for all compartment's bottom temperature	Each temperature zone sensor reading	Technical	High priority alarm
High-temperature alarm	If the temperature rises above 0.5 °C from the SP. It is applicable for all compartment's bottom temperature		Technical	High priority alarm
Low CO <sub>2</sub> level	When the CO <sub>2</sub> concentration drops by 1% from the SP, after 3 min the alarm will turn on	CO <sub>2</sub> sensor reading	Technical	High priority alarm
High CO <sub>2</sub> level	When the CO <sub>2</sub> concentration rises by 1% from the SP, after 3 min the alarm will turn on		Technical	High priority alarm
Low O <sub>2</sub> level	When the O <sub>2</sub> concentration drops by 1% from the SP, after 5 min the alarm will turn on	O <sub>2</sub> sensor reading	Technical	High priority alarm
High O <sub>2</sub> level	When the O <sub>2</sub> concentration rises by 1% from the SP, after 5 min the alarm will turn on		Technical	High priority alarm
Low incoming CO <sub>2</sub> pressure	If the pressure falls below 0.3 bar	Pressure sensor reading	Technical	High priority alarm
High internal CO <sub>2</sub> pressure	If the pressure rises above 0.7 bar	Pressure sensor reading	Technical	High priority alarm
Low incoming N <sub>2</sub> pressure	If the pressure falls below 0.3 bar	Pressure sensor reading	Technical	High priority alarm
High internal N <sub>2</sub> pressure	If the pressure rises above 0.7 bar	Pressure sensor reading	Technical	High priority alarm
UV alarm	If the UV lamp is malfunctioning	UV sensor reading	Technical	Informative alarm

## 15.8 Alarm verification

In the table below, there is a list of how and when to verify the functionality of the alarm system.

**Table 15.2** Every possible alarm in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators

Alarm name	How to verify an alarm	When to verify an alarm
High-temperature alarm	Decrease the setpoint value by 3.0 °C from the current setpoint	If you have a suspicion that alarms are malfunctioning
Low-temperature alarm	Put cold metal part in the middle of the compartment and close the lid	
High CO <sub>2</sub> level	Decrease the setpoint value by 3.0% from the current setpoint	
Low O <sub>2</sub> level	Increase the setpoint value by 3.0% from the current setpoint	
High O <sub>2</sub> level	Open the lid and leave it open for 5 min	
Low CO <sub>2</sub> level	Open the lid and leave it open for 3 min	
Low incoming CO <sub>2</sub> level	Disconnect the incoming CO <sub>2</sub> gas	
Low incoming N <sub>2</sub> pressure	Disconnect the incoming N <sub>2</sub> gas	

## 16 Changing the set points

### 16.1 The temperature set point

The temperature setpoint can be adjusted in the range between 24.9 °C to 40.0 °C.

 **The default temperature setpoint is 37.0 °C.**

To change the temperature setpoint, follow these instructions:

1. When the display shows the current temperature:




2. Press the rotatory button and rotate (↺) or (↻) to adjust the setpoint.
3. After changing the temperature value, press the rotatory button to save the setpoint.

If the display does not show the current temperature reading, rotatory button rotation (↺) or (↻) will toggle between the temperature, CO<sub>2</sub>, O<sub>2</sub> and mode readings.

## 16.2 The CO<sub>2</sub> gas concentration set point

The CO<sub>2</sub> concentration can be adjusted in the range between 1.9% to 9.9%.

 **The default CO<sub>2</sub> setpoint is 6.0%.**

To change the CO<sub>2</sub> gas concentration setpoint, follow these instructions:

1. When the display shows the CO<sub>2</sub> gas concentration:




2. Press the rotatory button and rotate (↶) or (↷) to adjust the setpoint.
3. After changing the value, press the rotatory key once more to save the value.

If the display does not show the current temperature reading, rotatory button rotation (↶) or (↷) will toggle between the temperature, CO<sub>2</sub>, O<sub>2</sub> and mode readings.

## 16.3 The O<sub>2</sub> gas concentration set point

The O<sub>2</sub> concentration can be adjusted in the range between 3.9% to 19.9%.

 **The default O<sub>2</sub> setpoint is 5.0%**

To change the O<sub>2</sub> gas concentration setpoint, follow these instructions:

1. When the display shows the O<sub>2</sub> gas concentration:



2. Press the rotatory button and rotate (↶) or (↷) to adjust the setpoint.
3. After changing the value, press the rotatory key once more to save the value.

If the display does not show the current temperature reading, rotatory button rotation (↶) or (↷) will toggle between the temperature, CO<sub>2</sub>, O<sub>2</sub> and mode readings.

## 16.4 The culture mode

The culture mode can be set for under “oil culture” or “open culture”. “Under oil” culture mode is used when the culture media has an oil or Paraffin overlay. “Open culture” mode is used when the culture media does not have any overlay.

 **The default setting is “Oil culture” mode.**

To change the culture mode, follow these instructions:

1. When the display shows the culture mode:



2. Press the rotatory button and rotate (↺) or (↻) to change the mode.
3. When the display shows the desired/correct mode, press the rotatory button. The mode is now set.

If the display does not show the current mode reading, rotatory button rotation (↺) or (↻) will toggle between the temperature, CO<sub>2</sub>, O<sub>2</sub> and mode readings.

Open culture is possible in a 4-well (or similar type of dish) in volumes not under 0.8 mL per well without oil overlay for up to a maximum of 4 hours. The Osmolality will change rapidly after that and reach over 300 mOsm/kg. In a more extended period risk of osmolality changes in media will increase rapidly.

In "Oil culture" mode, the lid temperature is kept 0.2 – 0.3 °C above the temperature setpoint. In "Open culture" mode, the lid temperature will be increased by an additional 1.0 °C above the temperature setpoint (thus making the lid 1.2 – 1.3 °C warmer than the insert surface).

### **Difference between open culture mode and oil culture mode**

The significant difference between open culture mode and oil culture mode is the amount of heat in the lid. Oil accumulates temperature, so higher lid temperature can be accumulated in oil and transferred in media, elevating temperature around the embryo.

Open culture mode is designed not for embryo culturing but (if there is a need) for media equilibration. Do not use open culture mode longer than 4 h. Media volume should be not less than 0.8 mL (in 4 well dishes). If the media stays longer without oil coverage, a high risk of media osmolality changes appears.


If you have any questions or uncertainty, consult Esco Medical or your local representative before using open culture mode in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator.

## 17 Surface temperatures and measuring temperature

In this section, the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators temperature controls system is described in more detail.

The Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators are equipped with 4 completely separate PID controllers for temperature measurement. Each controller is responsible for controlling the temperature of a separate area.

Each of the 4 available areas is equipped with its separate temperature sensor and heater, allowing the user to adjust the temperature in every area separately, thus achieving higher precision.

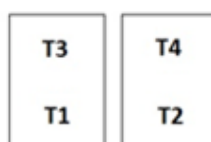

**Figure 17.1** Temperature zones in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators

Each area can be calibrated separately, using the item corresponding to the respective area in the menu. These items are placed in the menu and they are named: T1 CAL, T2 CAL, T3 CAL and T4 CAL.

An overview of the areas associated with the sensor names is shown in the table below:

**Table 17.1** Areas associated with sensors

Area	Bottom	Lid
Compartment 1	T1	T3
Compartment 2	T2	T4

To calibrate the temperature in a particular area, please find the corresponding sensor name and adjust it according to a measurement taken using a high-precision thermometer.



**Temperature calibration is done by adjusting the Tx (where x is the sensor number) according to a measurement done on the spot relevant to the dish placement.**



**After temperature adjustment, give it at least 15 minutes for the temperature to stabilize, use the thermometer to verify the correct temperature on each area.**

Be careful when changing the calibration settings – make sure that only the altered value corresponds to where the measurement is done. Give the system time to adjust.



**There is no crossover heating between the 2 compartments: this is a unique feature of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators. Lid temperature will, however, affect the bottom temperature in a compartment. The delta-T should always be 0.2 °C. Thus, if the bottom temperature 37.0 °C, the lid should be 37.2 °C.**

**Note:** how to calibrate the temperature at the T1 area can be found in this manual's 13.4.1 section.



**"T1" is used to adjust the bottom temperature of compartment 1. "T3" is used to adjust the temperature on the lid in the same compartment. Remember that the delta-T between the top and bottom should always be 0.2 °C.**



**Adjust according to a high precision measurement done with a suitable sensor placed in a dish with media and a mineral oil overlay. Place the dish on one of the designated spots indicated on the heating insert.**



**Proceed to validate if the lid temperature is precisely 0.2 °C higher than the bottom temperature.**



**Stick a suitable and calibrated sensor to the middle of the lid area and close the lid. Wait 15 minutes and record the temperature reading. Adjust the "T3" to the desired level, using the same procedure as described above. It may be necessary to do iterations before the zone is completely calibrated.**

The 2<sup>nd</sup> compartment is adjusted/calibrated in the same manner.


## 18 Pressure


### 18.1 CO<sub>2</sub> gas pressure

The CO<sub>2</sub> pressure can be read out in the CO<sub>2</sub> sub-menu:



The CO<sub>2</sub> pressure is shown in bar. External pressure must be between 0.4 – 0.6 bar (5.80 – 8.70 PSI) at all times. It cannot be adjusted on the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator; it must be done on the external gas regulator.

 **Remember there is a pressure alarm on the pressure limits: if the pressure falls below 0.3 bar or rises above 0.7 bar (4.40– 10.20 PSI).**


 **The internal pressure sensor cannot be calibrated by the user. Under normal circumstances, the pressure sensor is replaced every 2 years according to the maintenance plan.**


### 18.2 N<sub>2</sub> gas pressure

The N<sub>2</sub> pressure can be read out in the O<sub>2</sub> sub-menu:



The N<sub>2</sub> pressure is shown in bar. External pressure must be between 0.4 – 0.6 bar (5.80 – 8.70 PSI) at all times. It cannot be adjusted on the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator; it must be done on the external gas regulator.

 **Remember there is a pressure alarm on the pressure limits if the pressure falls below 0.3 bar or rises above 0.7 bar (4.40 – 10.20 PSI).**

 **The internal pressure sensor cannot be calibrated by the user. Under normal circumstances, the pressure sensor is replaced every 2 years according to the maintenance plan.**

## 19 Firmware

The firmware installed on your Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators is upgradeable. Whenever a critical update is available, it will be provided to our distributors around the world – they will make sure that your Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator runs with the newest available firmware. A service technician can do this during the scheduled annual service.

Please follow these steps to check the firmware which is currently installed on your unit:

1. In the menu locate the Service sub-menu “Serv” and rotate the rotatory button (↻). The display will show currently installed firmware version:



Ver 2.0 is only shown as an example. The current Mini MIRI® Dry multiroom IVF incubator firmware version is 1.4.2, and the Mini MIRI® Humidity multiroom IVF incubator firmware version is 1.4.1.

2. Rotate rotatory button (↻) to exit into the menu.

## 20 Safe sense function

There is a possibility to purchase the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator with an integrated SAFE Sense system.

BCSI created this system to provide pH monitoring within a closed environment (an incubator) to measure pH without disturbing the maintained optimal conditions.

Please read more about SAFE Sense software in the SAFE Sense User manual.

## 21 Cleaning instructions

### 21.1 Considerations about a sterile device

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are not sterile devices. They are not delivered sterile state and it is not possible to keep them sterile when in use.


However, their design was created with great care to make it easy for the user to keep the device sufficiently clean during use and not contaminate the key components.



The design features intended to provide cleanliness include:

- A circulated air system
- A HEPA filter continually cleans the incoming gas
- A HEPA/VOC filter, which continually cleans the air inside the system (do not apply for the Mini MIRI® Humidity multiroom IVF incubator)
- A removable heat optimization plate can be removed and cleaned (cannot be autoclaved!). As this serves as the main holding area for samples, this should be the highest priority to keep clean
- Compartments with sealed edges that can be cleaned
- Use of aluminum and PET parts that withstand cleaning well

## 21.2 Manufacturer recommended cleaning procedure

 **Always validate the cleaning procedures locally; for more guidance, consult either your manufacturer or the distributor.**

The routine cleaning procedure is recommended for regular processing and maintenance. The combination of standard cleaning procedures and disinfection procedures is recommended for event-related concerns such as media spills, visual accumulation of soil and/or other evidence of contamination. It is also recommended to clean and disinfect the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators immediately after any media spills.

### **Periodic cleaning of the device (with no embryos inside)**

Wearing gloves and good handling techniques are essential to successful cleaning.

1. It is recommended to clean the unit with an aqueous 70% isopropyl alcohol. Moisten a sterile wipe and clean all the device's internal and external surfaces by rubbing the surfaces' wipe.
2. After wiping, leave the device's lids open for some time to ensure that all alcohol fumes evaporate.
3. Finally, use purified or sterile water is used to wipe the surfaces of the device.
4. Inspect the device – if visually clean, consider it ready for use.

If the device is not visually clean, repeat the process from step 1.

## 21.3 Manufacturer recommended disinfection procedure

### Disinfection of the device (with no embryos inside)

The use of gloves and good handling techniques are essential for successful disinfection.

Proceed with the following steps (this procedure has been demonstrated during the on-site training program as part of the installation protocol):

1. Power off the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator (rear panel).
2. Open the lids.
3. Use the required disinfectant to disinfect the internal surface and a glass plate on the lid's top. Use sterile wipes to apply the disinfectant.
4. Wipe all internal surfaces and the top of the lid with three wipes at least. Repeat until the wipes are no longer discolored.
5. Change your gloves, and after 10 minutes of contact time, spray sterile water on the surfaces and wipe them with a sterile wipe.
6. Inspect the device – if it is visually clean, consider it ready for use. If the device is visually not clean, go to step 3 and repeat the procedure.
7. Turn on the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator (rear panel).

## 22 Heat optimization plates.

Insert the heat optimization plate. The heat optimization plate will ensure full contact with the dish. It generally means much more stable temperature conditions for the cells. The plate fits the compartment. The plate fits the compartment, and it is taken out for cleaning.



**Do not autoclave at the inserts. It will damage the inserts as high temperature bends them out of shape.**


Place the dish where it fits the pattern. The heat optimization plates can be applicable for Nunc®, Falcon®, Sparmed® and VitroLife® dishes.



**Use only the correct type of heat optimization plates for your dishes.**



Figure 22.1 Heat optimization plate

 **Never incubate without the plates in place and never use non-Esco Medical heating optimization plates. It may cause dangerous and unpredictable temperature conditions that may be harmful to the specimens.**

## 23 Humidification

### 23.1 Mini MIRI® Dry multiroom IVF incubator


The Mini MIRI® Dry multiroom IVF incubator is used for culture with mineral oil or Paraffin overlay. Set the culture mode for "Oil culture". It cannot be used for open culture mode.

The Mini MIRI® Dry multiroom IVF incubator must not be irrigated. Humidification of the Mini MIRI® Dry multiroom IVF incubator will damage the device – condensation will block internal pipes and damage electronic parts.

Open culture is possible in a 4-well (or similar type of dish) in volumes not under 0.8 mL per well without oil overlay for up to a maximum of 4 hours. The Osmolality will change rapidly after that and reach over 300 mOsm/kg. In a more extended period of time risk of osmolality changes in media will increase rapidly.

In "Oil culture" mode, the lid temperature is kept 0.2 – 0.3 °C above the temperature setpoint. In "Open culture" mode, the lid temperature will be increased by an additional 1.0 °C above the temperature setpoint (thus making the lid 1.2 – 1.3 °C warmer than the insert surface).

Please consult Esco Medical or your local representative before using "Open culture" mode in the Mini MIRI® Dry multiroom IVF incubator if you have any questions.

 **"Open culture" mode is designed not for embryo culturing but (if there is a need) for media equilibration.**

## 23.2 Mini MIRI® Humidity multiroom IVF incubator

The system has a built-in humidity sensor. And, the water bottle is closed on the side of the unit of easy control of water level and refilling.

The design runs a simulated humidity routine that will ensure that no evaporation occurs in all standard dishes if they are normally covered with the lid that comes with the dish.

The Mini MIRI® Humidity multiroom IVF incubator maintains humidity level circulating gas in the system through water in the humidification bottle. In other words, the Mini MIRI® Humidity multiroom IVF incubator does not actively control the humidification level in the system to reach certain humidity levels (despite that gas humidification is a continuous process).


👍 **Humidification bottle should be changed each month.**

👍 **Water in the humidification bottle must be changed at least once per week.**

👍 **One-third of the humidification bottle should be filled with sterile water for the Mini MIRI® Humidity to work properly and maintain the required humidity in the system.**

## 24 Temperature validation

The Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators are equipped with 2 PT-1000 Class-B sensors located in the center of the bottom of each compartment.

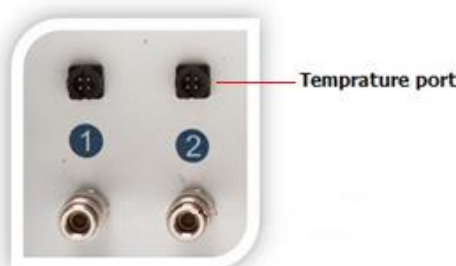

**Figure 24.1** PT-1000 Class B sensors

The sensors serve external validation purposes. They are entirely separate from the circuit of the unit.

The compartment's temperature conditions can be continuously logged through the external connectors on the unit's side without compromising its performance.

Any logging system that uses standard PT-1000 sensors may be used.

Esco Medical can supply an external logging system (MIRI® – GA12) for the sensors.

## 25 Gas level validation

Gas concentration in each compartment of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator can be validated by taking a gas sample from one of the 2 gas sample ports on the unit's side, using a suitable gas analyzer.



**Figure 25.1** Gas sample ports

Each sample port is directly connected to the corresponding compartment with the same number. The gas sample will be taken ONLY from the specific compartment.

 **An external automatic gas sampler can be connected to the ports for continual validation.**

 **Before any gas measurement, make sure that the lids had not been opened for at least 5 minutes.**

 **Taking out a large sample volume may affect gas regulation.**

 **Make sure that the gas analyzer is calibrated before use.**

## 26 Alarm switch for an external system

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator can be connected to an external monitoring system, ensuring maximum safety, especially during nights and weekends. The Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator is equipped with a 3.5 mm jack connector on the back that can be connected to a monitoring device.

Whenever an alarm goes off (that could be temperature alarm, gas alarms for CO<sub>2</sub> or O<sub>2</sub> levels, low-pressure or high-pressure alarms for CO<sub>2</sub> and N<sub>2</sub> gases) or if the power supply to the unit suddenly cut, the switch is indicating that the unit needs to be inspected by the user.

The connector can be connected either to a voltage source OR to a current source.

**⚠ Note that if a current source is attached to the 3.5 mm jack connector, the maximum current rating is between 0 – 1.0 Amp.**

**⚠ If a voltage source is attached, then the limitation is between 0 – 50V AC or DC.**

If there is no alarm, the switch within the unit will be in the "ON" position, as is illustrated below.

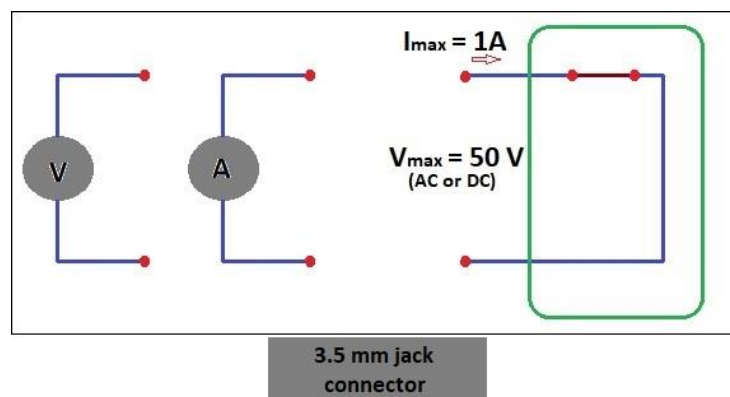


Figure 26.1 No alarm mode

Whenever the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator goes into an alarm mode, the switch will become an "open circuit". It means that no current can run through the system anymore.

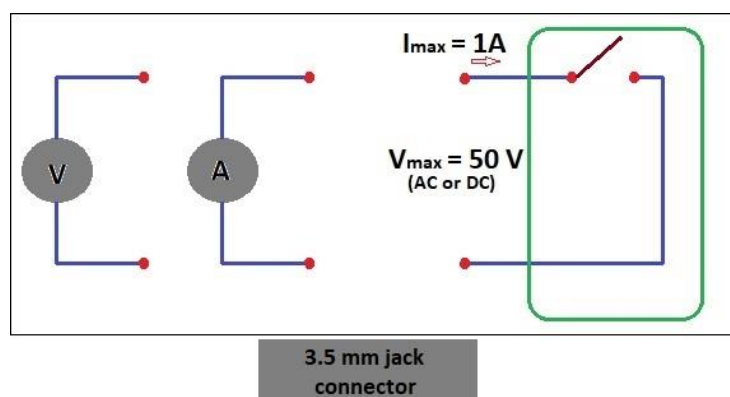



Figure 26.2 "Open circuit" alarm mode

 Whenever the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator's power cord is disconnected from the power source, this switch will automatically indicate an alarm! It is an extra safety feature intended to alert the personnel in case of a power cut in the laboratory.

## 27 Writing area on the compartment lids

Each compartment lid on the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator is made from white glass, optimized for writing text. The compartment's patient data or content can be noted down for easy reference during the incubation process.

The text can be wiped off with a cloth afterward. Use only a suitable non-toxic pen that allows the text to be erased later and will not damage the incubated samples.



**Figure 27.1** Area for patient information

## 28 Maintenance

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are designed to be user-friendly. Reliable and safe operation of this equipment is based on the following conditions:

1. Correct calibration of temperature and gas level, using high-precision equipment in the intervals prescribed based on clinical practice at the laboratory, where the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are used. The manufacturer recommends that the period between validation should be no longer than 14 days.
2. VOC/HEPA filters must be replaced every 3 months.
3. In-line HEPA filters must be replaced yearly during annual maintenance.
4. According to the clinical practice intervals, suitable cleaning is in the laboratory where the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are used. The manufacturer does not recommend periods longer than 14 days between cleaning.



**It is essential to perform the inspection and service at the intervals indicated in the MAINTENANCE section below. Failure to do so can have serious adverse outcomes, causing the unit to stop functioning as expected and cause damage to samples, patients or users.**



**Warranty is considered to be void if service and maintenance are not followed.**



**Warranty is considered void if service and maintenance procedures are done not by trained and authorized personnel.**

## 29 Emergency Procedures

### **Total loss of power to or on the unit:**

- Remove all the samples and place them in an alternative or backup device that is not affected by the problem.
- Without the power source, the internal temperature of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator will drop below 35 °C after being 10 minutes in an ambient environment of 20 °C.
- The CO<sub>2</sub> concentration will remain within 1% of the setpoint for 30 minutes if the lids remain closed.
- If a longer time to turn the power back on is needed, it may be useful to cover the unit with insulating blankets to slow the temperature drop.

### **If a single temperature alarm goes off:**

- Remove the samples from the affected compartment. They can be relocated to any of the other compartments, which happens to be unoccupied. All compartments are separate so that the remaining ones will function normally.

### **If multiple temperature alarms go off:**

- Remove the samples from the affected compartments. They can be relocated to any of the other compartments, which happens to be unoccupied. All compartments are separate so that the remaining ones will function normally.
- Alternatively, remove the samples from all the affected compartments and place them in an alternative or backup device that is not affected by the problem.

### **If the CO<sub>2</sub> level alarm goes off:**

There will be a 30-minute-long interval during which the user can assess if the condition is temporary or permanent. If the state is permanent, remove all the samples and place them in an alternative or backup device that is not affected by the problem. If the



condition is temporary and the CO<sub>2</sub> level is low, keep the lids shut. If the state is temporary and the CO<sub>2</sub> level is high, open a few lids to vent out some CO<sub>2</sub>.

**If the O<sub>2</sub> level alarm goes off:**

Usually, no emergency procedures are necessary in this case. If the condition is judged to be permanent, it may be advantageous to switch off O<sub>2</sub> regulation in the menu.

**If the CO<sub>2</sub> pressure alarm goes off:**

Inspect the external gas supply and gas supply lines. If the problem is external and not readily fixed, follow the guidelines under the section "CO<sub>2</sub> pressure alarm".

**If the O<sub>2</sub> pressure alarm goes off:**

Inspect the external gas supply and gas supply lines. If the problem is external and not readily fixed, follow the guidelines under the "O<sub>2</sub> pressure alarm" section.

**In case of a gas pressure alarm on the unit:**

Inspect the external gas supply and the gas supply lines. If the problem is external and not readily fixed, or if the problem is internal, follow the guidelines under gas level alarm.

## 30 User Troubleshooting

**Table 30.1** Heating system

Symptom	Cause	Action
No heating, the display is off	The unit is switched off at the back or not connected to the power	Switch the device on or connect the power
No heating	The setpoint for temperature is wrong	The temperature is more than 0.5 °C off the set temperature
		Check the desired temperature setpoint
Uneven heating	System not calibrated	Calibrate each zone according to the user manual, using a high precision thermometer

**Table 30.2** CO<sub>2</sub> gas regulator

Symptom	Cause	Action
No CO <sub>2</sub> gas regulation	The system is not powered	Check power mains
	The system is on standby or switched off	Switch the system on
	CO <sub>2</sub> gas regulator is off	Activate CO <sub>2</sub> gas regulator by setting "CO <sub>2</sub> " to "ON" in the menu
	No CO <sub>2</sub> or wrong gas attached to CO <sub>2</sub> gas input	Check the CO <sub>2</sub> gas supply, make sure that pressure is kept stable at 0.4 – 0.6 bar (5.80 – 8.70 PSI)
	The actual gas concentration is higher than the setpoint	Check CO <sub>2</sub> gas setpoint
Poor CO <sub>2</sub> gas regulation	Lid(s) are left open	Close lid(s)
	Seals missing on the lid(s)	Replace the seals on the lid(s)
"A CO <sub>2</sub> " is shown on the display	CO <sub>2</sub> gas concentration more than ±1 from the setpoint	Allow the system to stabilize by closing all lids
"CO <sub>2</sub> P" is shown on the display	No/wrong CO <sub>2</sub> gas pressure to the system	Check the CO <sub>2</sub> gas supply, make sure that pressure is kept stable at 0.4 – 0.6 bar (5.80 – 8.70 PSI)

**Table 30.3** O<sub>2</sub> gas regulator

Symptom	Cause	Action
No O <sub>2</sub> gas regulation	System not powered	Check mains
	The system is on standby or switched off	Switch the system on
	O <sub>2</sub> gas regulator is off	Activate the O <sub>2</sub> gas regulator by setting "O <sub>2</sub> " to "ON" in the menu
	No N <sub>2</sub> or wrong gas type attached to N <sub>2</sub> gas input	Check gas supply, make sure that pressure is kept stable at 0.4 – 0.6 bar (5.80 – 8.70 PSI)
	The actual gas concentration is higher than the setpoint	Check O <sub>2</sub> setpoint
Poor O <sub>2</sub> gas regulation	Lid(s) are left open	Close lid(s)
	Seals missing on the lid(s)	Replace the seals on the lid(s)
"A O <sub>2</sub> " is shown on the display	O <sub>2</sub> gas concentration more than $\pm 1\%$ from the setpoint	Allow the system to stabilize by closing all lids
"N <sub>2</sub> P" is shown on the display	No/wrong N <sub>2</sub> gas pressure to the system	Check N <sub>2</sub> gas supply and ensure that pressure is stable at 0.4 – 0.6 bar (5.80 – 8.70 PSI). If O <sub>2</sub> regulation is not needed, set the "O <sub>2</sub> " to "OFF" in the menu to deactivate O <sub>2</sub> gas regulation and abort the N <sub>2</sub> gas alarm

**Table 30.4** Datalogger

Symptom	Cause	Action
No data is sent to the PC	System not powered	Check mains
	The system is on standby or switched off	Switch the system on
	The data cable between multiroom IVF incubator and PC not correctly attached	Check connection. Use only the cable supplied with the unit
	Data logger software/USB driver not correctly installed	Please refer to the software installation guide

**Table 30.5** Display

Symptom	Cause	Action
Missing segment(s) in display	Failure in the PCB	Contact your Escó Medical Distributor to replace the PCB

**Table 30.6** Keyboard

Symptom	Cause	Action
The absent or erratic function of keys	Failure in the keys	Contact your Escó Medical Distributor to replace the keys

## 31 Specifications

**Table 31.1** Mini MIRI® Dry multiroom IVF incubator specifications

Technical specifications	Mini MIRI® Dry
Overall dimensions (W×D×H)	525 x 420 x 230 mm
Weight	22 kg
Material	Mild steel / Aluminum / PET / Stainless steel
Power supply	115V 60Hz OR 230V 50Hz
Power consumption	160 W
Temperature control range	24.9 °C – 40.0 °C
Gas consumptions (CO <sub>2</sub> ) <sup>3</sup>	< 2 liters per hour
Gas consumption (N <sub>2</sub> ) <sup>4</sup>	< 8 liters per hour
Premixed gas consumption	In purge < 50 liters per hour In normal run < 1 liter per hour
CO <sub>2</sub> range	1.9 % – 9.9%
O <sub>2</sub> range	3.9 % – 19.9%
Gas pressure CO <sub>2</sub> (input)	0.4 – 0.6 bar (5.80 – 8.70 PSI)
Gas pressure N <sub>2</sub> (input)	0.4 – 0.6 bar (5.80 – 8.70 PSI)
Alarms	Audible and visible for out-of-range temperature, gas concentration and gas pressure.
Shelf life	1 year
Use life	6 years

**Table 31.1** Mini MIRI® Humidity multiroom IVF incubator specifications

Technical specifications	Mini MIRI® Humidity
Overall dimensions (W×D×H)	525 x 420 x 230 mm
Weight	22 kg
Material	Mild steel / Aluminum / PET / Stainless steel
Power supply	115V 60Hz OR 230V 50Hz
Power consumption	160 W
Temperature control range	24.9 °C – 40.0 °C
Gas consumptions (CO <sub>2</sub> ) <sup>3</sup>	< 4 liters per hour
Gas consumption (N <sub>2</sub> ) <sup>4</sup>	< 12 liters per hour
Premixed gas consumption	In purge < 50 liters per hour In normal run < 1 liter per hour
CO <sub>2</sub> range	1.9 % – 9.9%
O <sub>2</sub> range	3.9 % – 19.9%
Gas pressure CO <sub>2</sub> (input)	0.4 – 0.6 bar (5.80 – 8.70 PSI)
Gas pressure N <sub>2</sub> (input)	0.4 – 0.6 bar (5.80 – 8.70 PSI)
Alarms	Audible and visible for out-of-range temperature, gas concentration and gas pressure.
Shelf life	1 year
Use life	6 years

<sup>3</sup> Under normal conditions (CO<sub>2</sub> set point reached at 6.0%, all lids closed)

<sup>4</sup> Under normal conditions (O<sub>2</sub> set point reached at 5.0%, all lids closed)

## 32 Electromagnetic compatibility

**Table 32.1** Electromagnetic emissions

Guidance and manufacturer's declaration – electromagnetic emissions		
The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are intended for use in the electromagnetic environment specified below. The customer or the user of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator should ensure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators do not use RF energy. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class A	The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are suitable for use in a hospital environment.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Class A	It is not suited for domestic establishments.

**Table 32.2** Electromagnetic immunity

Guidance and manufacturer's declaration – electromagnetic immunity			
The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are intended for use in the electromagnetic environment specified below. The customer or the user of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator should ensure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment- guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact  ±8 kV air	±6 kV contact  ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV for input/ output lines		
Surge IEC 61000-4-5	±1kV differential mode ±2kV common mode		
Voltage dips, short interruptions and voltage variations on power supply input lines  IEC 61000-4-11	<5 % 100V (>95% dip in 100V) for 0.5 cycle 40% 100V (60% dip in 100V) for 5 cycles 70% 100V (30% dip in 100V) for 25 cycles) dip in 100V) for 5 sec		

Power frequency (50/60 Hz) magnetic field  IEC 61000-4-8	3 A/m	Performance A	The power-frequency magnetic fields' level should be characteristic of a specific location in a commercial or hospital environment.
<b>Guidance and manufacturer's declaration – electromagnetic immunity</b>			
The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are intended for use in the electromagnetic environment specified below. The customer or the user of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator should ensure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment- guidance
Conducted RF IEC 61000-4-6  Radiated RF IEC 61000-4-3	10 Vrms 150kHz to 80 MHz in ISM bands  3 V/m 80 MHz to 2.5 GHz	3V/m from 80MHz to 2.5 GHz	<p>Portable and mobile RF communications equipment should be used no closer to any part of Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators, including cables, than the recommended separation distance calculated according to the equation, applicable to the transmitter's frequency.</p> <p><b>Recommended separation distance</b></p> <p><math>d = 0.35 P</math></p> <p><math>d = 0.35 P</math>, 80MHz to 800MHz  <math>d = 0.7 P</math>, 800MHz to 2.5GHz</p> <p>P is the maximum power output rating of the transmitter in watts (W) according to the transmitter manufacturer, <math>d</math> is the recommended separation distance in meters (m).</p> <p>As determined by an electromagnetic site survey, field strengths from fixed RF transmitters should be lower than the compliance level in each frequency range.</p> <p>Interference may occur in the vicinity of the equipment.</p>

**Table 32.3** Recommended separation distances

Recommended separation distances between portable and mobile RF communication equipment and Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators			
The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are intended to be used in an electromagnetic environment in which radiated RF disturbances are controlled. The customer, or the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator user, can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters). The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators are recommended below, according to the communications equipment's maximum output power.			
The rated maximum output power of the transmitter	Separation distance according to the frequency of the transmitter (m)		
	150 kHz to 80 MHz $d = 1.2\sqrt{P}$	80 MHz to 800 MHz $d = 1.2\sqrt{P}$	800 MHz to 2.5GHz $d = 2.3\sqrt{P}$
0.01 W	0.1m	0.1m	0.2m
0.1 W	0.4m	0.4m	0.7m
1 W	1.2m	1.2m	2.3m
10 W	3.7m	3.7m	7.4m
100 W	11.7m	11.7m	23.3m
For transmitters rated at a maximum output power not listed above, the recommended separation distance $d$ in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where $P$ is the maximum output power rating of the transmitter in watts (w), according to the transmitter manufacturer.			
<b>NOTE 1:</b> at 80 MHz and 800 MHz, the higher frequency range's separation distance applies.			
<b>NOTE 2:</b> these guidelines may not apply in all situations.			
Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Medical Devices may be affected by cellular telephones and other personal or household devices not intended for medical facilities. It is recommended to ensure that all equipment used near the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators product complies with the medical electromagnetic compatibility standard and checks before use that no interference is evident or possible. If the interference is suspected or probable, switching off the offending device is the specific solution as it is the usual practice in aircraft and medical facilities.

Medical electrical equipment must be treated with special precautions indicated by EMC and must be installed and put into service according to the EMC information provided. Portable and mobile RF communications equipment can affect medical electrical equipment.

## 33 The Validation guide

### 33.1 Product release criteria

The Esco Medical Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators undergoes strict quality and performance testing before being released for sale.

#### 33.1.1 Performance

Each component used in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator is tested during the manufacturing process to ensure a defect-free unit.

Before release, the incubator is tested per a release test having a duration of at least 24 hours, using high-performance thermometers and gas analyzers, along with real-time data logging to ensure that the unit lives up to expected performance standards.

**Pass I:** Internal sensor temperature variation from setpoint within  $\pm 0.1$  °C absolute.

**Pass II:** Internal sensor CO<sub>2</sub> concentration variation from setpoint within  $\pm 0.2\%$  absolute.

**Pass III:** Internal sensor N<sub>2</sub> concentration variation from setpoint within  $\pm 0.2\%$  absolute.

**Pass IV:** Gas flow CO<sub>2</sub> less than 2 l/h (for Mini MIRI® Humidity model – less than 4 l/h).

**Pass V:** Gas flow N<sub>2</sub> less than 8 l/h (for Mini MIRI® Humidity model – less than 12 l/h).

#### 33.1.2 Electrical safety

An electric safety test is also carried out using a high-performance medical safety tester with each unit to ensure that electric requirements for medical devices defined by the EN60601-1 3<sup>rd</sup> edition standards are met.

#### 33.1.3 Communication & data logging

Each unit is connected to a computer running the Mini MIRI® multiroom IVF incubator data logging software. Gas is supplied to the unit, and the system is activated. The data received by the PC program is analyzed to ensure communication between the multiroom IVF incubator and the PC.

#### 33.1.4 Gas concentration levels and consumption

A leak test is performed on each compartment. The maximum leakage allowed through the seals is 0.0 l/h.

The average CO<sub>2</sub> gas variation must stay within setpoint  $\pm 0.2\%$  absolute on all external sampling and internal sensor readings.



The gas flow under regular operation in the Mini MIRI® Dry multiroom IVF incubator is less than 2 liters per hour, whereas in the Mini MIRI® Humidity multiroom IVF incubator – 4 liters per hour. The average should be below 2 liters in the Mini MIRI® Dry multiroom IVF incubator, whereas in the Mini MIRI® Humidity multiroom IVF incubator – below 4 liters.

The average N<sub>2</sub> gas variation must stay within setpoint  $\pm 0.2\%$  absolute on all external and internal sensor readings.

The gas flow under regular operation in the Mini MIRI® Dry multiroom IVF incubator is less than 2 liters per hour, whereas in the Mini MIRI® Humidity multiroom IVF incubator – 8 liters per hour. The average should be below 2 liters in the Mini MIRI® Dry multiroom IVF incubator, whereas in the Mini MIRI® Humidity multiroom IVF incubator – below 12 liters.

### 33.1.5 Cosmetic

- No misalignment in the lids.
- Each lid should be opened and closed easily.
- The seals for the lids must be appropriately attached and aligned.
- No scratches or missing paint on the cabinet.
- Overall, the device must be presentable as a high-quality item.
- The heat optimization plates are checked for misalignment and shape. These are placed into the compartments to check for any mismatch due to the compartment and aluminum blocks' sizes.

## 34 Validation on-site

Even though at Esco Medical, we strive to do the most comprehensive tests before the device is shipped to the customer, there is no way to be sure that everything is still OK at the location when the device is set up.

Therefore, in keeping with established good medical device practice, we have set up a validation test regimen that must be completed before the device can be accepted into clinical use.

In the following, we describe these tests and the equipment necessary to perform them.

A test documentation form is also provided. A copy must be provided to Esco Medical for internal device tracking and device history record.

## 34.1 Mandatory equipment

 **All equipment must be of high quality and calibrated.**

- A thermometer with a suitable sensor for measuring in a droplet of media covered with Paraffin oil with a resolution minimum of 0.1 °C
- A thermometer with a suitable sensor for measuring on an aluminum surface with a resolution minimum of 0.1 °C
- A CO<sub>2</sub> analyzer with a range of 0.0 – 10.0%.
- An O<sub>2</sub> analyzer with a range of 0.0 – 20.0%.
- A Pressure tester with a range of 0.0 – 1.0 bar.
- A Multimeter.

## 34.2 Recommended additional equipment

 **All equipment must be of high quality and calibrated.**

- A VOC meter able to measure the most common volatile organic compounds at least at the ppm-level.
- A sample should be taken just above Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator with the laser particle counter. The reading should be noted down as the background particle level.


Recommended additional equipment can be used for further installation testing that will minimize the likelihood of problems on-site.

## 35 Testing

### 35.1 Gas supply CO<sub>2</sub>

For the regulation system to maintain the correct CO<sub>2</sub> concentration level in the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator compartments, the device must be connected to a stable source of 100% CO<sub>2</sub> 0.4 – 0.6 bar (5.80 – 8.70 PSI) of pressure.

Measure the CO<sub>2</sub> concentration in the gas supply by routing the gas line into a bottle without a lid and a suitably large opening. Set the pressure/flow so that the bottle is flushed continually with gas, without increasing pressure in the bottle (i.e., the amount of gas exiting the bottle should be equal to the gas volume entering the bottle).

 **Pressure build-up will affect the measured CO<sub>2</sub> concentration, as CO<sub>2</sub> concentration is pressure-dependent.**

Sample from the bottle near the bottom with the gas analyzer.

**PASS: CO<sub>2</sub> concentration measured must be between 98.0% - 100%.**



**Use of CO<sub>2</sub> gas with moisture will damage the flow sensors. Moisture level must be verified on the gas manufacturer's certificate: only 0.0 ppm v/v Max is permissible.**

### 35.1.1 About CO<sub>2</sub>

Carbon dioxide (CO<sub>2</sub>) is a colorless, odorless, non-combustible gas. Carbon dioxide above the triple point temperature of -56.6 °C and below the critical point temperature of 31.1 °C can exist in both a gaseous and a liquid state.

Bulk liquid carbon dioxide is commonly maintained as a refrigerated liquid and vapor at pressures between 1,230 kPa (approx. 12 bar) and 2,557 kPa (approx. 25 bar). Carbon dioxide may also exist as a white opaque solid with a temperature of -78.5 °C under atmospheric pressure.



**A high concentration of carbon dioxide (10.0% or more) can asphyxiate quickly without warning with no possibility of self-rescue regardless of the oxygen concentration.**

The User should make sure the CO<sub>2</sub> used is safe and moisture-free. Below is a list of some common component concentrations. Please note that the values given are NOT the proper amounts, only an example:

- Assay 99.9% v/v min.
- Moisture 50 ppm v/v max. (20 ppm w/w max).
- Ammonia 2.5 ppm v/v max.
- Oxygen 30 ppm v/v max.
- Oxides of Nitrogen (NO/NO<sub>2</sub>) 2.5 ppm v/v max each.
- The non-volatile residue (particulates) 10 ppm w/w max.
- The non-volatile organic residue (oil and grease) 5 ppm w/w max.
- Phosphine 0.3 ppm v/v max.
- Total volatile hydrocarbons (calculated as methane) 50 ppm v/v max. of which 20 ppm v/v.
- Acetaldehyde 0.2 ppm v/v max.
- Benzene 0.02 ppm v/v max.
- Carbon Monoxide 10 ppm v/v max.

- Methanol 10 ppm v/v max.
- Hydrogen Cyanide 0.5 ppm v/v max.
- Total Sulphur (as S) 0.1 ppm v/v max.

## 35.2 Gas supply N<sub>2</sub>

For the regulation to maintain the correct O<sub>2</sub> concentration levels in the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator compartments, the device must be connected to a stable source of 100% N<sub>2</sub> at 0.4 – 0.6 bar (5.80 – 8.70 PSI) of pressure.

Measure the N<sub>2</sub> concentration in the gas supply by routing the gas line into a bottle without a lid and a suitably large opening. Set the pressure/flow so that the bottle is flushed continually with gas, without increasing pressure in the bottle (i.e., the amount of gas exiting the bottle should be equal to the gas volume entering the bottle).

Sample from the bottle near the bottom with the gas analyzer.



**A gas analyzer that can measure 0% O<sub>2</sub> accurately can be used. 100% N<sub>2</sub> = 0 % O<sub>2</sub>.**

**PASS: N<sub>2</sub> concentration measured must be between 95.0%– 100%.**



**The use of N<sub>2</sub> gas with moisture will damage the flow sensors. Moisture level must be verified on the gas manufacturer's certificate: only 0.0 ppm v/v Max is permissible.**

### 35.2.1 About N<sub>2</sub>

Nitrogen makes up a significant portion of the earth's atmosphere with 78.08% by volume. Nitrogen is a colorless, odorless, tasteless, non-toxic, and almost inert gas. Nitrogen is principally shipped and used in either gaseous or liquid form.



**N<sub>2</sub> gas can act as a simple asphyxiant by displacing air.**

The User should make sure the N<sub>2</sub> used is safe and moisture-free. Below is a list of some standard component concentrations. Please note that the values given are NOT the proper amounts, only an example:

- Research Grade 99.9995%.
- Contaminant.
- Argon (Ar) 5.0 ppm.

- Carbon Dioxide (CO<sub>2</sub>) 1.0 ppm.
- Carbon Monoxide (CO) 1.0 ppm.
- Hydrogen (H<sub>2</sub>) 0.5 ppm.
- Methane 0.5 ppm.
- Oxygen (O<sub>2</sub>) 0.5 ppm.
- Water (H<sub>2</sub>O) 0.5 ppm.

### 35.3 CO<sub>2</sub> gas pressure check

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator requires a pressure of 0.4 – 0.6 bar (5.80 – 8.70 PSI) on the input CO<sub>2</sub> gas line. This gas pressure must be held stable at all times.

For safety, this unit has a built-in digital gas pressure sensor that monitors the incoming gas pressure and alerts the User if any drop is detected.

Remove the inlet gas line for the CO<sub>2</sub> gas. Attach the gas line to the gas pressure measuring device.

**PASS: The value must be 0.4 – 0.6 bar.**

Please refer to the “18.1 CO<sub>2</sub> gas pressure” section for more information.

### 35.4 N<sub>2</sub> gas pressure check

The Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator requires a pressure of 0.4 – 0.6 bar (5.80 – 8.70 PSI) on the input N<sub>2</sub> gas line. This gas pressure must be held stable at all times.

For safety, this unit has a built-in digital gas pressure sensor that monitors the incoming gas pressure and alerts the User if any drop is detected.

Remove the inlet gas line for the N<sub>2</sub> gas. Attach the gas line to the gas pressure measuring device.

**PASS: The value must be 0.4 – 0.6 bar.**

Please refer to the “18.2 N<sub>2</sub> gas pressure” section for more information.

### 35.5 Voltage supply

The voltage on-site must be verified.


Measure the UPS's output plug that the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator will be connected. Also, check that the UPS is attached to a properly grounded mains outlet.

Use a multimeter set for AC.

**PASS: 230V  $\pm$  10.0%**  
**115V  $\pm$  10.0%**

### 35.6 CO<sub>2</sub> gas concentration check

The CO<sub>2</sub> gas concentration is checked for deviation. The gas sample port on the back of the unit is used. Use sample port-2 for validation.

 **Remember not to open any lid at least 15 min before starting the test nor during the testing itself.**

Hook-up the gas analyzer inlet tube to the sample port. Make sure that the fit is perfect and that no air can enter or exit the system.


The gas analyzer must have a gas return port connected to the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator (i.e., another compartment). Only measure while the value on the gas analyzer stabilizes.

Please refer to the "13.4.2 CO<sub>2</sub> sub-menu" section for more information on how to perform the CO<sub>2</sub> gas calibration.

**PASS: CO<sub>2</sub> concentration measured must not deviate more than  $\pm$  0.2% from the setpoint.**

### 35.7 O<sub>2</sub> gas concentration check

The O<sub>2</sub> gas concentration is checked for deviation. The gas sample port on the back of the unit is used. Use sample port-2 for validation.

 **Remember not to open any lid at least 10 min before starting the test nor during the testing itself.**

Hook-up the gas analyzer inlet tube to the sample port. Make sure that the fit is perfect and that no air can enter or exit the system.

The gas analyzer must have a gas return port connected to the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator (i.e., another compartment). Only measure while the value on the gas analyzer stabilizes.

Please refer to the "13.4.3 O<sub>2</sub> sub-menu" section for more information on how to perform the CO<sub>2</sub> gas calibration.

**PASS: O<sub>2</sub> concentration measured must not deviate more than  $\pm 0.2\%$  from the setpoint.**

### 35.8 Temperature check: Compartment bottoms

The first part of the temperature check is performed using a thermometer with a sensor suitable for measuring temperature in a droplet of media covered with Paraffin oil, with a resolution of 0.1 °C as a minimum.

At least 6 dishes prepared in advance (with at least one microdroplet of media approximately 10 – 100 µL in each dish). The media should be covered with a layer of Paraffin oil. The dishes need not be equilibrated, as the pH will not be measured during the validation tests.

The dishes are placed with at least one dish in each compartment. The dishes should be placed on the corresponding size slot on the heat optimization plates.

Let the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator to warm the dishes and stabilize for at least 1-hour.

Open a compartment lid, remove the cover from the dish and place the sensor tip inside the droplet.

If the measuring device has a fast response time (less than 10 seconds), the quick droplet measurement method should give a useful result.

If the measuring device is slower, a method for retaining the sensor in the droplet spot must be found. Usually, taping the sensor led to a spot inside the compartment bottom is possible. Then close the lid and wait until the temperature has stabilized. Be careful when closing the lid so as not to dislocate the sensor placement in the droplet.

Place the thermometer sensor on each zone and verify the temperature.


If calibration is needed, please refer to the "13.4.1 Temperature sub-menu" section for more information on how to perform the temperature calibration.

**PASS: all temperatures measured on the bottom of the compartments where the dishes are located must not deviate more than  $\pm 0.1$  °C from the setpoint.**

### 35.9 Temperature check: Compartment lids

The second part of the temperature validation is performed using a thermometer with a suitable sensor for measuring temperature on an aluminum surface, with a resolution of 0.1 °C as a minimum.


Tape the sensor to the center of the lid and carefully close the lid. Ensure that the tape keeps the sensor in complete contact with the surface area of the aluminum.

 **Taping the inside of the lid is not an optimal procedure, as the tape will act as an insulator from the heat generated by the bottom heater. However, it is a usable compromise if the taped area's size is kept small and the tape used is strong, thin and light.**

Place the thermometer on each zone and verify the temperature.

**Pass: all temperatures measured on the compartments' lid must not deviate more than  $\pm 0.2$  °C from the setpoint.**

If calibration is needed, please refer to the "13.4.1 Temperature sub-menu" section for more information on how to perform the temperature calibration.

 **An iterative process may be needed if differences in the temperature levels are found and compensated through the calibration procedures. Bottom and lid temperatures will affect each other to some extent. There will be no crossover heat noticeable between compartments.**

### 35.10 6-hour stability test

Following the careful validation of the single parameter, a 6-hour (minimum duration) check must be initiated.

The device must be set up as closely as to the condition under which it will be running in clinical use.

If the preference of CO<sub>2</sub> setpoint is 6.0% or temperature is different from the default setting, an adjustment needs to be done before the test.



If the device will not be clinically operational with the O<sub>2</sub> regulation activated, but there is N<sub>2</sub> gas available, the test should be conducted with O<sub>2</sub> regulation switched on and with N<sub>2</sub> gas supply.

If the N<sub>2</sub> is not available, the test can be done without it.

Make sure that the Esco Medical data logger software is running.

Check that parameter are logged and give a meaningful reading. Let the device run without interfering for at least 6 hours. Analyze the results on the graphs.

**Pass I:** Internal sensor temperature variation from set point is within  $\pm 0.1$  °C absolute.


**Pass II** Internal sensor CO<sub>2</sub> concentration variation from setpoint within  $\pm 0.2\%$  absolute.

**Pass III:** Internal sensor N<sub>2</sub> concentration variation from setpoint within  $\pm 0.2\%$  absolute.

**Pass IV:** Gas flow CO<sub>2</sub> less than 2 l/h (for Mini MIRI® Humidity model – less than 4 l/h).

**Pass V:** Gas flow N<sub>2</sub> less than 8 l/h (for Mini MIRI® Humidity model – less than 12 l/h).

## 35.11 Cleaning

 **Always validate the cleaning procedures locally or consult the manufacturer or the distributor for more guidance.**


After the testing has been conducted successfully, it should be cleaned again before the device is taken into clinical use.

Inspect the unit for physical signs of dirt or dust. The unit should look generally tidy.

- Clean the unit externally with a lint-free cloth moistened with a 70% alcohol solution.
- Switch off the device and remove the mains lead.
- Remove all heat optimization plates and clean them with a lint-free cloth soaked with a 70% alcohol solution.
- Wipe the inside of the 2 compartments with a lint-free cloth moistened with a 70% alcohol solution.
- Wipe the lids in the same way.
- Let the lids remain open for 5 minutes.
- Wipe the 2 compartments and the heat optimization plates with a lint-free cloth soaked in sterilized water.
- Wipe the lids in the same way.
- Let the lids remain open for 10 minutes.
- Re-insert the heat optimization plates.
- Close the lids.

- Attach the power again and switch-on the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator.
- Let the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator run empty for at least 20 minutes before inserting a sample.

### 35.12 Test documentation form

 The "Installation report" form must be completed with the tests-passed status filled by installation personnel and submitted to Esco Medical before the device is taken into clinical use.

### 35.13 Recommended additional testing

#### 35.13.1 A VOC meter (applicable only for the Mini MIRI® Dry model)

A sample should be taken just above the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator with the VOC meter. The reading should be noted down as the background VOC level. Then a sample is taken from the gas sample port number – 2.

**Pass: 0.0 ppm VOC.**

 Ensure that the sample lines do not contain any VOC.

#### 35.13.2 A laser particle counter

A sample should be taken just above the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator with the laser particle counter. The reading should be noted down as the background particle level. Then a sample is taken from the gas sample port number – 2.

**Pass: 0.3-micron < 100 ppm.**

 Ensure that the sample lines do not contain any particles.

## 36 Clinical use

Congratulations! Your device is now ready for clinical use with the validation tests completed and the test report submitted to Esco Medical.

It should provide many years of stable service.

It is necessary to monitor the performance of the device continually.

Use the below scheme for in-use validation.



**Do not attempt to run the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator for clinical purposes without access to high-grade quality control validation equipment.**

**Table 36.1** Validation intervals

Task	Every day	Every week
Temperature check		X
CO <sub>2</sub> gas concentration check	X	
O <sub>2</sub> gas concentration check	X	
Check log for anomalies		X
CO <sub>2</sub> gas pressure check	X	
N <sub>2</sub> gas pressure check	X	

### 36.1 Temperature check

The temperature check is performed using a high-precision thermometer. Place the thermometer on each zone and verify the temperature. Calibrate if necessary.

Please refer to the "13.4.1 Temperature sub-menu" section for more information on how to perform the temperature calibration.

**PASS:**

- All temperatures measured on the bottom of the compartment in the locations where the dishes would be placed must not deviate more than  $\pm 0.2$  °C from the setpoint.
- All temperatures measured on the lid must not deviate more than  $\pm 0.5$  °C from the setpoint.

### 36.2 CO<sub>2</sub> gas concentration check

The CO<sub>2</sub> gas concentration is checked for deviations. The gas sample port on the side of the unit is used for this. Use sample port-2 for validation. It is essential to have a high-precision gas analyzer for CO<sub>2</sub> and O<sub>2</sub> available to do the test.

Please follow these simple rules while testing gas concentration:

- Check the CO<sub>2</sub> gas setpoint.
- Check the actual CO<sub>2</sub> gas concentration to ensure the setpoint is reached, and gas concentration is stabilized around the setpoint.
- Remember not to open any lids for at least 10 min, before starting the test or during the testing itself.

Please refer to the "13.4.2 CO<sub>2</sub> sub-menu" section for more information on how to perform the CO<sub>2</sub> gas calibration.

**PASS: CO<sub>2</sub> concentration measured must not deviate more than  $\pm 0.2\%$  from the setpoint.**

### 36.3 O<sub>2</sub> gas concentration check


The O<sub>2</sub> gas concentration is checked for deviations. The gas sample port on the side of the unit is used for this. Use sample port-2 for validation. It is essential to have a high-precision gas analyzer for CO<sub>2</sub> and O<sub>2</sub> available to do the test.

Please follow these simple rules while testing gas concentration:

- Check the O<sub>2</sub> gas setpoint.
- Check the actual O<sub>2</sub> gas concentration to ensure the setpoint is reached, and gas concentration is stabilized around the setpoint.
- Remember not to open any lids for at least 10 min, before starting the test or during the testing itself.

Please refer to the "13.4.3 O<sub>2</sub> sub-menu" section for more information on how to perform the O<sub>2</sub> gas calibration.

**PASS: O<sub>2</sub> concentration measured must not deviate more than  $\pm 0.2\%$  from the setpoint.**

 Gas analyzers use a small pump to draw outgas from the location being sampled. The pump capacity varies from brand to brand. The gas analyzer's ability to return the gas sample to the Mini MIRI® Dry and Mini MIRI® Humidity incubator (loop sampling) avoids negative pressure and ensures accuracy. Performance of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator will not be affected. The gas in the compartment is not under pressure, and the reading is just an artifact based on unsuitable measuring equipment. Contact Esco Medical or the local distributor for further guidance.

### 36.4 CO<sub>2</sub> gas pressure check

The Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator requires a pressure of 0.4 – 0.6 bar on the input CO<sub>2</sub> gas line. This gas pressure must be held stable at any time.

For safety reasons, this unit has a built-in digital gas pressure sensor control that monitors the incoming gas pressure and alerts the user if any drop is detected.

It is recommended to check the CO<sub>2</sub> gas pressure in the menu by inspecting the value for an item called 'CO<sub>2</sub> P' (CO<sub>2</sub> pressure).

**PASS:** The value must be 0.4 – 0.6 bar.

Please refer to the "18.1 CO<sub>2</sub> gas pressure" section for more information.

### 36.5 N<sub>2</sub> gas pressure check

The Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator requires a pressure of 0.4 – 0.6 bar on the input N<sub>2</sub> gas line. This gas pressure must be held stable at any time.

For safety reasons, this unit has a built-in digital gas pressure sensor control that monitors the incoming gas pressure and alerts the user if any drop is detected.

It is recommended to check the N<sub>2</sub> gas pressure in the menu by inspecting the value for an item called 'N<sub>2</sub> P' (N<sub>2</sub> pressure).

**PASS:** The value must be 0.4 – 0.6 bar.

Please refer to the "18.2 N<sub>2</sub> gas pressure" section for more information.

## 37 The Maintenance guide

Your Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator from Esco Medical contains high precision quality components. These components are chosen to ensure the high durability and performance of the equipment.

However, continual validation of the performance is necessary.

User validation should be done as a minimum according to instructions given in the "33 The Validation guide" section.

If problems are encountered, contact Esco Medical or your local representative.

However, to sustain the high-performance level and avoid system errors, the owner is responsible for having a certified technician who performs components replacements according to table 37.1.

These components must be replaced in the time intervals specified below. Failure to follow these instructions may, in the worst-case scenario, result in damage to the specimens in the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator.



**Warranty void if service intervals are not followed according to table 37.1.**



**Warranty void if non-original parts are used or non-trained and non-authorized personnel carry out the servicing.**

The table below shows time intervals in which components must be replaced:

**Table 37.1** Service interval plan

Component name	Every month	Every 3 month	Every year	Every 2 years	Every 3 years	Every 4 years
VOC/HEPA filter capsule <sup>5</sup>		X				
Humidification Bottle <sup>6</sup>	X					
HEPA in-line filter for CO <sub>2</sub> gas			X			
HEPA in-line filter for N <sub>2</sub> gas			X			
O <sub>2</sub> sensor			X			
CO <sub>2</sub> sensor						X
UV light <sup>5</sup>			X			
Cooling fan					X	
Internal gas pump <sup>5</sup>				X		
Pump module <sup>6</sup>				X		
Proportional valves					X	
Gas lines					X	
Flow sensors				X		
Pressure regulators						X
Internal 0.2μ filter for CO <sub>2</sub>			X			
Internal 0.2μ filter for N <sub>2</sub>			X			
A firmware update (if a new version has been released)			X			

### 37.1 VOC/HEPA filter capsule (applicable only for the Mini MIRI® model)

The VOC/HEPA filter capsule is placed on the Mini MIRI® Dry multiroom IVF incubator unit's back for easy replacement. In addition to the active carbon component, this capsule also has an integrated HEPA filter inside, enabling it to remove particles and volatile organic compounds from the air being re-circulated to the compartments. Because of the carbon component's lifespan, all VOC filters' lifetime is limited, and they must be replaced

<sup>5</sup> Only for the Mini MIRI® Dry model

<sup>6</sup> Only for the Mini MIRI® Humidity model

often. According to table 37.1, the VOC filter installed in the Mini MIRI® Dry multiroom IVF incubator must be replaced every 3 months.

Please follow these safety precautions when changing the VOC filter:

- Always use the original filter (contact Esco Medical or your local distributor for more details or ordering).
- Change filter every 3 months.
- Failure to change the filter on time will result in low/no air-cleaning within the system.
- Warranty void if wrong/non-original filter is used.

Please refer to the “12.1 Installation of a new filter capsule” section for the replacement instructions.

### 37.2 Humidification bottle (applicable only for the Mini MIRI® Humidity model)

A humidification bottle contains water that is used to maintain the humidity in the chamber. It should be changed each month.

Water in the humidification bottle must be changed at least once per week.

### 37.3 In-line HEPA filter for CO<sub>2</sub> gas

The round-shape in-line 0.2μ HEPA filter for CO<sub>2</sub> gas removes any particles found in the incoming CO<sub>2</sub> gas. Failure to use the HEPA filter may cause damage to the high precision flow sensor, calculate the amount of CO<sub>2</sub> gas entering the system, and disturb the CO<sub>2</sub> regulator system.

Please follow these safety precautions when changing the filter:

- Always use the original filter (contact Esco Medical or your local distributor for more details or ordering).
- Change the filter once every year.
- Failure to change the filter on time will result in low/no cleaning of incoming CO<sub>2</sub> gas.
- Warranty void if wrong/ non-original filter is used.

Please refer to the service manual for replacement instructions.

### 37.4 In-line HEPA filter for N<sub>2</sub> gas

The round-shape in-line 0.2μ HEPA filter for N<sub>2</sub> gas removes any particles found in the incoming N<sub>2</sub> gas. Failure to use the HEPA filter may cause damage to the high precision flow sensor, calculate the amount of N<sub>2</sub> gas entering the system, and disturb the N<sub>2</sub> regulator system.

Please follow these safety precautions when changing the filter:


- Always use the original filter (contact Esco Medical or your local distributor for more details or ordering).
- Change the filter once every year.
- Failure to change the filter on time will result in low/no cleaning of incoming N<sub>2</sub> gas.
- Warranty void if wrong/ non-original filter is used.

Please refer to the service manual for replacement instructions.

### 37.5 O<sub>2</sub> sensor

The oxygen regulation uses the Oxygen sensor to keep the O<sub>2</sub> gas concentration at a desired level inside the chambers. The lifetime of this sensor is limited due to its construction. From the day the sensor is unpacked, a chemical process is activated within the sensor core. The chemical reaction is entirely harmless to its surroundings. Still, it is necessary for measuring the amount of oxygen with very high precision needed in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators.

After 1-year, the chemical process in the sensor core stops and the sensor must be replaced. Therefore, it is essential to replace this sensor **WITHIN year from the date it was unpacked and installed.**

 **Oxygen sensors must be replaced at least once every year from the date they were installed in the unit. Irrespective of the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator being used or not.**

In the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator "Installation report" form, the User will see when this sensor was installed. This date must be used to calculate the date for the following O<sub>2</sub> sensor change.



Please follow these safety precautions when changing sensor:

- Always use an original O<sub>2</sub> sensor (contact Esco Medical or your local distributor for more details or ordering).
- Change the O<sub>2</sub> sensor within 1 year from the date of the previous sensor installation.
- Failure to change the oxygen sensor on time will result in low/no regulation of O<sub>2</sub> concentration.
- Warranty void if wrong/ non-original sensor is used.

Please refer to the service manual for replacement instructions.

### 37.6 CO<sub>2</sub> sensor

The CO<sub>2</sub> regulation uses the CO<sub>2</sub> sensor to keep the gas concentration at the chambers' desired level.

This sensor's lifetime is more than 6 years, but for safety reasons, Esco Medical recommends the sensor to be replaced once every 4-years.

Please follow these safety precautions when changing the sensor:

- Always use an original CO<sub>2</sub> sensor (contact Esco Medical or your local distributor for more details or ordering).
- Change the CO<sub>2</sub> sensor within 4 years from the date of installation.
- Failure to change the CO<sub>2</sub> sensor on time can result in low/no CO<sub>2</sub> gas concentration regulation.
- Warranty void if wrong/ non-original sensor is used.

Please refer to the service manual for replacement instructions.

### 37.7 UV light (applicable only for the Mini MIRI® model)

For safety reasons and to clean the re-circulating air, this equipment has a 254 nm UV light installed. The UV-C light has a limited lifetime and must be replaced every year, according to table 37.1.



Figure 37.1 UV light warning



**Exposure to UV-C radiation may cause severe damage to your skin and eyes. Always power off before removing any cover.**

Please follow these safety precautions when changing the UV-C light:

- Always use an original UV-C light bulb (contact Esco Medical or your local distributor for more details or ordering).
- Change UV-C light bulb within 1 year from date of installation.
- Failure to change the UV-light bulb on time can result in contamination build-up.
- Warranty void if wrong/ non-original UV-light bulb is used.

Please refer to the service manual for replacement instructions.

### 37.8 Cooling fan

The cooling fan is responsible for cooling down the electronics installed in the unit. A breakdown of the cooling fan will stress the components due to temperature rise within the system. It may cause the electronics to drift, resulting in low temperature and gas regulation.

To avoid this, Esco Medical recommends that the cooling fan be replaced once every 3 years.

Please follow these safety precautions when changing the cooling fan:

- Always use an original fan (contact Esco Medical or your local distributor for more details or ordering).
- Change the fan within 3 years from the date of installation.
- Failure to change the fan may cause the electronics to drift, resulting in low temperature and gas regulations.
- Warranty void if wrong/ non-original fan is used.

Please refer to the service manual for replacement instructions.

### 37.9 Internal gas pump (applicable only for the Mini MIRI® model)

The internal gas pump is used to transport the mixed gas through the VOC/HEPA filter, UV light and the chambers. In time the performance of this pump can be affected, causing a longer recovery time.

Therefore, this pump must be replaced once every 2 years to maintain the fast recovery time after lid openings.

Please follow these safety precautions when changing the internal gas pump:

- Always use an original gas pump (contact Esco Medical or your local distributor for more details or ordering).
- Change the gas pump within 2 years from the date of installation.
- Failure to change the pump may cause slow recovery times or breakdowns.
- Warranty void if wrong/ non-original pump is used.

Please refer to the service manual for replacement instructions.

### 37.10 Pump module (applicable only for the Mini MIRI® Humidity model)

The pump is used to transport the mixed gas through the chambers. In time the performance of this pump can be affected, causing a longer recovery time.

Therefore, this pump must be replaced once every 2 years to maintain the fast recovery time after lid openings.

Please follow these safety precautions when changing the internal gas pump:

- Always use an original gas pump (contact Esco Medical or your local distributor for more details or ordering).
- Change the gas pump within 2 years from the date of installation.
- Failure to change the pump may cause slow recovery times or breakdowns.
- Warranty void if wrong/ non-original pump is used.

Please refer to the service manual for replacement instructions.

### 37.11 Proportional valves

The internal valves make gas regulation possible. If the proportional valves are worn, gas regulation may be affected. It may cause more prolonged recovery time, incorrect gas concentration or breakdown. Therefore, these proportional valves must be replaced once every 3 years to maintain system safety and stability.

Please follow these safety precautions when changing valves:

- Always use original proportional valves (contact Esco Medical or your local distributor for more details or ordering).
- Change the valves within 3 years from the date of installation.
- Failure to change the valves may cause slow recovery times or breakdowns.

- Warranty void if wrong/ non-original valves are used.

Please refer to the service manual for replacement instructions.

### 37.12 Gas lines

The internal gas lines are used to transport mixed gas through the VOC/HEPA filter, UV light and the chambers. Over time, the lines' performance can be affected, causing more extended recovery time due to clogging.

Therefore, the gas lines must be replaced once every 3 years to maintain the fast recovery time after lid opening.

Please follow these safety precautions when changing gas lines:

- Always use original gas lines (contact Esco Medical or your local distributor for more details or ordering).
- Change the gas lines within 3 years from the date of installation.
- Failure to change the gas lines may cause slow recovery times or breakdowns.
- Warranty void if wrong/ non-original gas lines are used.

Please refer to the service manual for replacement instructions.

### 37.13 Flow sensors

The flow sensors are used by the CO<sub>2</sub>/N<sub>2</sub> regulations and for logging the unit's gas consumption.

This sensor's lifetime is more than 3 years, but Esco Medical recommends the sensor be replaced once every 2<sup>nd</sup> year for safety reasons.

Please follow these safety precautions when changing sensors:

- Always use an original flow sensor (contact Esco Medical or your local distributor for more details or ordering).
- Change flow sensors within 2 years from the date of installation.
- Failure to change the flow sensors on time may result in low/no CO<sub>2</sub> and O<sub>2</sub> gas concentration regulation.
- Warranty void if wrong/ non-original sensors are used.

Please refer to the service manual for replacement instructions.

### 37.14 Pressure regulators

The internal pressure regulators protect the system from too high external gas pressures that would damage the gas circuit's sensitive parts. If the pressure regulators are worn, they may begin to drift and not offer the protection they are supposed to. It could cause breakdowns or leaks in the internal gas circuit. Therefore, the regulators must be replaced once every 4<sup>th</sup> year to maintain the system safe and stable.

Please follow these safety precautions when changing regulators:

- Always use original pressure regulators (contact Esco Medical or your local distributor for more details or ordering).
- Change the regulators within 4 years from the date of installation.
- Failure to change the regulators may cause breakdowns.
- Warranty void if wrong/ non-original regulators are used.

Please refer to the service manual for replacement instructions.

### 37.15 Internal 0.2μ filter for CO<sub>2</sub> gas

The round-shape in line 0.2μ HEPA filter for CO<sub>2</sub> gas removes any particles found in the incoming CO<sub>2</sub> gas. Failure to use the HEPA filter may cause damage to the high precision flow sensor, calculate the amount of CO<sub>2</sub> gas entering the system, and disturb the CO<sub>2</sub> regulator system.

Please follow these safety precautions when changing the filter:

- Always use the original filter (contact Esco Medical or your local distributor for more details or ordering).
- Change the filter once every year.
- Failure to change the filter on time will result in low/no cleaning of incoming CO<sub>2</sub> gas.
- Warranty void if wrong/ non-original filter is used.

Please refer to the service manual for replacement instructions.

### 37.16 Internal 0.2μ filter for N<sub>2</sub> gas

The round-shape in line 0.2μ HEPA filter for N<sub>2</sub> gas removes any particles found in the incoming N<sub>2</sub> gas. Failure to use the HEPA filter may cause damage to the high precision flow sensor, calculate the amount of N<sub>2</sub> gas entering the system, and disturb the N<sub>2</sub> regulator system.

Please follow these safety precautions when changing the filter:

- Always use the original filter (contact Esco Medical or your local distributor for more details or ordering).
- Change the filter once every year.
- Failure to change the filter on time will result in low/no cleaning of incoming N<sub>2</sub> gas.
- Warranty void if wrong/ non-original filter is used.

Please refer to the service manual for replacement instructions.

### 37.17 Firmware update

If Esco Medical has released a newer version of the firmware, this should be installed on the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubators during the yearly scheduled service.

Please refer to the service manual for instructions on how to update the firmware.

## 38 The Installation guide

This section describes when and how to install the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator in the IVF clinic.

### 38.1 Responsibilities

All technicians or embryologists installing the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator must identify problems and perform any necessary calibrations, adjustments and maintenance.

Installation personnel performing MEA (Mouse Embryo Assay) must be thoroughly familiar with the MEA and all functions of the instrument, calibration and testing procedures, and instruments used in the instrument's testing. MEA test is a supplemental installation test and is not mandatory.

All individuals who will perform installation, repair and/or maintenance of the instrument must be trained by Esco Medical or at a qualified training center. Experienced service technicians or embryologists conduct training to ensure that the installation personnel clearly understand the instrument's functions, performance, testing, and maintenance.

Installation personnel must be updated regarding alterations or additions to this document and the "Installation report" form.

## 38.2 Before installation

2 – 3 weeks before the installation due, the user/owner at the clinic is contacted via e-mail to plan the exact time to perform the installation. When a convenient time has been determined, travel and accommodation arrangements can be made.

The released Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator must be sent 1 – 3 weeks before installation, depending on the clinic location. Check with shippers about local customs regulations and delays that could arise from that.

The clinic must be informed about the site requirements before installation and should have signed the customer requirement checklist:

1. The lab must have an idle sturdy and stable lab bench for standing operation.
2. The Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator weight is approx. 22 kg.
3. The required space for placement is 1.0 m x 0.6 m.
4. Temperature control should be able to maintain a stable temperature, never exceeding 30 °C.
5. Humidity control to avoid condensation.
6. Uninterrupted power supply (UPS) with 115 or 230 V, minimum 120 W.
7. Proper grounding.
8. CO<sub>2</sub> gas outlet with 0.6 – 1.0 atm above ambient.
9. N<sub>2</sub> gas outlet with 0.6 – 1.0 atm above ambient if the clinic uses reduced oxygen levels
10. Tubes that fit 4 mm hose nipple and HEPA filter.
11. Access to a PC with USB for the data logging.

## 38.3 Preparing for installation

- Print out the installation test from the validation manual. Make sure it is the latest and current version only.
- Fill out the following blank boxes in the form: the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator serial number (S/N) and customer.
- The service tool kit is checked for content before every installation trip to ensure it contains the necessary tools.
- Always bring the latest versions of firmware and data logging software. Bring these files on a labeled memory stick to the service site.

## 38.4 Bring the following to the installation site

- "Installation report" form.

- Service manual for the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubators.
- Updated service tool kit.
- Memory stick with the latest released firmware & software.
- High precision thermometer with a resolution not less than 0.1 °C.
- Calibrated gas analyzer with precision at least 0.1% for CO<sub>2</sub> and O<sub>2</sub> and the possibility of returning gas samples to the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator.
- Extension cable for USB connection.

### 38.5 Installation procedure at the site


1. Follow the guidelines in the safety instructions and warnings section ("2 Safety warning" section).
2. Connect the power cable to the UPS.
3. Connect the power cable to the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator.
4. Connect the gas lines.
5. Set gas pressure on the external gas regulator at 0.4 – 0.6 bar (5.80 – 8.70 PSI).
6. Switch on the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator on the back.
7. Observe for standard functionality.
8. Let the unit warm up and stabilize for 30 min.
9. Follow the guidelines in the "33 The Validation guide".
10. Complete user training and finish reading instructions.
11. After a burn-in phase of 24-hours, the unit is ready for use IF the testing is successful.


### 38.6 User training

1. Mains switch on/off.
2. Explain the Mini MIRI® Dry and Mini MIRI® Humidity incubators essential function and incubation with a multi-room facility to store the samples.
3. Explain temperature control in the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator (direct heat transfer with heated lids).
4. Gas regulation on/off.
5. Setpoint for temperature, CO<sub>2</sub> and O<sub>2</sub>.
6. Explain how N<sub>2</sub> is used to suppress the O<sub>2</sub> level.
7. Alarm turn off procedure (temperature, CO<sub>2</sub>, O<sub>2</sub>) and revert times.
8. Insertion and removal of heating inserts from the Mini MIRI® Dry and Mini MIRI® Humidity multiroom IVF incubator.
9. How to toggle the "Oil culture" and "Open culture" modes, and when which mode should be used.



10. Emergency procedures (can be found in the "29 Emergency Procedures" section).
11. Explain how to clean the device and heating inserts.
12. External measurement and calibration of temperature.
13. External measurement and calibration of gas concentration (InControl or Galaxy).
14. How to add and remove a sample.
15. Demonstrate how to replace the VOC-HEPA filter (can be found in the "12.1 Installation of new filter capsule" section). Not applicable in the Mini MIRI® Humidity multiroom IVF incubator.
16. Datalogger functionality, how to establish a connection and re-connection.


 **Use the User manual section as much as possible to get the User well acquainted with it.**

 **The user/owner is informed that the first VOC filter change is 3 months after installation and 3-month intervals. The first service check is under normal circumstances after 1 year.**

### 38.7 After the installation

When the installation trip is finished, a copy of the original "Installation report" form must be sent to Esco Medical Ltd. It will be saved with the device records. According to the ISO procedure and Medical Device Directive, a paper copy of the completed and signed installation test form is stored in the unique device's device history record. The date of installation is written in the instrument overview file. The date of installation is also written in the service schedule.

Suppose the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator user or owner make inquiries about a written "Installation report". The completed and signed "Installation report" form must be sent to the clinic. Any deviations/complaints/suggestions from the Installation visit are reported in the CAPA system. If a critical error has occurred, information about this will be reported directly to QC or QA.

 **If the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator fails any of the "Installation report" form acceptance criteria, or it in any way suffer from a severe error and incubation parameters are compromised, the Mini MIRI® Dry or Mini MIRI® Humidity multiroom IVF incubator must be taken out of service until it is repaired/ exchanged, or anew test approves the Mini MIRI® or Mini MIRI® Humidity multiroom IVF incubator. The User and owner must be informed about this and arrangements to solve the problems must be initiated.**