

MIRI[®] GA

Gas and Temperature Validation Unit

NEXT GENERATION Desktop CO₂ / O₂ & Temperature Validation Unit

MIRI[®] GA is a tabletop device intended to make external incubator validation easier and safer. It is capable of monitoring the temperature (PT1000 connector) & gas concentration, flow and pressure. It can validate up to 6 chambers simultaneously 24 hours a day. The MIRI[®] GA can connect to any brand of incubator and is a perfect accessory to MIRI[®] TL and MIRI[®] Multiroom Incubators.



Key Features



Constantly validate up to 6 x CO₂ / O₂ incubators




Controllable flow rate



Monitor up to 6 x PT1000 sensors



6 ports for sequential gas samples

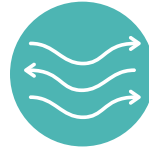

Gas feedback returns sampled gas to incubator or exhaust



Audible and visual alarms



As validation unit, MIRI[®] GA can provide

- Parameter data logging
- Adjustable alarm thresholds

Connected to a PC, the logging software will document parameters and optionally send e-mail alarms if set up.



Want to know more about this product? Scan the QR Code

Esco Micro Pte. Ltd.

21 Changi South Street 1 • Singapore 486 777 • Tel +65 6542 0833 • Fax +65 6542 6920
medical@escoglobal.com • esco-medical.com

Esco Global Offices: Bangladesh | China | Denmark | Germany | Hong Kong | Indonesia | Italy | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam