


# CultureCoin®

Culture Dish for MIRI® TL

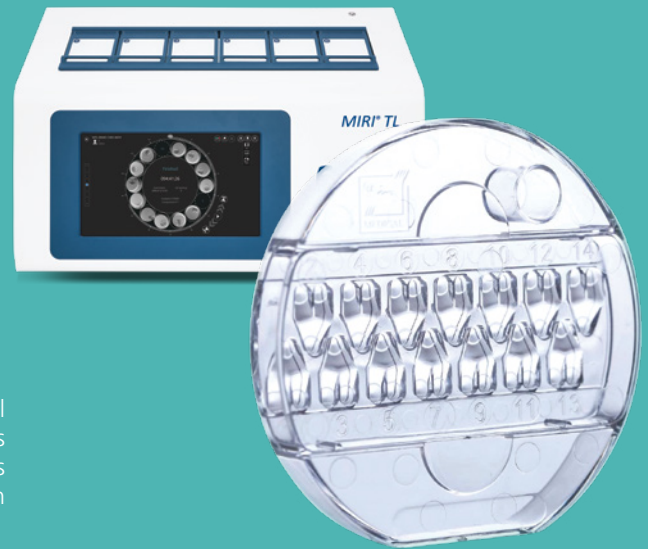

*A culture dish, exclusively designed for the MIRI® TL*

The Esco Medical CultureCoin® is a sterile, single use, culture vessel specifically developed for use in IVF. The dish holds up to 14 embryos in the centrally placed line of culture wells. The culture well features a 300µ, optically clear area where the embryo is situated and is an area optimized for microscopy.

## Key Features

### Numbered culture wells (1-14)

For easy location of each well by the user. Each culture well has two washing wells for cell manipulations, washing media or replacement media.

### Mounted lid

located inside the circumference of the dish covering the culture area completely.

### pH reservoir

For pH measuring without affecting the cultured samples.

### Flat bottom

Ideal for heat transfer from below.



Ergonomic design for easy, safe, and secure handling of embryos.





Gamma-sterilized.



Oxygen plasma treated surface for the effective prevention of bubble formation.

## Parts of the CultureCoin®

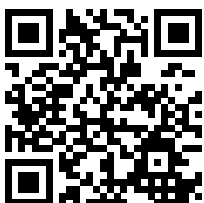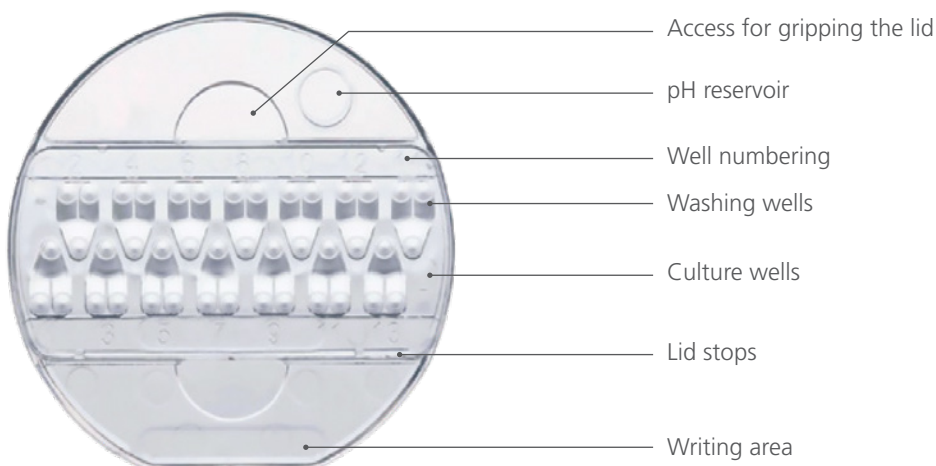

**Want to know more about this product? Scan the QR Code**

**Esco Micro Pte. Ltd.**

21 Changi South Street 1 • Singapore 486 777 • Tel +65 6542 0833 • Fax +65 6542 6920  
medical@escoglobal.com • esco-medical.com

**Esco Global Offices:** Bangladesh | China | Denmark | Germany | Hong Kong | Indonesia | Italy | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam