

# MIRI®

## Multiroom Incubator

*VIP treatment for high success rate embryos in IVF*  
The MIRI® is a revolution, in form and functionality, of CO<sub>2</sub> incubators for *In Vitro* Fertilization (IVF). With 6 chambers, the MIRI® is a Multiroom Incubator that allows users to access their cultures in one chamber without affecting the neighbouring chambers. Thus, the harmful effects of fluctuations in temperature and gas caused by frequent incubator access are avoided. Built specifically to equip IVF laboratories and clinics to provide the best standard of care, it boasts a unique set of features that cannot be found elsewhere.



### Key Features



#### Heated Lid

Prevents condensation and enhances temperature regulation and recovery.



#### Heated bottom

Provides direct heat transfer to the cultures through the optimization plate for stable heat regulation.




#### Removable Heating Optimization Plate

With wide selection of inserts.



#### Multiroom System

Completely individual chambers for easier calibration, faster recovery, less disruption, and prevents cross-contamination.

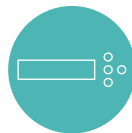

#### Control Panel Buttons and LED Display

A display that can be easily seen from a distance and a control panel that allows easy and intuitive operation.




#### Optional SAFE Sens Integration

For continuous pH monitoring.



#### Mute Button

Temporarily mutes alarm messages and sound for five (5) minutes.

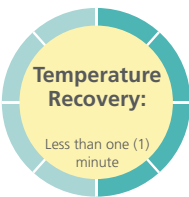

#### Fast Recovery

- 1 minute temperature recovery.
- 3 minutes gas recovery.

#### Built-in pH meter

For accurate validation.

#### Solid Validation System

- Six (6) PT1000 sensors and Gas ports for validation outputs
- SAFE Sens Integration for pH monitoring (optional).
- External Data Logging.
- Alarm relay contact
- Total capacity of up to 48 standard culture dishes.
- Separate CO<sub>2</sub> and O<sub>2</sub> regulation, expensive mixed gases not required!
- Air is continuously cleaned by HEPA/VOC filters, and UV light.

*This equipment is a CE-marked device and is in conformity with the essential requirements of the Medical Devices Directive 93/42/EEC as amended by Directive 2007/147/EC.*



**Want to know more about this product? Scan the QR Code**

**Esco Micro Pte. Ltd.**

21 Changi South Street 1 • Singapore 486 777 • Tel +65 6542 0833 • Fax +65 6542 6920  
medical@escoglobal.com • esco-medical.com

**Esco Global Offices:** Bangladesh | China | Denmark | Germany | Hong Kong | Indonesia | Italy | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam