

CultureCoin®

Culture Dish for MIRI® TL

A culture dish, exclusively designed for the MIRI® TL

The Esco Medical CultureCoin® is a sterile, single use, culture vessel specifically developed for use in IVF. The dish holds up to 14 embryos in the centrally placed line of culture wells. The culture well features a 300µ, optically clear area where the embryo is situated and is an area optimized for microscopy.

Key Features

Numbered culture wells (1-14)

For easy location of each well by the user. Each culture well has two washing wells for cell manipulations, washing media or replacement media.

Mounted lid

located inside the circumference of the dish covering the culture area completely.

pH reservoir

For pH measuring without affecting the cultured samples.

Flat bottom

Ideal for heat transfer from below.



Ergonomic design for easy, safe, and secure handling of embryos.





Gamma-sterilized.



Oxygen plasma treated surface for the effective prevention of bubble formation.

Parts of the CultureCoin®

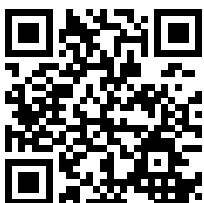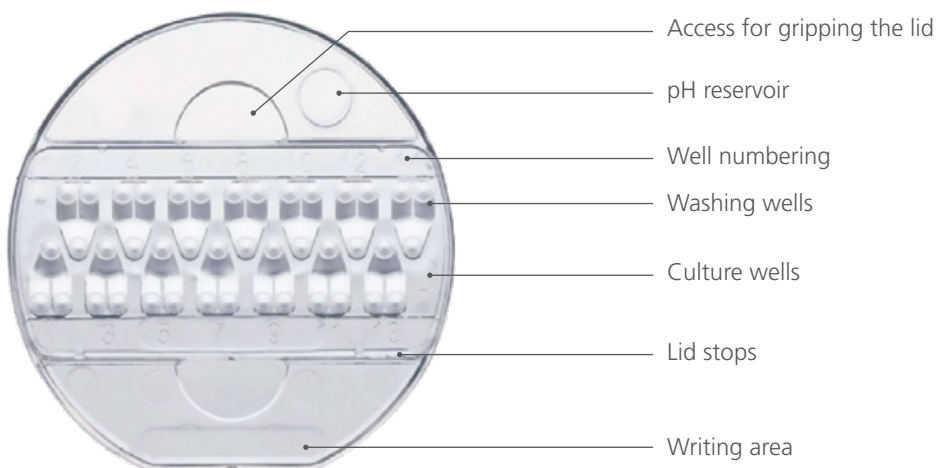

Want to know more about this product? Scan the QR Code

Esco Micro Pte. Ltd.

21 Changi South Street 1 • Singapore 486 777 • Tel +65 6542 0833 • Fax +65 6542 6920
medical@escoglobal.com • esco-medical.com

Esco Global Offices: Bangladesh | China | Denmark | Germany | Hong Kong | Indonesia | Italy | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam