

Mini MIRI®

Dry and Humidity Incubators

A compact incubator ensuring unique environment for human embryo culture

Built on the strong and reliable MIRI® Multiroom, the Mini MIRI® is an incubator that provides a stable culture environment. It has two chambers that prevent cross-contamination while VOC/HEPA filtration cleans the incoming airstream. The compact design and direct heat regulation further translate to faster temperature and gas recovery.



Key Features



Ergonomic Rotatory Key

This is used to access the menu, to toggle among settings and to set parameter values.



Built-In Gas Mixer

Premixed gas is not required (but possible to use).



Heated Lid

Prevents condensation and improves temperature regulation and recovery enhancing uniformity within each chamber.



Dual Chamber System


This is ideal because any disruption (e.g. temperature drop after opening lid) has zero impact on the rest of the system. Also, calibration is much simpler since there is no crossover of heat from adjacent chamber. While the chambers operate independently, setpoints are common across all chambers rather than individual.

Mini MIRI® comes in two variants to fit your culture's environmental needs.

Mini MIRI® Dry



Airflow Diagram

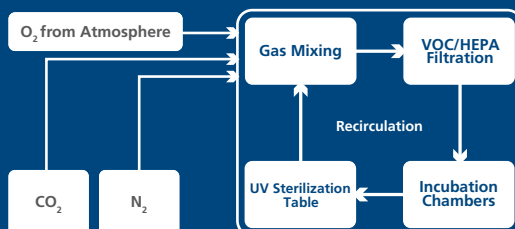

High quality airstream via VOC/HEPA filter + UV

The filter module can be easily replaced once used.


The gas in the MIRI® is continuously recirculated through a VOC/HEPA filter. A UV-C light (254 nm) sterilizes the airstream before passing through the filter.

Mini MIRI® Humidity

Passive humification system.



Airflow Diagram

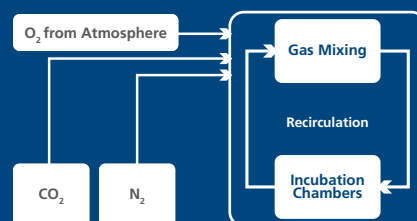


This equipment is a CE-marked device and is in conformity with the essential requirements of the medical devices EU regulation 2017/745.

Want to know
more about
this product?

Scan the
QR Code

www.esco-medical.com

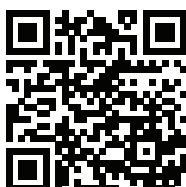

Esco Medical, Aps

Kringelld 10, 8250 Egaa, Denmark
medical@escolifesciences.com
+65 6542 0833

Esco Global Offices: Bangladesh | China | Denmark | Germany | Hong Kong | Indonesia | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam

9010450_Mini MIRI Flyer_A4_vG_101725

Esco can accept no responsibility for possible errors in catalogues, brochures and other printed materials. Esco reserves the right to alter its products and specifications without notice. All trademarks and logotypes in this material are the property of Esco and the respective companies.

Esco Micro, Pte. Ltd. (Headquarters)

19 Changi South Street 1, Singapore 486 779
Tel +65 6542 0833