



Enabling better fertility outcomes with first-in-class solutions

Fertility

Product Catalogue

table of contents

About Esco Medical04

MIRI® TL Incubator06

CultureCoin®07

Multiroom Incubator08

 MIRI®08

 MIRI® II-1209

Mini MIRI® Dry and Humidity Incubators11

Inserts for MIRI®, MIRI® II-12 and Mini MIRI®12

CelCulture® CO₂ Incubator13

Esco Multi-Zone ART Workstation14

Esco Workstation with MIRI® Chambers15

Stereomicroscope16

MIRI® AVT22

Versati™ Tabletop Centrifuge23

Quality Assurance and Validation Units24

 MIRI® GA24

 SAFE Sens25

Traceability Tool27

 MIRI® Evidence27

MIRI® GA
Gas and Temperature
Validation Unit



MIRI® AVT
Anti-Vibration Table



Esco Multi-Zone
ART Workstation



MIRI® Laminar Flow



Mini MIRI®
Dry and Humidity Incubators



CellCulture®
CO₂ Incubator



MIRI® Evidence
RFID Witnessing System



MIRI® Humidity
Multiroom Incubator



MIRI®
Multiroom Incubator



Versati™
Tabletop Centrifuge



MIRI® II - 12
Multiroom Incubator



CultureCoin®



MIRI® TL6
Time-Lapse Incubator



MIRI® TL12
Time-Lapse Incubator



About Esco



Welcome to Esco

Since the establishment of Esco in 1978, we never stopped developing, providing, and delivering innovative solutions. From one, we have progressed into four business units with a worldwide presence, namely Esco Scientific, Esco Healthcare, Esco Medical, and Esco Aster—remaining true to our tagline “World-class. Worldwide.”

Last 2020, we shifted from Esco Group of Companies to **Esco Lifesciences Group**, carrying a new tagline **“Improving lives through science.”** The transformation of the company name and brand signifies Esco’s vigor in keeping up, responsive, and adaptive with the fastchanging world while keeping focused on its mission to deliver enabling technologies and provide service all over the world—and improve lives through science.

At Esco Medical, life has begun

Esco Medical is one of the divisions of the Esco Lifesciences Group, apart from Esco Scientific, Esco Healthcare and Esco Aster. Esco Medical provides innovative technological solutions for fertility clinics and laboratories.

The slightest deviation, usually considered as insignificant, often result in non-optimal conditions for embryo growth and lowered pregnancy success. In Esco, we understand that even the smallest details affect the *In vitro* Fertilization process. Thus, Esco Medical’s primary focus is to provide fertility technologies and solutions to help the world’s leading IVF centers to improve, standardize and automate their processes in order to achieve better clinical outcomes and patient satisfaction.

Esco Medical is the leading manufacturer and innovator of high-quality equipment such as Time-Lapse Incubator, Multiroom Embryo Incubators, IVF Workstation, CO₂ Incubator, Anti-Vibration Table, and Gas Analyser. Most of our medical products are designed in Denmark and made in the EU.

*I promise...
I'll raise you like flying a kite.
Will hold you with both hands and slowly let go.
The day you fly, I'll watch you chase the sky.*



MIRI® Time-Lapse Incubator



MIRI® TL is a Time-Lapse incubator that monitors embryo development. The MIRI® TL, optimized for clinical and IVF procedures, is designed to support existing work and quality assurance routines. This value-added treatment provides the most unique incubation environment with the market's most secure and safest procedures. It lessens disturbance and minimizes stressful factors that may be introduced when taking the dishes out of the incubator. This incubation system also ensures predictability in the daily handling and currently offers the market's lowest cost of ownership.



Unique Incubation Environment

- Has independent multi-chamber system
- Gas recirculation through HEPA/VOC filters and UV light.
- Built-in gas mixer. Premixed gas is not required.

MIRI® TL6: 6 Individual chambers

MIRI® TL12: 12 Individual chambers


Gas recovery: less than three (3) minutes

Temperature recovery: less than one (1) minute



Unprecedented Faster Recovery

- Excellent recovery time for both temperature and gas parameters.
- Opening one chamber will have no impact on the rest of the system.
- Heated upper lid and bottom plate for excellent temperature regulation and uniformity.

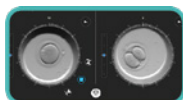

2 Temperature Mode Options:

- **Single:** Uniform set points for all 6/12 (six/twelve) chambers.
- **Multi:** Individual set points for each chamber.

Sophisticated Annotation Tools


- Freedom to personalize instrument and parameter settings.
- Do a side-by-side comparison and compare actual timings to ideal.



Quality checking an easy breeze!

- Has 6/12 temperature sensors to ensure constant temperature stability.
- Independent PT1000 sensor and gas sample port for external validation for each chamber.
- Built-in pH measuring system.
- Data logging system.

This equipment is a CE-marked device and is in conformity with the essential requirements of the medical devices EU regulation 2017/745.

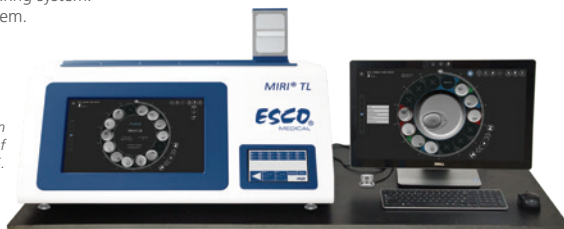

Embryo Analysis and Evaluation System

The MIRI® TL Viewer Software is a simple yet sophisticated information-providing tool that can help embryologists process the data generated. You can review, annotate and compare the morphokinetic parameters of each embryo to select or deselect embryos for transfer and export data for retrospective analysis.



Navigation through the stacked timeline is easy and intuitive as the revolver shows the videos of the 14 wells of one single CultureCoin®. You can play the individual videos, annotate and compare each single embryo.

Shown on the image is a magnified view of embryo #3

General Specifications

Specifications	TL6	TL12
Overall Dimensions (W x D x H)	805 x 590 x 375 mm (31.7 x 23.2 x 14.8")	950 x 685 x 375 mm (37.4 x 27.0 x 14.8")
Compartment Dimensions	120 x 90 x 26 mm (4.7 x 3.5 x 1")	
Temperature Control Range	28.7 - 41.0 °C	
Gas Consumption (CO ₂) *	< 2 L/h	
Gas Consumption (N ₂) **	< 5 L/h	
CO ₂ Control Range	2.9% - 9.9%	
O ₂ Control Range	2.0% - 20.0%	
Input Gas Pressure	0.4 – 0.6 bar (5.80 – 8.70 PSI)	
Built-in Microscope	Zeiss 20x, objective has numerical aperture of 0.35, specialized for 640 nm illumination	
Embryo Illumination	0.064s per image, using 1W single red LED (635nm)	
Camera Resolution	1920 x 1200. Monochrome, 12-bit, IDS system	
Optics Tube Ratio	3.00 px/µm	
Imaging Focal Planes	5, 10 and 20min intervals in 3, 5 and 7 focal planes	
Number of pixels in stored image	670 x 670	860 x 860

* Under normal condition (CO₂ set point reached at 6.0%, all lids closed).

** Under normal condition (O₂ set point reached at 5.0%, all lids closed).

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
Unit		
2070091	MRI-TL-MN-6C-8	MIRI® Time-Lapse Incubator, Mini, 6 Chambers, 230 V, 50/60 Hz
2070092	MRI-TL-MN-6C-9	MIRI® Time-Lapse Incubator, Mini, 6 Chambers, 115 V, 50/60 Hz
2070098	MRI-TL-MN-6C-SS-8	MIRI® Time-Lapse Incubator, 6 chambers with SAFE Sens, 230 V, 50/60 Hz
2070099	MRI-TL-MN-6C-SS-9	MIRI® Time-Lapse Incubator, 6 chambers with SAFE Sens, 115 V, 50/60 Hz
2070100	MRI-TL-12C-8	MIRI® Time-Lapse Incubator, 12 Chambers, 230 V, 50/60 Hz
2070101	MRI-TL-12C-9	MIRI® Time-Lapse Incubator, 12 Chambers, 115 V, 50/60 Hz
2070114	MRI-TL-SS-12C8	MIRI® Time-Lapse Incubator, 12 Chambers, 230V 50/60Hz, with SAFE Sens
2070115	MRI-TL-12C-SS-9	MIRI® Time-Lapse Incubator, 12 Chambers, 110V 50/60Hz, with SAFE Sens
MIRI® TL Viewer Software		
2070042	MRI-VIEWER	MIRI® Time-Lapse Viewer
1320095	MRI-SERVER	MIRI® Time-Lapse Server
Accessories		
1320011	MRA-1007	HEPA + VOC filter (to be replaced every 3 months)
1320088	MRI-CC	CultureCoin® for Time-Lapse of 14 embryos (25 pcs. per pack)
1320045	MRI-GA	MIRI® GA CO ₂ /O ₂ & Temperature Validation Unit, 115V/ 230V



CultureCoin® for MIRI® TL

- Holds up to 14 embryos with individual numbered wells (1-14).
- For single and separated culture where each embryo are cultured in its own environment.
- Ergonomic design for easy handling and location of embryos.
- Separate well for pH measurements.
- Corona plasma treated surface for the effective prevention of bubble formation.
- Packed in 1 dish pouches and delivered in boxes of 25 pcs.

General Specifications

Overall dimensions (Diameter x Height)	Ø 71 x 10 mm
CultureCoin® weight in total	13.8 grams
Material	Styrene Methyl Methacrylate (SMMA)
Incubation Temperature Range	28.7 - 40.0 °C
Incubation CO ₂ Range	1.9 - 10.0%
Incubation O ₂ Range	4.9 – 20.0%
Sterilization Method	Gamma Irradiation
Lifetime	2 years
Biocompatibility Tests	<ol style="list-style-type: none"> 1. Mouse Embryo Assay (MEA) test with thawed 1-cell mouse embryos. Acceptance criteria: at least 80% of embryos developed to the blastocyst stage. 2. Limulus Amebocyte Lysate (LAL) test. Acceptance criteria: < 20 EU/device.

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
1320088	MRI-CC	CultureCoin® for Time-Lapse of 14 embryos

Multiroom Incubator

MIRI® Incubation System

The Top-of-the-Line Features of the MIRI® Incubation System

- **Heated Lid**
Prevents condensation. Enhances temperature regulation.
- **Direct Heat Transfer**
Less than one minute temperature recovery.
- **Completely Independent Chambers**
Any disruption (e.g., temperature drop after opening the lid) has zero impact on the rest of the system.
- **A Complete Incubation Environment**
 - Has a built-in gas mixer. Premixed gas is not required.
 - Built-in pH measuring system and data logging system.



MIRI® Multiroom Incubator

The MIRI® is a revolution, in form and functionality, of benchtop incubators for *In Vitro* Fertilization (IVF). With 6 chambers, the MIRI® is a Multiroom Incubator that allows users to access their cultures in one chamber without affecting the neighbouring chambers. Thus, the harmful effects of fluctuations in temperature and gas caused by frequent incubator access are avoided. Built specifically to equip IVF laboratories and clinics to provide the best standard of care, it boasts a unique set of features that cannot be found elsewhere.

Key Features

Fast Recovery


- <1 minute temperature recovery.
- <3 minutes gas recovery.

Built-in pH meter

For accurate validation.

Solid Validation System

- Six (6) PT1000 sensors and Gas ports for validation outputs
- SAFE Sens Integration for pH monitoring (optional).
- External Data Logging.
- Alarm relay contact
- Total capacity of up to 48 standard culture dishes.
- Separate CO₂ and O₂ regulation, expensive mixed gases not required!
- Air is continuously cleaned by HEPA/VOC filters, and UV light. (For MIRI® Dry only)

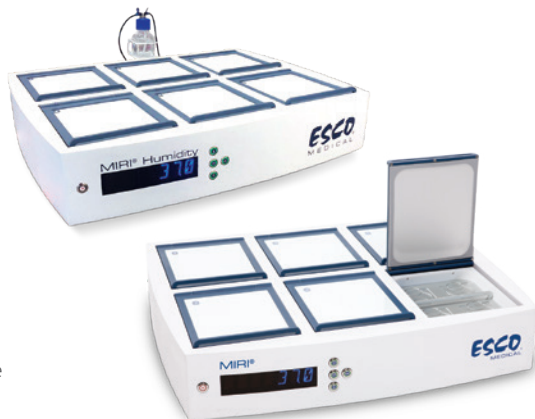

This equipment is a CE-marked device and is in conformity with the essential requirements of the medical devices EU regulation 2017/745.

Stacking Frames



MIRI® Stacking Frame with drawer



MRA-1014 - Stacking frame for 2 units

MIRI® II-12 Multiroom Incubator



The MIRI® II-12 is an incubator that provides unique features for every IVF laboratories and clinics. The chambers are specially designed to accommodate one patient ensuring personal space for each embryo. Having an excellent footprint, MIRI® II-12 is made to perfectly fit every IVF lab.

Independent Chambers

Each chamber is specially designed for one patient. Hence, there is no disturbance to other chambers even when a lid is opened/closed.

Excellent footprint

With its compact size, it can perfectly fit in every IVF Lab.

Low gas consumption

The built-in gas mixer and efficient recirculation system allows you to save more.



This equipment is a CE-marked device and is in conformity with the essential requirements of the medical devices EU regulation 2017/745.

Just a fitting solution...

MIRI® II-12 comes with specific inserts matching the type of dishes used in the laboratories.

Heating optimization plates


Each chamber contains a heating optimization plate to facilitate heat transfer directly to the culture dishes.

- Has inserts to fit various dish sizes



The dishes fit into the inserts so that the heat is directly transferred to the media.

Airflow Diagram

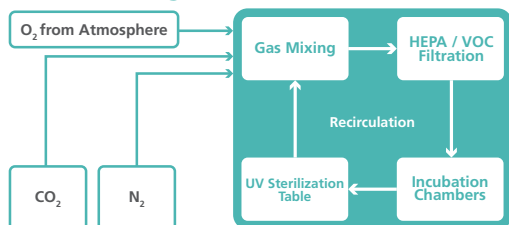

Provides total control of the gas phase environment. The built-in gas mixer and the high-performance CO₂ and O₂ sensors allow accurate control of gas composition in the chambers.

General Specifications



MIRI® Multiroom Incubator

Model	MIRI®	MIRI® Humidity
Overall Dimensions (W x D x H)	700 x 585 x 165 mm (27.6 x 23.0 x 6.5")	700 x 645 x 280 mm (27.6 x 25.4 x 11.0")
Compartment Dimensions	200 x 176 x 25 mm (7.9 x 6.9 x 1")	
Power Supply	115 / 230V, 50/60 Hz	
Power Consumption	300 W	
Temperature Control Range	24.9 – 40.0 °C	
*CO ₂ Gas Consumption	<2 L/h	<4 L/h
**N ₂ Gas Consumption	<12 L/h	
CO ₂ Control Range	2.0 – 9.9%	
O ₂ Control Range	5.0 – 20.0%	
Input Gas Pressure	0.4 – 0.6 bar (5.80 – 8.70 PSI)	
Net Weight	40 kg (88.2 lbs)	
Shipping Weight	50 kg (110.2 lbs)	
Shipping Dimension	860 x 760 x 460 mm (33.9" X 30" x 18.1")	

* Under normal condition (CO₂ set point reached at 5.0%, all lids closed)

** Under normal condition (O₂ set point reached at 5.0%, all lids closed)

MIRI® II-12 Multiroom Incubator

Overall Dimensions (W x D x H)	740 x 575 x 215 mm (29.1 x 22.6 x 8.5")
Compartment Dimensions	120 x 90 x 26 mm (4.7 x 3.5 x 1")
Power Supply	115 / 230V, 50/60 Hz
Power Consumption	500 W
Temperature Control Range	25.0 – 40.0 °C
*CO ₂ Gas Consumption	<2 L/h
**N ₂ Gas Consumption	<12 L/h
CO ₂ Control Range	3.0 – 10.0%
O ₂ Control Range	5.0 – 10.0%
Input Gas Pressure	0.4 – 0.6 bar (5.80 – 8.70 PSI)
Net Weight	47 kg
Shipping weight	63 kg (138.9 lbs)
Shipping dimension	900 x 720 x 460 mm (34.6 x 27.6 x 14")

* Under normal condition (CO₂ set point reached at 5.0%, all lids closed)

** Under normal condition (O₂ set point reached at 5.0%, all lids closed)

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
MIRI® Multiroom Incubator		
2070047	MRI-6A10-8	MIRI® Incubator, 230V, 50/60Hz
2070048	MRI-6A10-9	MIRI® Incubator, 115V, 50/60Hz
2070086	MRI-6A10-SS-8	MIRI® Incubator, with SAFE Sens for pH measurement, 230V, 50/60Hz
2070087	MRI-6A10-SS-9	MIRI® Incubator, with SAFE Sens for pH measurement, 115V, 50/60Hz
2070183	MRI-6A10-H-8	MIRI® Humidity 6 Chamber, 230V, 50/60Hz
2070184	MRI-6A10-H-9	MIRI® Humidity 6 Chamber, 115V, 50/60Hz
2070185	MRI-6A10-H-SS-8	MIRI® Humidity 6 chamber with SAFE Sens for pH monitoring, 230V, 50/60Hz
2070186	MRI-6A10-H-SS-9	MIRI® Humidity 6 chamber with SAFE Sens for pH monitoring, 115V, 50/60Hz
1320045	MRI-GA	MIRI® GA CO ₂ / O ₂ & Temperature Validation Unit, 115V / 230V
MIRI® II-12 Multiroom Incubator		
2070164	MRI2-12C-8	MIRI® II-12 Incubator with 12 chambers, 230V, 50/60Hz
2070165	MRI2-12C-9	MIRI® II-12 Incubator with 12 chambers, 115V, 50/60Hz
2070166	MRI2-12C-SS-8	MIRI® II-12 Incubator with 12 chambers, with SAFE Sens, 230V, 50/60Hz
2070167	MRI2-12C-SS-9	MIRI® II-12 Incubator with 12 chambers, with SAFE Sens, 115V, 50/60Hz

Mini MIRI® Incubator



Built on the strong and reliable MIRI® Multiroom, the Mini MIRI® is an incubator that provides a stable culture environment. It has two chambers that prevent cross-contamination while HEPA/VOC filtration cleans the incoming airstream. The compact design and direct heat regulation further translate to faster temperature and gas recovery.

Comes in two models:



Mini MIRI® Humidity

- The water bottle is located on the side of the unit for refilling and easy control of water level.
- Passive humidification system



Mini MIRI® Dry

- Has a built-in gas mixer. Premixed gas is not required
- Comes with a UV module and HEPA/VOC filter.

This equipment is a CE-marked device and is in conformity with the essential requirements of the medical devices EU regulation 2017/745.

General Specifications

Model	Mini MIRI® Dry	Mini MIRI® Humidity
Overall Dimensions (W x D x H)	525 x 420 x 230 mm (20.7 x 16.5 x 9.1")	
Compartment Dimensions	200 x 176 x 25 mm (7.9 x 6.9 x 1")	
Power Supply	115 / 230V, 50/60 Hz	
Temperature Control Range	24.9 – 40.0 °C	
CO ₂ Gas Consumption	<2 L/hr	< 4 L/h
N ₂ Gas Consumption	<8 L/h	<12 L/h
Input Gas Pressure	0.4 – 0.6 bar (5.80 – 8.70 PSI)	
CO ₂ Control Range	1.9 % – 9.9%	
O ₂ Control Range	3.9 % – 19.9%	
Net weight	22 kg (48.5 lbs)	
Shipping weight	30 kg (66.1 lbs)	
Shipping Dimensions	700 x 400 x 460 mm (27.6 x 15.7 x 18.1")	

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
2070155	MRI-MINI-H-8	Mini MIRI® Humidity, 230V, 50/60Hz
2070156	MRI-MINI-H-9	Mini MIRI® Humidity, 115V, 50/60Hz
2070157	MRI-MINI-H-SS-8	Mini MIRI® Humidity with SAFE Sens for pH monitoring, 230V, 50/60 Hz
2070158	MRI-MINI-H-SS-9	Mini MIRI® Humidity with SAFE Sens for pH monitoring, 115V, 50/60 Hz
2070143	MRI-MINI-D-8	Mini MIRI® Dry, without Humidification, 230V, 50/60Hz
2070144	MRI-MINI-D-9	Mini MIRI® Dry, without Humidification, 115V, 50/60Hz
2070145	MRI-MINI-SS-D-8	Mini MIRI® Dry with SAFE Sens for pH monitoring, 230V, 50/60 Hz
2070146	MRI-MINI-SS-D-9	Mini MIRI® Dry with SAFE Sens for pH monitoring, 115V, 50/60 Hz

Inserts for MIRI[®], MIRI[®] II-12, and Mini MIRI[®]



Extensive list of inserts for MIRI[®], MIRI[®] II-12, and Mini MIRI[®]

When ordering a MIRI[®], MIRI[®] II-12 or a Mini MIRI[®], you just have to pick the right insert/s that fit the dishes used in your laboratory. You have the freedom and flexibility to choose -- no limitations. The MIRI[®], MIRI[®] II-12 and Mini MIRI[®] can easily fit-in to your existing work routine.

All inserts are optimized for the direct transfer of heat to the dishes and are totally removable for easy cleaning. This is to ensure optimal conditions for your embryos.



Nunc™



Falcon®



Vitrolife



LifeGlobal®
GPS Dishes



Nipro™



SparMED Oosafe®

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
For MIRI[®] and Mini MIRI[®]		
1320003	MRA-FD	Insert for Falcon® Dishes
1320004	MRA-ND	Insert for Nunc™ Dishes
1320070	MRA-VD	Insert for Vitrolife Dishes
1320099	MRA-NID	Insert for Nipro™ Dishes
1320100	MRA-LD	Insert for LifeGlobal® GPS Dishes
1320101	MRA-PD	Insert Without Footprint for Plain Dishes
1320118	MRA-OD	Insert for SparMED Oosafe®
For MIRI[®] II-12		
1320429	MRA2-FD	Insert for Falcon® Dishes
1320430	MRA2-ND	Insert for Nunc™ Dishes
1320431	MRA2-VD	Insert for Vitrolife Dishes
1320433	MRA2-LD	Insert for LifeGlobal® GPS Dishes
1320436	MRA2-OD	Insert for SparMED Oosafe®
1320434	MRA2-PD	Insert for Plain Dish

CelCulture® CO₂ Incubators



The CO₂ Incubator has a vital role in providing an optimal environment in embryo development during IVF and other ART procedures. Sleek, reliable and intuitive, the Esco CelCulture® CO₂ incubator is packed with outstanding features such as rapid parameter recovery, ISO Class 5 Cleanliness, ISOCIDE™ antimicrobial coating, optional Inner Door Kit that reduces contamination risk, and other accessories for specialized applications.



CelCulture® CO₂ Incubators
available in 3 sizes, 50 L, 170 L, and 240 L.

CelCulture® CO₂ Incubators

CelCulture® is equipped with 90°C Moist Heat Decontamination System evaluated by HPA-UK. It utilizes ULPA filter to keep the chamber at ISO Class 5 cleanliness which ensures that all contaminants are filtered and clean air is recirculated.

Key Features

- Wider temperature range, from (Ambient +5) temperature to 60°C.
- Complete contamination control methods to protect your precious samples.
- All gas inputs are filtered via 0.2micron in-line filter and ULPA filtration system.
- 90°C moist heat decontamination cycle, validated by HPA-UK.

ISOCIDE™



Voyager Software Kit

PC-based software for remote monitoring, data logging and programming.



Floor Stand with Adjustable Feet

Nominal range of 180 mm to 250 mm (7.1" to 9.8")



Floor Stand with Casters

Support stand raises the incubator to a height of 700 mm (27.6") above the floor.



Roller Base

With casters for mobility of your incubators.

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
IR Sensor Model with Stainless Steel Chamber		
2170257	CCL-050B-8-IVF	Celculture® Incubator, 50 L, IR sensor, CO ₂ control, Moist Heat Decon, with Sealed Inner Door Kit for 50 L (2 Glass Doors with Latches), Factory Installed, 230 VAC, 50/60 Hz
2170272	CCL-170B-8-IVF	CelCulture® Incubator 170 L IR Sensor, CO ₂ Control ULPA, Moist Heat Decon, with Sealed Inner Door Kit for 170 L (4 Glass Doors with Latches), Factory Installed, 230 VAC 50/60 Hz
2170278	CCL-240B-8-IVF	CelCulture® Incubator 240 L IR Sensor, CO ₂ Control, ULPA, Moist Heat Decon, with Sealed Inner Door Kit for 240 L (4/6 Glass Doors with Latches), Factory Installed, 230 VAC 50/60 Hz
2170258	CCL-050B-9-IVF	Celculture® Incubator, 50 L, IR sensor, CO ₂ control, Moist Heat Decon, with Sealed Inner Door Kit for 50 L (2 Glass Doors with Latches), Factory Installed, 115 VAC, 50/60 Hz
2170273	CCL-170B-9-IVF	CelCulture® Incubator 170 L IR Sensor, CO ₂ Control ULPA, Moist Heat Decon, with Sealed Inner Door Kit for 170 L (4 Glass Doors with Latches), Factory Installed, 115 VAC 50/60 Hz
2170279	CCL-240B-9-IVF	CelCulture® Incubator 240 L, IR Sensor, CO ₂ Control, ULPA, Moist Heat Decon with Sealed Inner Door Kit for 240 L (4/6 Glass Doors with Latches), Factory Installed, 115 VAC 50/60Hz
Suppressed O₂ Model with Stainless Steel Chamber		
2170260	CCL-050T-8-IVF	Celculture® Incubator, 50 L, IR sensor, CO ₂ & O ₂ control, Moist Heat Decon, with Sealed Inner Door Kit for 50 L (2 Glass Doors with Latches), Factory Installed, 230 VAC, 50/60 Hz
2170275	CCL-170T-8-IVF	CelCulture® Incubator 170 L IR Sensor, CO ₂ & O ₂ Control ULPA, Moist Heat Decon, with Sealed Inner Door Kit for 170 L (4 Glass Doors with Latches), Factory Installed, 230 AC 50/60 Hz
2170281	CCL-240T-8-IVF	Celculture® Incubator, 240 L, IR sensor, CO ₂ & O ₂ control, Moist Heat Decon, with Sealed Inner Door Kit for 240 L (4/6 Glass Doors with Latches), Factory Installed, 230 VAC 50/60 Hz
2170261	CCL-050T-9-IVF	Celculture® Incubator, 50 L, IR sensor, CO ₂ & O ₂ control, Moist Heat Decon, with Sealed Inner Door Kit for 50 L (2 Glass Doors with Latches), Factory Installed, 115 VAC, 50/60 Hz
2170276	CCL-170T-9-IVF	CelCulture® Incubator 170 L IR Sensor, CO ₂ & O ₂ Control ULPA, Moist Heat Decon, with Sealed Inner Door Kit for 170 L (4 Glass Doors with Latches), Factory Installed, 115 VAC 50/60 Hz
2170282	CCL-240T-9-IVF	Celculture® Incubator, 240 L, IR sensor, CO ₂ & O ₂ control, Moist Heat Decon, with Sealed Inner Door Kit for 240 L (4/6 Glass Doors with Latches), Factory Installed, 115 VAC 50/60 Hz

Esco Multi-Zone ART Workstation



The Esco Multi-Zone ART Workstation is the most advanced workstation in its class. It is designed for use in applications that require a high level of control over environmental conditions. Applications can range from animal embryo culture in research to human embryo manipulation done in fertility laboratories.


Key Features



Multi-Zone Heating System

1 set point, 10 independent zones** with their own heating elements and sensors allow excellent uniformity. The heating system will automatically prioritize power distribution to ensure effective temperature control with fast recovery. In models that feature MIRI® chambers, temperature regulation of the chambers will always be on regardless of mode.

- Accuracy: $\pm 0.2\text{ }^{\circ}\text{C}$
- Uniformity: $\pm 0.2\text{ }^{\circ}\text{C}$

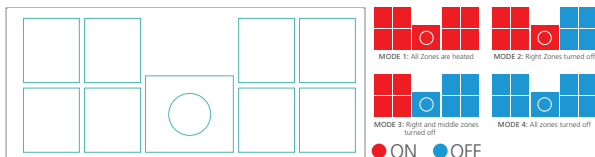

*When any of the heating zones are OFF, the monitor shall not display real-time temperature as there are no controlled heating to give uniformity across the OFF zone.

**The provided example of heating zone overlay is applicable for MAW-4D_model. For different heating zone layouts, please read pages 6 to 10 and pages 13 to 16. The heating modes are the same for all the models.



Humidification System

It allows digital control of premixed Gas which enables complete environmental control for the samples. Gas outlet port is located on the table surface and the plastic cover will enclose the humidified gas which creates small "incubator".



Heated Glass Stage

The heated glass stage has its own independent heating zone further to enhance the temperature control and recovery in this zone.

Stainless Steel Tabletop

The main material used in the tabletop surface is stainless steel. while aluminum is used for covering the bottom of the tabletop.



Microscope Integration Provision

The integrated stereoscope in the work chamber allows users to maintain culture dishes at steady temperature during observation and manipulation. Fewer movements will also reduce risk of accidents.



Surveillance System*

Provides the user with real-time information of zone performance and other work area parameters such as gas pressure and flow rate.

*When any of the heating zones are OFF, the monitor shall not display real-time temperature as there are no controlled heating to give uniformity across the OFF zone.

Esco Multi-Zone Workstation with MIRI® Chambers



The workstation is now even better with its integrated MIRI® chambers to further secure your embryos while inside the workstation. The MIRI® incubator is popular for its top-notch features such as stable culture environment and faster parameter recovery.

Multi-Zone Heating System

The independent zones have its own heating elements and sensors.

MIRI® Chambers

Your specimens are more secured than ever with the integrated MIRI® chambers, known for its stable and precise temperature output.



Support Stand Options

More options to choose from to meet your requirements.



Fit in the same MIRI® Inserts



Nunc™



Falcon®



Vitrolife



LifeGlobal®
GPS Dishes



Nipro™



SparMED Oosafe®

Accessories




Support Stand:


- With leveling feet
- With Caster Wheels
- Motorized stand with Caster Wheels



UV Kit

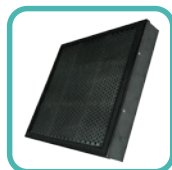

Carbon Pre-filter

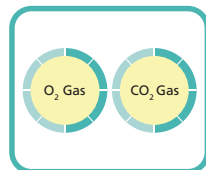

Gas Mixer for MAW

Stereomicroscope



A Stereomicroscope is a type of optical microscope that allows the user to see a three-dimensional view of a specimen. This Esco MS-I Stereomicroscope is very useful in assessing oocyte and embryo morphology and viability. It has a built-in camera mount where you can place your camera set for a better viewing and documentation of images. Moreover, additional objective lenses (0,3x, 0,4x, 0,75x) can be screwed on the microscope head for additional magnification.



Key Features

PK-Mount and Camera Set

Allows mounting a camera set that will help the user to have a better look on their slide.

Adjustable Interpupillary Distance

Allows adjustment of the distance between two eyepieces depending on their preference.

Black and White Work Board

This enhances the user's viewing of specimen/samples.



Up Light Source

This is one of the adaptations in an inverted designed microscope that

is useful in an IVF laboratory as it is designed for inspection of gametes and embryos.

The Esco Steremicroscope and Esco IVF Workstations

The Esco Multi-Zone ART Workstations (MAW) are IVF workstations with controlled temperature and humidity working environment. These workstations are intended for working with gametes and / or embryos at / or near body temperature in a multi-zone heated set-up during in vitro fertilization (IVF) / assisted reproductive technology (ART) treatments.

The Esco Multi-Zone ART Workstations are even better with its integrated MIRI® chambers to further secure the embryos while inside the workstation. The MIRI® incubator is popular for its top-notch features such as stable culture environment and faster parameter recovery. Also, these workstations have a provision for one's microscope of choice like the Esco Stereomicroscope.






Esco Multi-Zone ART Workstations



Available in a variety of sizes and configurations to meet the needs of the laboratory


MAW-3D8

Width: 3 Feet
Microscope: Single
Basic Configuration
1 user
For small laboratories

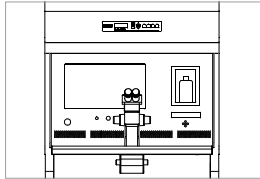

MAW-6D_ MONO


Width: 6ft
Microscope: Single
1 user
More space for other work.

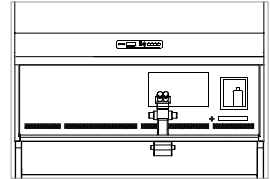

(Front View)

MAW-4D_


Width: 4ft
Microscope: Single
Basic Configuration
1 user
For small Laboratories

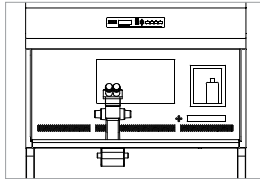

(Front View)

MAW-6D_-MP


Width: 6ft
Microscope: Single
Stereomicroscope,
1 Inverted microscope set-up

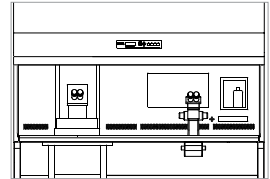

(Front View)

MAW-6D_-DUAL

Width: 6ft
Microscope: Dual
2 users
For efficient use of space

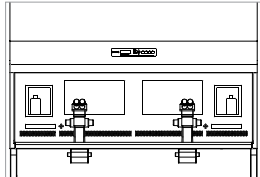


(Front View)




MAW-4D_-MC

Width: 4ft
Microscope: Single
1 user
For small Laboratories
Basic Configuration with
2 MIRI® Chambers

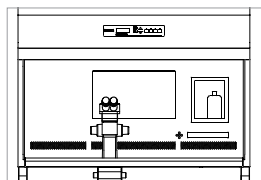


(Front View)

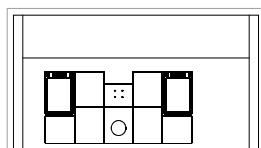


(Worktop View)

MAW-6D_- DUAL-MC

Width: 6ft
Microscope: Dual
2 users
Basic Configuration with 3
MIRI® Chambers
For efficient use of space

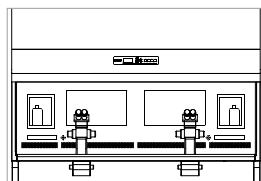


(Front View)

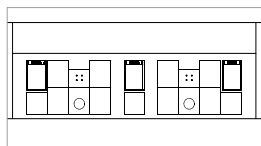


(Worktop View)

MAW-6D_- MONO-MC

Width: 6ft
Microscope: Single
1 user
More space for other work
Basic Configuration with 2
MIRI® Chambers

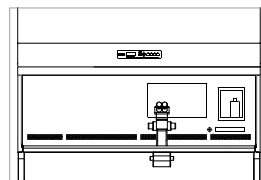


(Front View)

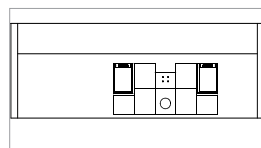


(Worktop View)

MAW-6D_-MP- MC

Width: 6ft
Microscope: Single
Stereomicroscope,
1 Inverted
microscope set-up
2 users
Multi-purpose
Embryo inspection &
micromanipulation
Basic Configuration with 2
MIRI® Chambers

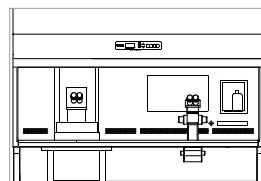

(Front View)

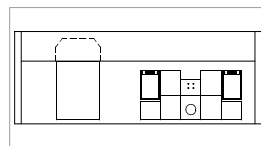

(Worktop View)



General Specifications



Esco Multi-Zone ART Workstations

Model	MAW-3D8	MAW-4D_	MAW-6D_-MONO	MAW-6D_-DUAL	MAW-6D_-MP
Nominal Size	0.9 meter (3')	1.2 meter (4")	1.8 meter (6")		
Internal Work area dimensions (W x D x H)*	950 x 500 x 710 mm (37.4" x 19.7" x 28.0")	1260 x 500 x 710 mm (49.6" x 19.7" x 28.0")	1870 x 500 x 710 mm (73.6" x 19.7" x 28.0")		1870 x 490 x 780 mm (73.6" x 19.3" x 30.7")
Laminar Air Velocity	Average of 0.21m/s or 41 fpm (± 20%)				
Filter Efficiency	>99.999% for particle size between 0.1 to 0.3 microns per IEST-RP-CC001.3 / H14 per EN 1822				
Pre-filter	Disposable and non-washable polyester fibers with 85% arrestance / EU3 rated				
Noise level (per NSF 49)**	<48 dBa	47 dBa	52 dBa		
Set of (9+1) Heating Zones	(7+1) Heating Zone	1 set	1 set	2 sets	1 set
Surveillance System with data logger	1 set	1 set	1 set	2 sets	1 set
Humidification System***	1 set	1 set	1 set	2 sets	1 set
PT 1000 ports	5 ports	5 ports	5 ports	10 ports	5 ports
Microscope	Required, Not included (see microscope ordering information)				
	Position for 1 stereo microscope			Position for 2 stereo microscopes	Position for 1 stereo microscope and 1 inverted microscope
Transmitted Light Source	1 set	1 set	1 set	2 sets	1 set
Shipping weight	195 kg	250 kg	350 kg		

* The actual MAW-6D_-MP's work area dimension will be customized to fit the inverted microscope.

** Noise reading in open field condition / anechoic chamber. Noise reading in normal room varies by room size, layout, and background noise, but may reach roughly 3-4 dBA above these values.

***1 set of humidification system includes one (1) water bottle with tubing, one (1) sample carry tray and one (1) plastic cover.

Esco Multi-Zone ART Workstation with MIRI® Chambers

Model	MAW-4D_-MC	MAW-6D_-MONO-MC	MAW6D_-DUAL-MC	MAW-6D_-MP-MC
Nominal Size	1.2 meter (4')	1.8 meter (6')	1.8 meter (6')	1.8 meter (6')
Internal Work area dimensions (W x D x H)*	1260 x 500 x 710 mm (49.6" x 19.7" x 28.0")	1870 x 500 x 710 mm (73.6" x 19.7" x 28.0")	1870 x 500 x 710 mm (73.6" x 19.7" x 28.0")	1870 x 490 x 780 mm (73.6" x 19.3" x 30.7")
External dimensions w/o support stand (W x D x H)	1340 x 640 x 1400 mm (52.8" x 25.2" x 55.1")	1950 x 640 x 1400 mm (76.8" x 25.2" x 55.1")	1950 x 640 x 1400 mm (76.8" x 25.2" x 55.1")	1950 x 640 x 1460 mm (76.8" x 25.2" x 57.5")
External Dimensions w/ SML-8 type support stand (W x D x H)	1340 x 640 x 2160 mm (52.8" x 25.2" x 85.0")	1950 x 640 x 2160 mm (76.8" x 25.2" x 85.0")	1950 x 640 x 2160 mm (76.8" x 25.2" x 85.0")	1950 x 640 x 2220 mm (76.8" x 25.2" x 87.4")
Laminar Air Velocity	Average of 0.21m/s or 41 fpm (± 20%)			
ULPA Filter Efficiency	>99.999% for particle size between 0.1 to 0.3 microns per IEST-RP-CC001.3 / H14 per EN 1822			
Pre-Filter	Disposable, non-washable polyester fibre, 85% arrestance, EU3 rated			
Set of (7+1) Heating Zones	1	1	2	1
Number of MIRI® Chambers	2	2	3	2
Included Accessories	1 x water bottle for HS-1, including tubing 2 x sample carry tray	1 x water bottle for HS-1, including tubing 2 x sample carry tray	2 x water bottle for HS-1, including tubing 3 x sample carry tray	1 x water bottle for HS-1, including tubing 2 x sample carry tray
Microscope	Required, Not included (see microscope ordering information)			
	Position for 1 stereo microscope	Position for 1 stereo microscope	Position for 2 stereo microscopes	Position for 1 microscope and 1 inverted microscope
Support stand	Required, Not included (see support stand ordering information)			

Stereomicroscope

Overall Dimensions (W x D x H)	250 x 250 x 240 mm (9.8" x 9.8" x 9.4")
Net Weight	4.5kg
Type of Light Source	Makes use of a light source integrated into a laminar table Note that an external LED ring could be used, mounted on a microscope (voltage is 110 – 220 V)
Camera Compatibility	Microscope can come with or without camera. Most important camera needs to be C-mount- ed. The recommended camera is IDS model UI-3880LE-C-HQ.

What we offer: IDS (from Germany) Camera USB 3.1 Gen 1, uEye, LE, 1/1.8" Sony IMX178 CMOS Color Sensor, 3088x2076 Pixels, Rolling Shutter, USB Type-C Connector, CS/C-Mount, HQ-IR Filter, Dimension 47 mm x 46 mm.

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
Esco Multi-Zone ART Workstation Basic Configuration and Multi-Purpose		
2070177	MAW-3D8	Esco Multi-Zone ART Workstation, 3ft 3ft (0.91m), 230V 50/60Hz
2070017	MAW-4D8	Esco Multi-Zone ART Workstation, 4ft (1.2m), 230V 50/60Hz
2070025	MAW-4D9	Esco Multi-Zone ART Workstation, 4ft (1.2m), 115V 50/60Hz
2070018	MAW-6D8-MONO	Esco Multi-Zone ART Workstation, 6ft (1.8m), 230V 50/60Hz
2070026	MAW-6D9-MONO	Esco Multi-Zone ART Workstation, 6ft (1.8m), 115V 50/60Hz
2070050	MAW-6D8-DUAL	Esco Multi-Zone ART Workstation, Double Heated Zone, 6ft (1.8m), 230V 50/60Hz
2070039	MAW-6D9-DUAL	Esco Multi-Zone ART Workstation, Double Heated Zone, 6ft (1.8m), 115V 50/60Hz
2070036	MAW-6D8-MP	Esco Multi-Zone ART Workstation, Multi-Purpose, 6ft (1.8m), 230V 50/60Hz
2070038	MAW-6D9-MP	Esco Multi-Zone ART Workstation, Multi-Purpose, 6ft (1.8m), 115V 50/60Hz
Esco ART Workstation with MIRI® Chambers		
2070068	MAW-4D8-MC	Esco Multi-Zone ART Workstation, MIRI® Chambers, 4ft (1.2m), 230V 50/60Hz
2070070	MAW-4D9-MC	Esco Multi-Zone ART Workstation, MIRI® Chambers, 4ft (1.2m), 115V 50/60Hz
2070069	MAW-6D8-MONO-MC	Esco Multi-Zone ART Workstation, MIRI® Chambers, 6ft (1.8m), 230V 50/60Hz
2070071	MAW-6D9-MONO-MC	Esco Multi-Zone ART Workstation, MIRI® Chambers, 6ft (1.8m), 115V 50/60Hz
2070075	MAW-6D8-DUAL-MC	Esco Multi-Zone ART Workstation, Double Heated Zone, MIRI® Chambers, 6ft (1.8m), 230V 50/60Hz
2070074	MAW-6D9-DUAL-MC	Esco Multi-Zone ART Workstation, Double Heated Zone, MIRI® Chambers, 6ft (1.8m), 115V 50/60Hz
2070072	MAW-6D8-MP-MC	Esco Multi-Zone ART Workstation, Multi-Purpose, MIRI® Chambers, 6ft (1.8m), 230V 50/60Hz
2070073	MAW-6D9-MP-MC	Esco Multi-Zone ART Workstation, Multi-Purpose, MIRI® Chambers, 6ft (1.8m), 115V 50/60Hz
ART Workstation with MIRI® Chambers and Trigas Mixer		
2070187	MAW-4D8-MC-G	Esco Multi-Zone ART Workstation, MIRI® Chambers, Gas Mixer, 4ft (1.2m), 230V 50/60Hz
2070188	MAW-4D9-MC-G	Esco Multi-Zone ART Workstation, MIRI® Chambers, Gas Mixer, 4ft (1.2m), 115V 50/60Hz
2070189	MAW-6D8-MONO-MC-G	Esco Multi-Zone ART Workstation, MIRI® Chambers, Gas Mixer, 6ft (1.8m), 230V 50/60Hz
2070190	MAW-6D9-MONO-MC-G	Esco Multi-Zone ART Workstation, MIRI® Chambers, Gas Mixer, 6ft (1.8m), 115V 50/60Hz
2070191	MAW-6D8-DUAL-MC-G	Esco Multi-Zone ART Workstation, Double Heated Zone, MIRI® Chambers, Gas Mixer, 6ft (1.8m), 230V 50/60Hz
2070192	MAW-6D9-DUAL-MC-G	Esco Multi-Zone ART Workstation, Double Heated Zone, MIRI® Chambers, Gas Mixer, 6ft (1.8m), 115V 50/60Hz
2070193	MAW-6D8-MP-MC-G	Esco Multi-Zone ART Workstation, Multi-Purpose, MIRI® Chambers, Gas Mixer, 6ft (1.8m), 230V 50/60Hz
2070194	MAW-6D9-MP-MC-G	Esco Multi-Zone ART Workstation, Multi-Purpose, MIRI® Chambers, Gas Mixer, 6ft (1.8m), 115V 50/60Hz
Stereomicroscope		
1320375	MS – I	Stereomicroscope

MIRI® Anti-Vibration Table



The MIRI® AVT (Anti-Vibration Table) features an anti-vibration mechanism for passive dampening of the microscope. This is mainly used for micromanipulation procedures like Intra-Cytoplasmic Sperm Injection (ICSI) procedures. Exclusively designed in Denmark and made in E.U., the stainless steel table and sturdy frame add mass to the anti-vibration table. MIRI® AVT is constructed to be easy-to-use and almost maintenance-free.

Key Features

- Anti-vibration mechanism for passive dampening
- Sturdy frame
- Stainless steel table
- Range (HZ) vibration could be eliminated while using this AVT: 5.5-50Hz

General Specifications

Overall Dimensions	1200 x 800 x 785 mm (47.2 x 31.5 x 31")
Net weight	117 kg (258 lbs)
Material	Powder-painted mild steel, Stainless steel
Float Dimension (W x D)	540 x 340 mm (21.3 x 13.4")

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
1320484	MIRI® AVT	Anti-Vibration Table MIRI® AVT



Versati™ Tabletop Centrifuge



Versati™ Tabletop Centrifuge Overview



Versati™ Tabletop centrifuge stands out among the same-level products with its versatility, running features, and easy handling. It can be used with high-capacity and low-to-high-speed general-purpose centrifuge applications. It is suitable for sperm preparation for human and animal IVF

Key Features

- Compact Design
- Incredible Flexibility
- Overspeed Protection

Model: TCV-1500

- Maintenance-free brushless motor
- Superior safety (Automatic rotor recognition)
- Audible and visible alarms
- Up to 1500 ml capacity

Accessories



Fixed-angle Rotor

- Mainly for hospital use in centrifugation of urine and blood
- Rotor and adapters are autoclavable (maximum of 121°C, 20 minutes).
- Can be used for Sperm Preparation for Human and Animal IVF

Ordering Information

ITEM CODE	MODEL CODE	DESCRIPTION
2220005	TCV-1500-8	Tabletop Centrifuge Ventilated 230 VAC, 50/60 Hz
2220006	TCV-1500-9	Tabletop Centrifuge Ventilated 120 VAC, 50/60 Hz

Guide to Models

TCV - 1500 - 8

Model	Code	Centrifuge Capacity	Code	Electrical Supply	Code
Tabletop Centrifuge Ventilated	TCV	Maximum Sample Capacity (ml)	1500	230 VAC 50/60 Hz	8
				120 VAC 50/60 Hz	9
Tabletop Centrifuge Refrigerated	TCR			230 VAC 50/60 Hz	8
				120 VAC 50/60 Hz	9

TCV / TCR High Speed Tabletop Centrifuge for up to 16,000 rpm

Ideal for multi-purpose centrifugation at high-speed.

Model Code	Item Code	Description
TCV-1500-8	2220005	Tabletop Centrifuge Ventilated 230 VAC, 50/60 Hz
TCV-1500-9	2220006	Tabletop Centrifuge Ventilated 120 VAC, 50/60 Hz

Contact Esco or your local Sales Representative for ordering information

Quality Assurance and Validation Units



MIRI® GA Gas and Temperature Validation Unit

MIRI® GA is a tabletop device intended to make external incubator validation easier and safer. It is capable of monitoring the temperature (PT1000 connector) & gas concentration, flow and pressure. It can validate up to 6 chambers simultaneously 24 hours a day. It also has an adjustable flow rate which gives it the ability to properly sample small volume incubation chambers. Moreover, MIRI® GA comes with a full Data Logger software which is helpful in monitoring each parameter. The MIRI® GA can connect to any brand of incubator and is a perfect accessory to MIRI® TL, MIRI® Multiroom Incubator, and Mini MIRI® Dry Incubators.

Key Features

- Constantly validate up to 6 x CO₂ / O₂ incubators
- CO₂ / O₂ incubators controllable flow rate Monitor up to 6 x
- PT1000 sensors
- 6 ports for sequential gas samples
- Gas feedback returns sampled gas to incubator or exhaust

General Specifications

Input ports	6 x PT1000 ports for temperature monitoring 6 x gas sampling ports
Output ports	1 x gas feedback port, 1 x USB port
Shipping dimensions and weight	460 x 450 x 250 mm (18.1" x 17.7" x 9.8"), 10 kg (22 lbs)

Ordering Information*

ITEM CODE	MODEL CODE	DESCRIPTION
1320045	MRI-GA	MIRI® GA CO ₂ / O ₂ / Temp validation Unit, 115/230V, 50/60Hz

*Includes data logger software, 1pc PT1000 cable, 1pc Gas connection tube, 1pc Gas feedback tube

Accessories

ITEM CODE	MODEL CODE	DESCRIPTION
1320063	MRA-1101	1pc PT1000 cable
1320064	MRA-1102	Set of 6pcs PT1000 cables
1320065	MRA-1103	1pc Gas connection tube
1320066	MRA-1104	Set of 6pcs Gas connection tubes



Continuous pH Monitoring

Worry-free pH monitoring of your Esco incubators

SAFE Sens (Sterile, Automated Fluoroscopic Evaluation) offers fast, effective, and non-invasive continuous pH monitoring product for in vitro fertilization (IVF) procedures.

The SAFE Sens technology employs an optical fluorescent measurement technology, used in combination with disposable sensors, which accurately and reliably monitors the pH of small volumes of fluids such as the media used in IVF.

Order SAFE Sens with your new MIRI®, Mini MIRI®, MIRI® TL*

Key Features



Continuous pH measurement

- Reading and recording every 30 minutes (default setting - adjustable).
- Single use sensor probe for up to seven (7) days of pH readings.



Easy-to-implement

- Easy-to-align (no buffers, no hassles).
- Easy-to-use and maintain.




Data Logging System

- Data Logging and user alarms.
- Each TrakStation® can be connected to multiple incubators.



Compact and Efficient

No more unnecessary openings of your incubator for spot pH measurement.

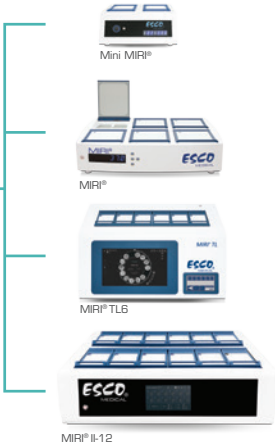

*Factory-installed. No option yet to install on-site for previously ordered MIRI®.

**SAFE Sens is a trademark brand of Blood Cell Storage, Inc. (BCSI)

Ordering Guide for SAFE Sens

Step1: Choose the correct electrical rating for each equipment.

ITEM CODE	MODEL CODE	DESCRIPTION
MIRI® TL6		
2070098	MRI-TL-MN-6C-SS-8	MIRI® Time-Lapse Incubator, Mini, 6 chambers with SAFE Sens, 230 V, 50/60 Hz
2070099	MRI-TL-MN-6C-SS-9	MIRI® Time-Lapse Incubator, Mini, 6 chambers with SAFE Sens, 115 V, 50/60 Hz
MIRI® TL12		
2070100	MRI-TL-12C-8	MIRI® Time-Lapse Incubator, 12 Chambers, 230 V, 50/60 Hz
2070101	MRI-TL-12C-9	MIRI® Time-Lapse Incubator, 12 Chambers, 115 V, 50/60 Hz
MIRI®		
2070086	MRI-6A10-SS-8	MIRI® Incubator, with SAFE Sens for pH measurement, 230V, 50/60Hz
2070087	MRI-6A10-SS-9	MIRI® Incubator, with SAFE Sens for pH measurement, 115V, 50/60Hz
MIRI® Humidity		
2070185	MRI-6A10-H-SS-8	MIRI® Humidity 6 chamber with SAFESens for pH monitoring, 230V
2070186	MRI-6A10-H-SS-9	MIRI® Humidity 6 chamber with SAFESens for pH monitoring, 115V
Mini MIRI® Dry		
2070145	MRI-MINI-SS-D-8	Mini MIRI® Dry with SAFE Sens for pH monitoring, 230V, 50/60 Hz
2070146	MRI-MINI-SS-D-9	Mini MIRI® Dry with SAFE Sens for pH monitoring, 115V, 50/60 Hz
Mini MIRI® Humidity		
2070157	MRI-MINI-H-SS-8	Mini MIRI® Humidity with SAFE Sens for pH monitoring, 230V, 50/60 Hz
2070158	MRI-MINI-H-SS-9	Mini MIRI® Humidity with SAFE Sens for pH monitoring, 115V, 50/60 Hz
MIRI® II-12		
2070166	MRI2-12C-SS-8	MIRI® II-12 Incubator with 12 chambers, with SAFE Sens, 230V, 50/60Hz
2070167	MRI2-12C-SS-9	MIRI® II-12 Incubator with 12 chambers, with SAFE Sens, 115V, 50/60Hz

Step 2: Order the SAFE Sens Accessories.*

ITEM CODE	MODEL CODE	DESCRIPTION
1081277	MRA-SS-SV2	SAFE Sens SV2 Sensor, Pack of 10 pieces (shelf-life 12 months)
1081278	MRA-SS-QC2	
1320191	MRA-SS-TS	SAFE Sens TrakStation, a tablet with SAFE Sens Software, for pH monitoring.

***Notes:** (1) One QC2 alignment tool can be used on all incubators. If incubators are located in separate rooms, you may have to order more than one QC2 tool.
(2) QC2 Alignment tool and SV2 sensors have an expiration date of one (1) year.
(3) The MIRI® with SAFE Sens automatically comes with free one (1) pack of SV2 sensors, which is to be used for Site Standardization. Please determine how many additional packs you need for routine pH testing . Put into consideration that the sv2 sensor has one year expiration date. *See page 5 for stock-keeping matrix*
(4) One TrakStation can connect up to eight (8) incubators by using a USB 3.0 Hub. Determine how many PC tablets you need.
(5) The Esco data logger should be installed separately on another Windows PC/ tablet. Please see page 5 for the minimum requirements.

Step 3: Choose Heat Optimization Plate. Only one chamber has the SAFE Sens installed so order only one (1) plate that has a hole for the SAFE Sens. See guide below.

No Heat Optimization Plate needed for MIRI® TL



Note: Heat plates should only be ordered for MIRI® and Mini MIRI®.

Heating Plates with Hole for SAFE Sens

ITEM CODE	MODEL CODE	DESCRIPTION
For MIRI® and Mini MIRI®		
1320220	MRA-ND-SS	Insert for Nunc™ Dishes, with hole for SAFE Sens
1320221	MRA-VD-SS	Insert for Vitrolife Dishes, with hole for SAFE Sens
1320222	MRA-NID-SS	Insert for Nipro™ Dishes, with hole for SAFE Sens
1320223	MRA-LD-SS	Insert for LifeGlobal® GPS Dishes, with hole for SAFE Sens
1320224	MRA-PD-SS	Insert Without Footprint for Plain Dishes, with hole for SAFE Sens
1320225	MRA-OD-SS	Insert for SparMED Oosafe®, with hole for SAFE Sens
1320219	MRA-FD-SS	Insert for Falcon® Dishes, with hole for SAFE Sens
For MIRI® II-12		
1320417	MRA2-FD-SS	Insert for Falcon® Dishes, with hole for SAFE Sens
1320418	MRA2-ND-SS	Insert for Nunc™ Dishes, with hole for SAFE Sens
1320419	MRA2-VD-SS	Insert for Vitrolife Dishes, with hole for SAFE Sens
1320421	MRA2-LD-SS	Insert for LifeGlobal® GPS Dishes, with hole for SAFE Sens
1320437	MRA2-OD-SS	Insert for SparMED Oosafe®, with hole for SAFE Sens
1320422	MRA2-PD-SS	Insert for Plain Dish with hole for SAFE Sens

Traceability Tool



The ultimate traceability tool for fertility clinics

MIRI® Evidence is redefining the way traceability in laboratory procedures is managed and stored in fertility clinics, eliminating errors, preventing system mix-ups, saving time, and helping clinics to comply with regulations. This tool provides a single, secure platform that ensures all patients are scheduled on time and documented properly - saving time and improving outcomes across the board.

The MIRI® Evidence workflow proposes tools to improve scheduling and documentation procedures, refining both clinical and administrative processes in a fertility clinic. Tasks or procedures are spontaneously arranged and scheduled as soon as the patient starts her treatment. Any changes to her treatment initiate an update of the scheduled task. The built-in validation engine warrants that all data is registered. If data is missing, a task is automatically created for the person or group of persons responsible for maintaining such information.

Benefits of Using the MIRI® Evidence



• Tracking

The scanners used with MIRI® Evidence are 2D imagers with no laser from a market-leading supplier and specially made for healthcare. The scanner can be used to scan wristbands, capture images, and photograph ID cards and passports.



• Workflow

The MIRI® Evidence workflow offers tools to enhance scheduling and documentation procedures, improving both clinical and administrative workflows in a fertility clinic. Once a patient's workflow has been started, the system automatically allocates the tasks. The combination of events determines the type of treatment.



• Compliance


MIRI® Evidence Tracking has been made to comply with the EDQM guidelines and EU regulations. The workflow and tracking system assists the fertility clinic in documenting all procedures in the laboratory and office, thus excluding double-witnessing and saving time.



• Mix-up prevention

Oocytes, embryos, and semen samples are tracked and validated in and out of the containers eliminating the need for human double-witnessing. The validation guarantees that the correct sample is being processed and prevents patient samples from being mixed up.

ESCO LIFESCIENCES GROUP

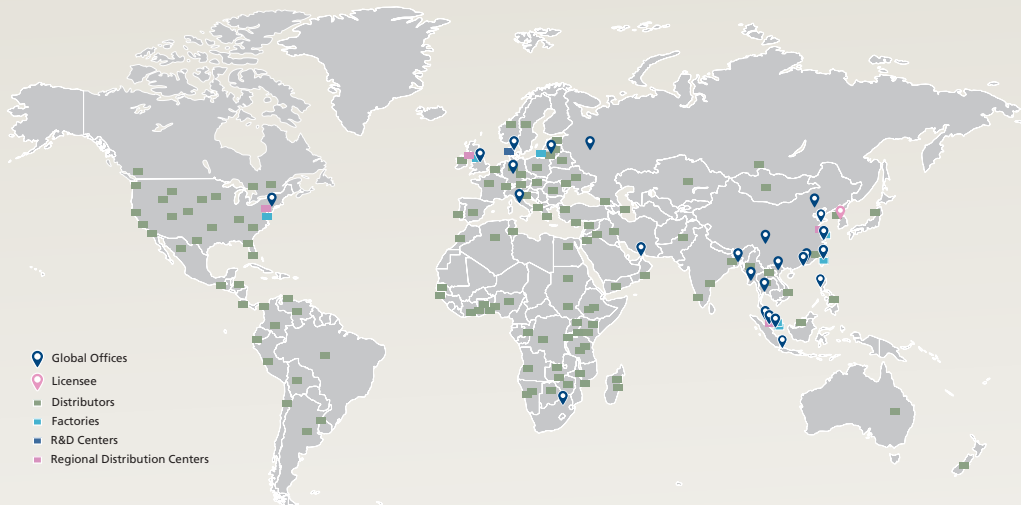

Esco Medical Products

MIRI® TL6 Time-Lapse Incubator
MIRI® TL12 Time-Lapse Incubator
MIRI® Multiroom Incubator
MIRI® Humidity Multiroom Incubator
MIRI® II-12 Multiroom Incubator
Mini MIRI® Dry Incubator
Mini MIRI® Humidity Incubator
Esco Multi-Zone ART Workstation
MIRI® Laminar Flow

MIRI® Evidence RFID Witnessing System
CelCulture® CO₂ Incubator
MIRI® GA (Gas and Temperature Validation Unit)
MIRI® AVT
Versati™ Tabletop Centrifuge
CultureCoin®

Infertility is a problem that has a significant social, psychological, and economic impact on afflicted individuals and couples. It is a global concern that knows no race or creed. It has been estimated that 1 in 6 couples struggle with infertility at least once in their lifetime.

Esco Medical is one of the divisions of the Esco Lifesciences Group. We provide innovative technological solutions for fertility clinics and laboratories. We aim to become the leading manufacturer of high-quality equipment such as long-term embryo incubators, ART workstations, anti-vibration tables, and time-lapse incubators.

Our products are designed with the Silent Embryo Hypothesis as a guiding principle. The Silent Embryo Hypothesis states that the less disturbed an embryo can remain, the better its developmental potential will be. Most of our products are designed in Denmark and made in the EU. Our primary focus is to increase pregnancy success rates and patient satisfaction.

ESCO MEDICAL

19 Changi South Street 1, Singapore 486779
Tel +65 6542 0833 • medical@escolifesciences.com
www.esco-medical.com

Esco Global Offices: Bangladesh | China | Denmark | Germany | Hong Kong | Indonesia | Lithuania | Malaysia | Myanmar | Philippines | Russia | Singapore | South Africa | South Korea | Taiwan | Thailand | UAE | UK | USA | Vietnam



Designed in Denmark



Made in the E.U.

